

AM 502
DIFFERENTIAL
AMPLIFIER

INSTRUCTION MANUAL

BEFORE READING

PLEASE CHECK FOR CHANGE INFORMATION AT THE REAR OF THIS MANUAL.

THIS MANUAL REPRINTED MARCH 1978

TEKTRONIX®

AM 502
DIFFERENTIAL
AMPLIFIER

INSTRUCTION MANUAL

Tektronix, Inc. P.O. Box 500 Beaverton, Oregon 97077

Serial Number

WARRANTY

This Tektronix instrument is warranted against defective materials and workmanship for one year. Any questions with respect to the warranty should be taken up with your Tektronix Field Engineer or representative.

All requests for repairs and replacement parts should be directed to the Tektronix Field Office or representative in your area. This will assure you the fastest possible service. Please include the instrument type number or part number and serial number with all requests for parts or service.

Copyright © 1973, 1978 by Tektronix, Inc. All rights reserved. Contents of this publication may not be reproduced in any form without the permission of Tektronix, Inc.

Products of Tektronix, Inc. and its subsidiaries are covered by U.S. and foreign patents and/or pending patents.

TEKTRONIX, TEK, SCOPE-MOBILE, TELEQUIPMENT, and are registered trademarks of Tektronix, Inc. Printed in U.S.A. Specification and price change privileges are reserved.

TABLE OF CONTENTS


		Pag
SECTION 1	OPERATING INSTRUCTIONS	1-1
SECTION 2	SPECIFICATION AND PERFORMANCE CHECK	2-1

WARNING

The remaining portion of this Table of Contents lists servicing instructions that expose personnel to hazardous voltages. These instructions are for qualified service personnel only.

SECTION 3	ADJUSTMENT	3-1
SECTION 4	MAINTENANCE AND INTERFACING INFORMATION	4-1
SECTION 5	CIRCUIT DESCRIPTION	5-1
SECTION 6	OPTIONS	6-1
SECTION 7	REPLACEABLE ELECTRICAL PARTS	7-1
SECTION 8	DIAGRAM AND CIRCUIT BOARD ILLUSTRATION	8-1
SECTION 9	REPLACEABLE MECHANICAL PARTS AND EXPLODED VIEW	9-1

CHANGE INFORMATION

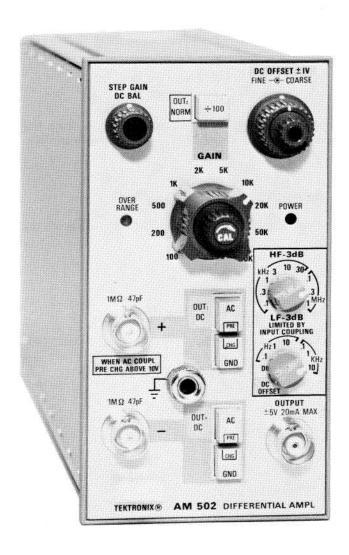

1582-1

AM 502 Differential Amplifier plug-in module.

OPERATING INSTRUCTIONS

INTRODUCTION

Description

The AM 502 is a dc-coupled differential amplifier with excellent common-mode rejection capabilities and high gain for low voltage measurements. The dc offset capability permits nulling up to 1 volt dc, so that low level, low frequency signals impressed on a dc level can be amplified without the degradation often introduced by ac input coupling. High and low frequency —3 dB points can be selected at the front panel to suit the application. Signal inputs and outputs are available at the rear connector as well as at the front panel. A front panel lamp indicates most overrange conditions of excessive input signal, excessive gain, or excessive offset.

The input circuit can be represented by approximately 1 megohm to ground paralleled by approximately 47 picofarads. The input impedance can be raised to about 200 megohms with the removal of a jumper.

Overload protection is provided by fuses in series with the input which open when large amounts of current flow due to overloading conditions.

A STEP GAIN DC BALANCE control is provided to adjust for zero shift as the gain switch is changed from one position to another.

The DC OFFSET COARSE and FINE controls give offset up to ± 1 volt dc potential difference at the input connectors. The amplifier internal bias is changed to accomplish the offset. The LF -3 dB switch must be in the DC OFFSET position to actuate these controls.

The HF -3 dB switch is used to reduce the upper bandwidth limit as necessary to improve the signal-to-noise ratio when using the AM 502 in low-frequency applications. The LF -3 dB control increases the lower bandwidth frequency.

Use of the pre-charging feature prevents surge currents, due to charging the ac coupling capacitor, from damaging the circuit under test.

Installation and Removal



Turn the power module off before inserting the plugin; otherwise, damage may occur to the plugin circuitry. Because of the high current drawn by the AM 502, it is also recommended that the power module be turned off before removing the AM 502. Refer to Fig. 1-1. Check to see that the plastic barriers on the interconnecting jack of the selected power module compartment match the cut-outs in the AM 502 circuit board edge connector.

Align the AM 502 chassis with the upper and lower guides of the selected compartment. Push the module in and press firmly to seat the circuit board in the interconnecting jack.

To remove the AM 502, pull on the release latch located in the lower left corner, until the interconnecting jack disengages and the AM 502 will slide out.

Controls and Connectors

Refer to Fig. 1-2. Even though the AM 502 is fully calibrated and ready to use, the functions and actions of the controls and connectors should be reviewed before attempting to use it. Pull the Power switch on the power module to apply power to the AM 502. The POWER indicator light indicates when power is applied to the AM 502.

OPERATING CONSIDERATIONS

Overheating

The AM 502 is designed to operate at an ambient temperature from 0° to +50°C. However, when operating several power supplies in a multi-plug-in power module, especially at low output voltages, or when operating close to other heat-producing equipment, internal temperature may exceed safe limits and actuate a thermal cutout in the power module. Refer to the power module instruction manual for more complete information.

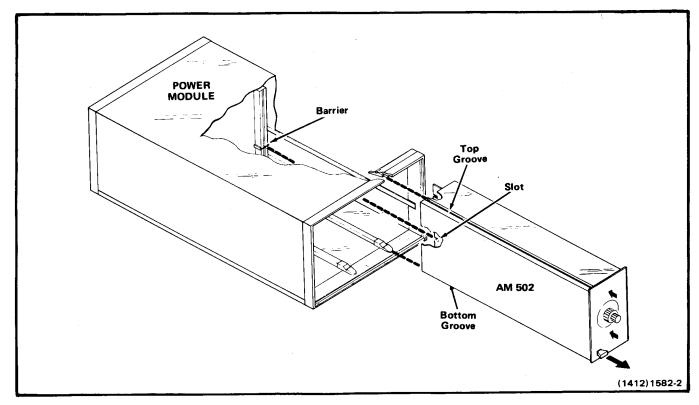


Fig. 1-1. Plug-in installation and removal.

Input Connections

Unshielded test leads can be used to connect the AM 502 to a signal source when a high-level, low-frequency signal is monitored at a low impedance point. However, when any of these factors is missing, it becomes increasingly important to use shielded signal cables. In all cases, the signal-transporting leads should be kept as short as practical.

When making single-ended input measurements (conventional amplifier operation), be sure to establish a common ground connection between the device under test and the AM 502. The shield of a coaxial cable is normally used for this purpose. See Fig. 1-3 for reference. Push the GND button for the input not connected to the device under test.

In some cases differential measurements require no common chassis ground connection, and therefore are less susceptible to interference by ground-loop currents.

Consider the change in the source operating characteristics due to loading by the signal input cables. The circuit at the input connectors can be represented by approximately 1 megohm to ground paralleled by approximately 47 picofarads. Two feet of 50 ohm coaxial cable increases the parallel capacitance by about 60 picofarads, which could be excessive in many situations. To minimize these effects, use a higher impedance cable or an attenuator probe.

Probes


Attenuator probes decrease the resistive-capacitive loading of a signal source. They also extend the measurement range of the AM 502 to substantially higher voltages. Some measurements require a higher resistance input to the AM 502, with very little source loading or signal attenuation. In such situations use a FET probe or the high-impedance input provision of the AM 502. Contact your Tektronix Representative for further information on probes.

High Impedance Input

To raise the internal input impedance of the AM 502 to about 200 megohms, remove the P40 plug (Fig. 3-1). Make certain the attenuator is in the NORM mode. Signal source impedance now becomes an important factor. For example, a 100 picoampere gate current through 10 megohms produces a one-millivolt offset. This offset may result in significant error when small voltages are measured.

Input Overloading

When measuring unknown dc voltages, push the $\div 100$ pushbutton in, and start with the 100 position on the GAIN switch. Increase the GAIN switch setting and finally release the $\div 100$ pushbutton until a suitable output signal is obtained. If the input circuit of the AM 502 is overdriven, large amounts of current will flow, opening the protective fuses.

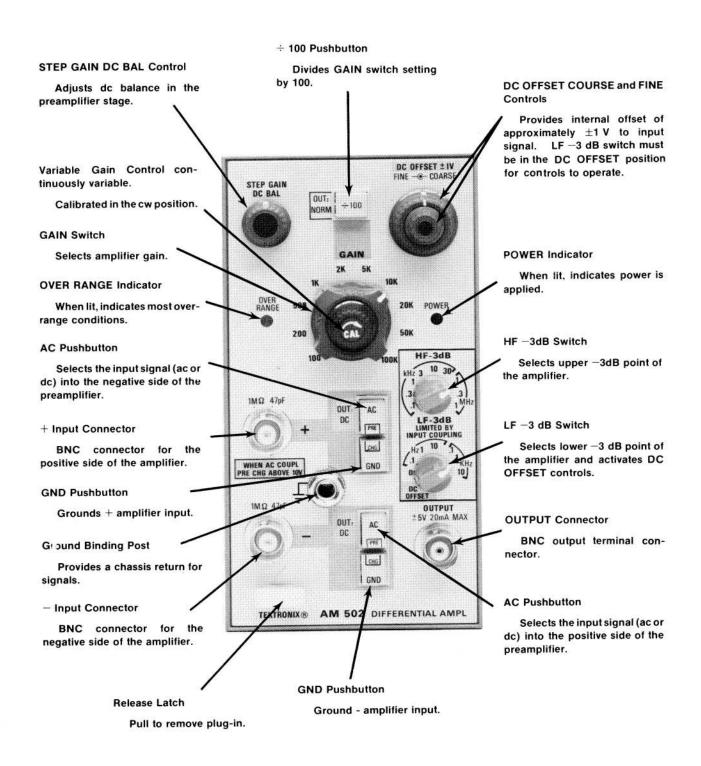

1582-9

Fig. 1-2. AM 502 controls and connectors.

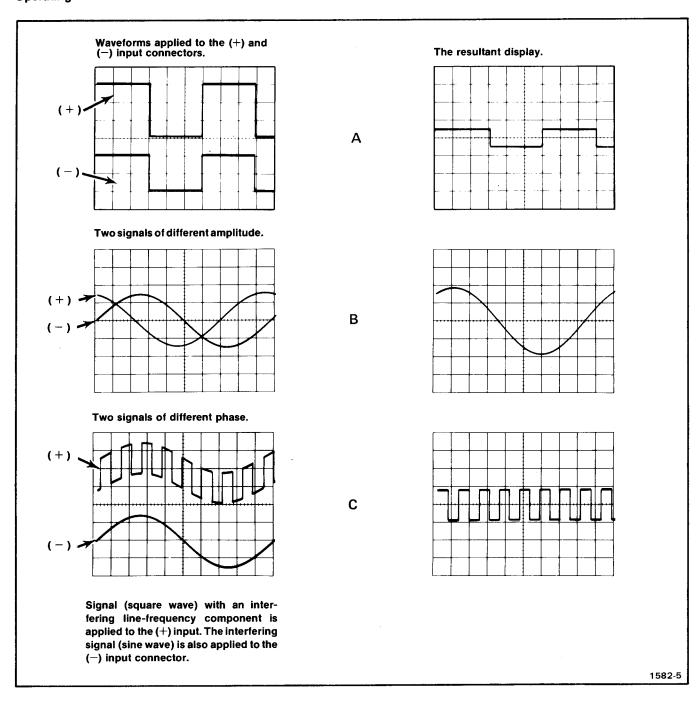


Fig. 1-3. Waveforms showing differential rejection of a common-mode signal. Resultant waveforms show the difference between the two signals.

Output Connections

Make output connections using a bnc to dual binding post connector, or a coaxial cable with at least one bnc connector. To prevent current limiting in the output stage, do not load the output with less than 250 ohms. Output current is limited to 20 milliamperes.

Step Gain DC Balance

If this control is misadjusted, the dc output level will shift as the GAIN switch position is changed. Push both GND buttons and place the GAIN switch in the 100 position. Rotate the GAIN switch from stop to stop while adjusting the STEP GAIN DC BAL control for no dc shift at the OUTPUT terminal.

DC Offset Coarse and Fine

Use these controls to offset up to ± 1 volt dc potential difference at the input connectors. The amplifier internal bias is changed to accomplish the offset. The differential rejection capabilities of the AM 502 are not affected. The LF -3dB switch must be in the DC OFFSET position to activate these controls.

HF and LF Bandwidth Reduction

Use the HF-3dB switch to reduce the upper bandwidth limit, as necessary, to improve the signal-to-noise ratio when using the AM 502 in low-frequency applications. The LF-3dB control increases the lower bandwidth frequency. Use this control to reduce dc drift, when raising the lower bandwidth does not undesirably reduce the bandwidth for the signal being measured.

Pre-Charging

Use of this feature prevents surge currents, due to charging the ac coupling capacitor in the AM 502, from damaging the circuit under test. Before connecting the AM 502 to a signal containing a dc component, push the AC and GND pushbuttons. Connect the input to the circuit under test. Wait about one second for the coupling capacitor to charge. Release the GND pushbutton, and the coupling capacitor is charged to the value of the dc voltage to be measured.

Differential Operation

A differential measurement is made by connecting each of the two inputs to selected points in the test circuit. The input to the amplifier will then be the difference in voltage of the two selected points. Consideration should be given to the proper connection method used between the AM 502 and the circuit under test; otherwise improper measurement results may occur. See Fig. 1-4 for reference.

Differential voltage measurements are made by applying the signals to the +input and -input connectors. Set the input coupling switches to the same position, AC or DC, depending on the measurement being made. In differential measurements, only the voltage difference between the two signals is amplified. Common mode signals (common in amplitude, phase, and frequency) are rejected. See Fig. 1-3 A, B, and C for reference.

Single-ended measurements often yield unsatisfactory results because of the interference resulting from ground-loop currents between the AM 502 and the device under test. In other cases, it may be desirable to eliminate a dc voltage by means other than the use of a blocking

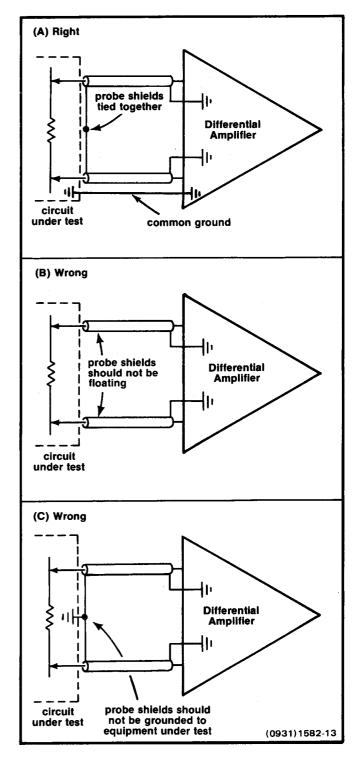


Fig. 1-4. Connecting a differential amplifier to a signal source.

capacitor, which could limit the low-frequency response. The limitations of single-ended measurements are effectively eliminated by using differential measurements.

DC Offset Operation

By using the FINE and COARSE DC OFFSET controls, it is possible to use the AM 502 differentially in a slide-back mode, to observe a small signal whose dc potential difference may be considerable. The offset voltage is continuously adjustable from plus 1 volt to minus 1 volt, and is internally available for all of the GAIN switch positions. The LF —3dB switch must be in the DC OFFSET position to activate the offset circuit.

Common Mode Rejection Ratio

The ability of the AM 502 to reject common-mode signals is indicated by the common-mode rejection ratio (cmrr). For example, assume that a signal consists of unwanted 60 hertz at 10 volts peak-to-peak (common mode connected to both inputs), plus a desired signal at 1 millivolt peak-to-peak (differentially connected to one input). The AM 502 gain is set at 200. The output of the AM 502 shows the desired signal at an amplitude of 0.2 volt (1 millivolt X 200), and the 60 hertz signal is viewed at an amplitude of 0.02 volts. The cmrr in this application is 100,000 to 1. This figure is calculated by multiplying the value of the common-mode signal (10 volt) by the gain of the amplifier (200) for a product of 2000 volts. This product is divided by the observed 60 hertz voltage at the output (0.02 volt) and result is the cmrr, 100,000 to 1. It would, of course be impossible to view the 1 millivolt signal superimposed on the 60 hertz signal by using single input methods.

Degradation of Common Mode Rejection

There are a number of factors that degrade common-mode rejection ratio (cmrr). The principal requirement for maximum rejection is for the common-mode signal to arrive at the input FET gates at the same phase and amplitude. A difference of only 0.01% in the attenuation ratios of the input attenuators will reduce the rejection ratio to 10,000 to 1. Also, any difference in source impedance will degrade the rejection ratio. Figures 1-5 and 1-6 show common-mode rejection degradation due to differences in source impedance. The frequency of the common-mode signal also affects the common-mode rejection ratio. Generally, as the frequency of the input signal increases, the common-mode rejection ratio is more difficult to maintain.

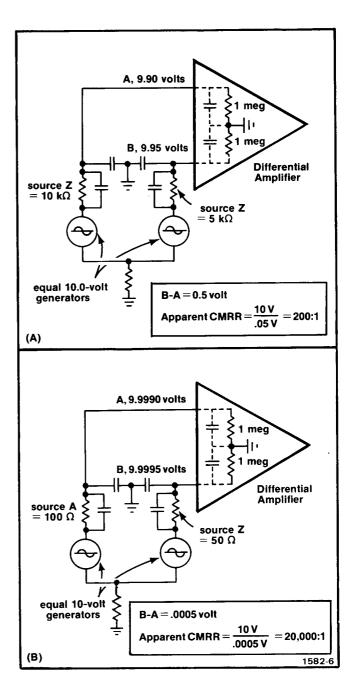


Fig. 1-5. Relationship of test point source impedance to the amplifier input impedance and the apparent CMRR caused by (A) large difference between test-point impedances and (B) low impedance test points.

1-6 REV. C MAR 1978

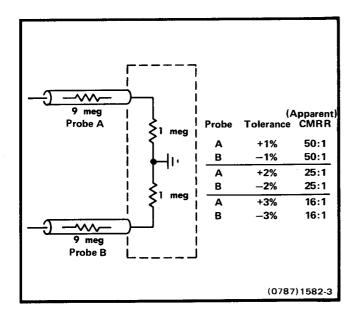


Fig. 1-6. Simplified input circuit and table showing the change in apparent CMRR due to 10X probes that are within 1, 2, and 3% of their attenuation values (with matched 1 megohm resistors).

The high frequency cmrr will also depend upon the signal source impedance, since various shunt capacitances between the source and the input gate must charge and discharge through that impedance.

Outside influences such as magnetic fields can also degrade the performance, particularly when low level signals are involved. Magnetic interference may be reduced by using identical signal transporting leads to the two inputs. Twist the two leads together over as much of their length as possible. Low-frequency measurements can be similarly protected by using a shielded cable that contains a twisted pair of conductors.

SPECIFICATION AND PERFORMANCE CHECK

SPECIFICATION

Performance Conditions

The electrical characteristics are valid only if the AM 502 has been calibrated at an ambient temperature between +20°C and +30°C and is operating at an ambient temperature between 0°C and +50°C unless otherwise noted

Items listed in the Performance Requirements column of the Electrical Characteristics are verified by completing the Performance Check in this manual. Items listed in the Supplemental Information column are not verified in this manual; they are either explanatory notes or performance characteristics for which no limits are specified.

Table 2-1
ELECTRICAL CHARACTERISTICS

Characteristics	Performance Requirements	Supplemental Information	
Gain			
Normal Mode	100 to 100,000 within 2%, 10 steps in 1, 2, 5 sequence.		
÷100 Mode	1 to 1000 within 2%.		
Variable Range		Continuously variable uncalibrated gain to at least 2.5 times the calibrated gain setting.	
Frequency Response			
Direct Coupled	dc to 1 MHz	With Gain control set to 20K or less, selectable HF -3dB points are within +1dB, -2dB. Upper -3dB point reduces to 500 KHz at 50K gain, and 250 KHz at 100K gain.	
AC Coupled	2 Hz or less, at lower frequency -3dB point to 1 MHz.		
Bandwidth			
HF -3dB	100 Hz to 1 MHz, 9 steps in 1-3 sequence.	Reference approximately 5.0 V peak-to-peak output at 1 KHz.	
Accuracy	+1dB, -2dB	Reading range -2dB to -5dB.	
LF -3dB	0.1 Hz to 10 kHz 6 steps in 1-10 sequence.		
Accuracy	+1dB, -2dB	Reading range -2dB to -5dB.	

Table 2-1 (Con't)

Characteristics	Performance Requirements	Supplemental Information
Output		
Voltage Swing	±5 V.	
Current		±20 mA
R _o		5 Ω or less
Minimum Load Impedance		250 Ω
Common Mode		
Normal Mode Range	±5 V.	
Rejection Range	100 dB, dc to 50 kHz.	Direct coupled with inputs shorted together.
÷100 Mode Range	±50 V.	
Rejection Range	50 dB, dc to 50 kHz.	Direct coupled with less than 100 V peak-peak sine-wave input
DC Offset Range		At least + and - 1 V.
Maximum Safe Input Voltage		
Direct Coupled Normal Mode		15 V (dc plus peak ac) to 5 MHz or less.
÷100 Mode		350 V (dc plus peak ac) to 5 MHz or less.
AC Coupled		350 V (dc plus peak ac) with coupling capacitor precharged.
Maximum Input Gate Current (each input)		50 pA at 25°C.
Noise (Tangentially Measured)	25 μV or less.	With 25 Ω or less source resistance, from 5 Hz to 1 MHz.
Voltage Drift with	,	
Short Term		10 μV (peak-to-peak) per hour
Long Term		20 μV (peak-to-peak) per hour.
Voltage Drift with Temperature		100 μV per °C.
Input R and C		Resistance, 1 MΩ. Capacitance, Approximately 47 p

Table 2-2
ENVIRONMENTAL

Characteristics	Information	
Temperature		
Operating	0°C to +50°C	
Storage	-40°C to +75°C	
Altitude		
Operating	To 15,000 feet; maximun operating temperature decreased by 1° C/1000 feet from 5000 to 15,000 feet.	
Storage	To 50,000 feet.	
Vibration		
Operating and Non-operating	With the instrument complete and operating, vibration frequency swept from 10 to 55 to 10 Hz at 1 minute per sweep. Vibrate 15 minutes in each of the three major axes at 0.015 inch total displacement. Hold 10 minutes at any major resonanace, or in none, at 55 Hz. Total time 75 minutes.	
Shock		
Operating and Non-operating	30 g's 1/2 sine, 11 ms duration 2 shocks in each direction along 3 major axes, for a total of 18 shocks.	

Table 2-3
PHYSICAL CHARACTERISTICS

Characteristics	Information
Overall Dimensions (measured at maximum points)	
Height	5.0 inches 12.7 centimeter
Width	2.5 inches 6.35 centimeter
Length	11.8 inches 30.0 centimeter
Net Weight (Instrument Only)	1.8 lbs. .82 kilograms

PERFORMANCE CHECK

Introduction

This procedure checks the electrical characteristics of the AM 502 that appear in the Specification section of this manual. If the instrument fails to meet the requirements given in this performance check, the calibration procedure should be performed. This procedure can also be used by an incoming inspection facility to determine acceptability of performance.

The electrical characteristics in Section 2 are valid only if the AM 502 is calibrated at an ambient temperature of $+20^{\circ}$ C to $+30^{\circ}$ C and operated at an ambient temperature of 0° C to $+50^{\circ}$ C. Forced air circulation is required for ambient temperature above $+40^{\circ}$ C.

Tolerances that are specified in this performance check procedure apply to the instrument under test and do not include test equipment error.

Test Equipment Required

The following test equipment, or equivalent, is required to perform the performance check. Test equipment characteristics listed are the minimum required to verify the performance of the equipment under test. Substitute equipment must meet or exceed the stated requirements. All test equipment is assumed to be operating within tolerances.

Special test devices are used where necessary to facilitate the procedure. Most of these are available from Tektronix, Inc. and can be ordered through your local Tektronix Field Office or representative.

Table 2-4
LIST OF TEST EQUIPMENT REQUIREMENTS

Description	Performance Requirements	Applications	Example
Oscilloscope	Bandwidth, dc to 2 MHz; Minimum deflection factor, 20 mV/div; sweep rate, .2 ms/div to 1 µs/div; accuracy, within 3%.	Used throughout procedure to provide display.	TEKTRONIX 5110. 5A13N, 5B10N.
Power Module	Three compartments or more.	All tests	TEKTRONIX TM 503 or TM 504.
Calibration Gen- erator	Amplitude calibration, 5 mV to 10 V; accuracy, ±0.25% into 1 MΩ; output, square wave at approximately 1 kHz.	Amplifier gain check	TEKTRONIX PG 506 Calibration Generator.*
Function Generator	Waveforms, sine and square; voltage amplitude, 10 V p-p; frequency range, 1 kHz to 50 kHz; accuracy, within 3%.	Signal generation for cross neutral- ization and high frequency cmrr check.	TEKTRONIX FG 501 Function Generator.*
Termination	Impedance, 50 Ω; accuracy, within 2%; connectors, bnc.	Output termination for signal gen-erator.	Tektronix Part No. 011-0049-01
Attenuator, 10X (4 required)	Impedance, 50 Ω; accuracy, within 2%; connectors, bnc.	Output attenuation for signal gen-erator. Noise check.	Tektronix Part No. 011-0059-02
Coaxial cable (2 required)	Impedance, 50 Ω; length, 36 inches; connectors, bnc.	Provides signal interconnection.	Tektronix Part No. 012-0057-01

^{*} Requires TM 500-Series power module

Preliminary Procedure

- 1. Ensure that all test equipment and the AM 502 under test are suitably adapted to the line voltage to be applied. Refer to the installation section of the power module manual.
- 2. Ensure that all test equipment is suitably adapted to the applied line voltage.
- 3. Install the AM 502 into the power module, and if applicable, install the TM 500 series test equipment into the test equipment power module.
- 4. Connect the equipment under test and the test equipment to a suitable line voltage source. Turn all equipment on and allow at least 20 minutes for the equipment to stabilize.

Initial Control Settings

Set the following controls during warm-up time:

AM 502

STEP GAIN DC BAL midrange ± 100 pushbutton out DC OFFSET **FINE and COARSE** midrange 100 **GAIN** fully clockwise (cal) CAL HF -3dB 1 MHz DC OFFSET LF -3dB pushbutton out +AC+ GND pushbutton out pushbutton out - AC - GND pushbutton in

Oscilloscope

Intensity, Focus

Set for well-defined trace and normal brightness.

Vertical Amplifier

Volts/Div + Input

Variable fully clockwise (cal)

1 V

dc

Time Base Plug-In

Time/Div .5 ms Variable (cal in)

Triggering + Slope

Mode p-p auto

Coupling ac
Source internal

Position Set so trace starts at left

side of graticule.

selected

Display Mode main sweep

Magnifier X1

Calibration Generator

Amplitude .5 V

PERFORMANCE CHECK PROCEDURE

NOTE

The oscilloscope vertical amplifier system gain, the 50-ohm termination, and 50-ohm attenuator are required to be calibrated within 0.5% accuracy prior to proceeding with Steps 1 and 2. The PG 506 calibration generator may be used to set the system to 0.5% accuracy.

- 1. Check Amplifier Gain. Gain Accuracy is within 2%
- a. Connect the calibration generator to the vertical amplifier input, using a 50-ohm cable.
- b. Set the calibration generator for a five-division, square-wave display; set the Deflection Error control for zero per cent error. Disconnect the 50-ohm cable from the vertical amplifier input.
- c. Connect the calibration generator to the AM 502 + input connector through a 50-ohm termination, and a 50-ohm 10X attenuator, using a 50-ohm cable.
- d. Connect a 50-ohm cable from the AM 502 OUTPUT connector to the vertical amplifier input connector.
- e. Adjust the AM 502 DC OFFSET controls (FINE and COARSE) until the OVER RANGE indicator light goes out. Position the crt display to the center of the graticule with the vertical amplifier position control.

Specification and Performance Check—AM 502

f. Check — using the AM 502 GAIN and calibration generator amplitude settings given in Table 2-5, check the vertical deflection within the given limits.

NOTE

The OVER RANGE indicator light must remain off during all switch settings. Adjust the DC OFFSET controls as required during the check procedure.

NOTE

Install both plug-in side covers and insert plugin into power module for checking gain below 5K setting to minimize noise on the display.

Table 2-5
AMPLIFIER GAIN ACCURACY

AM 502 GAIN Setting	AM 502 HF -3 dB Switch Setting	Calibration Generator Amplitude Setting	Vertical Deflection in Divisions
100	1 MHz	1 V	4.9 - 5.1
200	1 MHz	.5 V	4.9 - 5.1
500	1 MHz	.2 V	4.9 - 5.1
1K	1 MHz	.1 V	4.9 - 5.1
2K	1 MHz	50 mV	4.9 - 5.1
5K	1 MHz	20 mV	4.9 - 5.1
10K	10 kHz	10 mV	4.9 - 5.1
20K	10 kHz	5 mV	4.9 - 5.1
50K	10 kHz	2 mV	4.9 - 5.1
100K	10 kHz	1 mV	4.9 - 5.1

2. Check ÷100 Amplifier Gain Accuracy. Gain Accuracy is within 2%

- a. Disconnect all cables and repeat Step 1 parts a,b,c, and d.
- b. Set the AM 502 GAIN switch to 10K and the HF -3 dB switch to 1 MHz; push in the \div 100 pushbutton.
- c. Set the calibration generator amplitude control to 0.5 V. Adjust the AM 502 DC OFFSET control (FINE and COARSE) until the OVER RANGE indicator light goes out.
- d. Position the crt display to the center of the graticule with the vertical amplifier position control.
- e. Check—the crt display amplitude for 4.9 to 5.1 divisions.

- f. Disconnect the 10X attenuator from the AM 502 \pm input connector and connect it to the \pm input connector. Push in the \pm input GND pushbutton and push and release the \pm input GND pushbutton.
 - g. Check-repeat parts d and e of this step.
- h. Disconnect the 50-ohm termination, 10X attenuator, and cable.

3. Check Common Mode Rejection. Rejection Ratio is 100 dB, dc to 50 kHz

- a. Connect a dual-input connector cable between the AM 502 + input and input connector.
- b. Set the AM 502 \div 100 pushbutton out; the + input and input AC and GND pushbuttons out, and the GAIN switch to 100.
- c. Set the AM 502 HF -3 dB switch to .3 MHz and the LF -3 dB switch to DC.
- d. Set the vertical amplifier deflection factor for 10 mV/div.
- e. Connect a 50-ohm cable from the function generator to the center connector of the dual-input connector. Set the function generator for a 10 V, 50 kHz sine-wave output signal.
- f. Check the crt display for one division of vertical deflection or less.

4. Check ÷100 Common Mode Rejection. Rejection Ratio is 50 dB, dc to 50 kHz

- a. Disconnect the 50-ohm cable from the function generator and connect it to the calibration generator output. Press in the $\div 100$ pushbutton.
- b. Set the calibration generator for a 100 volt, square-wave output signal.
- c. Check—crt display for 3.1 divisions or less of vertical deflection.
- d. Disconnect the dual-input connector and 50-ohm cable from the AM 502 input connectors.

5. Check HF -3 dB Bandwidth. Accuracy is +1 dB, -2 dB

- a. Set the AM 502 HF -3 dB switch to 1 MHz, and the input pushbutton to GND.
- b. Set the vertical amplifier deflection factor to 1 V/div, and the time-base sweep rate to 1 ms/div.
- c. Set the function generator controls for a 50 mV, 1 kHz sine-wave output signal. (Use appropriate attenuation to eliminate input overdrive condition.)
- d. Connect a 50-ohm cable from the function generator to the AM 502 + input connector. Adjust the vertical deflection amplitude for a five-division display.
- e. Set the function generator output frequency to 1 MHz.
- f. Check—the amplitude of the crt display for 3.15 to 3.85 divisions.

NOTE

The specification in part f of this step must be met before proceeding with part g of this step.

g. Check—the remaining settings of the HF -3 dB switch, using Table 2-6 as reference. (Change time-base sweep rate as needed for lower frequencies.)

Table 2-6
HF -3 dB BANDWIDTH ACCURACY

AM 502 HF -3 dB Switch Setting	Function Generator Output Frequency	Vertical Deflection in Divisions
.3 MHz	300 kHz	3.15 - 3.85
.1 MHz	100 kHz	3.15 - 3.85
30 kHz	30 kHz	3.15 - 3.85
10 kHz	10 kHz	3.15 - 3.85
3 kHz	3 kHz	3.15 - 3.85
1 kHz	1 kHz	3.15 - 3.85
.3 kHz	300 Hz	3.15 - 3.85
.1 kHz	100 Hz	3.15 - 3.85

i. Set the HF -3 dB switch to 1 MHz.

6. Check LF -3 dB Bandwidth. Accuracy is +1 dB, -2 dB

a. Check—the settings of the LF -3 dB switch, using Table 2-7 as reference. Adjust the time-base sweep rate to obtain an appropriate display.

Table 2-7

LF -3 dB BANDWIDTH ACCURACY

AM 502 HF -3 dB Switch Setting	Function Generator Output Frequency	Vertical Deflection in Divisions	
10 kHz	10 kHz	3.15 - 3.85	
1 kHz	1 kHz	3.15 - 3.85	
.1 kHz	.1 kHz	3.15 - 3.85	
10 Hz	10 Hz	3.15 - 3.85	
1 Hz	1 Hz	3.15 - 3.85	

NOTE

The components used in the .1 Hz position are also used in the other positions of the switch; therefore, the tolerance of the .1 Hz position is checked.

- b. Press in the AM 502 + input AC pushbutton. Set the function generator controls for a 2 kHz output signal.
- c. Check—the amplitude of the crt display for 3.15 to 3.85 divisions.
- d. Disconnect the 50-ohm cable from the AM 502 \pm input connector and connect it to the \pm input connector.
- e. Press in the AM 502 + input GND pushbutton and input AC pushbutton. Press to release the input GND pushbutton.
- f. Check—the amplitude of the crt display for 3.15 to 3.85 divisions.
- g. Set the LF $-3\,\mathrm{dB}$ switch to the DC OFFSET position. Press in the AM 502 input GND pushbutton.
- h. Disconnect the 50-ohm cable from the AM 502 input connector.

7. Check Overall Noise (Tangentially Measured)

- a. Set the AM 502 GAIN control to 100K; press in to release the \pm input AC and GND pushbutton, and the \pm input AC pushbutton.
- b. Set the vertical amplifier deflection factor to 5 V/div. Set the time-base sweep rate to 10 μ s/div, and the trigger source switch to external.
- c. Connect a 50-ohm termination to the AM 502 \pm input connector; connect four series-connected 10X attenuators to the 50-ohm termination.
- d. Connect a 50-ohm cable from the calibration generator fast-rise output connector to the end of the attenuator string. Set the pulse duration control for 1 ms.
- e. Adjust the AM 502 DC OFFSET controls (FINE and COARSE) until the OVER RANGE indicator light goes out.
- f. Adjust the calibration generator pulse amplitude control and observe two noise bands as shown in Fig. 2-1A (remove one attenuator if necessary, to produce the desired display).
- g. Decrease the calibration generator pulse amplitude until the noise bands just merge. See Fig. 2-1B.
- h. Remove three of the attenuators and connect the signal through the 50-ohm attenuator (including the 50-ohm termination), to the vertical amplifier input and measure the pulse amplitude. Calculate the tangentially measured display noise as follows:

Noise (in
$$\mu$$
V) = $\frac{\text{Signal level (measured in part h)}}{\text{Attenuation Removed}}$

Typical figures are:

$$\frac{12 \text{ mV}}{10^{-3}}$$
 = 12 x 10⁻⁶ = 12 μ V of noise

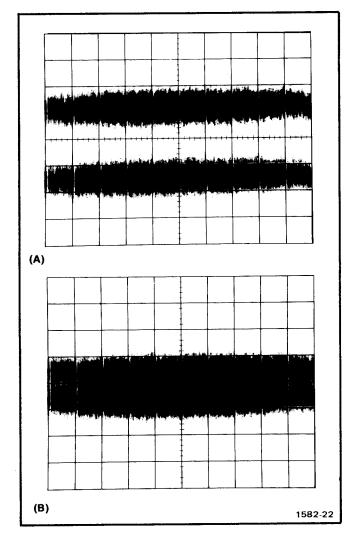


Fig. 2-1. Typical display of (A) two noise bands and (B) noise bands merged.

i. Disconnect all cables and equipment.

This completes the Performance Check of the AM 502 Differential Amplifier.

WARNING

THE FOLLOWING SERVICING INSTRUCTIONS ARE FOR USE BY QUALIFIED PERSONNEL ONLY. TO AVOID PERSONAL INJURY, DO NOT PERFORM ANY SERVICING OTHER THAN THAT CONTAINED IN OPERATING INSTRUCTIONS UNLESS YOU ARE QUALIFIED TO DO SO.

ADJUSTMENT

Introduction

This adjustment procedure is to be used to restore the AM 502 to the original performance specification. Adjustment need not be performed unless the instrument fails to meet the Performance Requirements of the Electrical Characteristics listed in the Specification section, or the Performance Check cannot be completed satisfactorily.

Completion of all adjustment steps in this procedure ensures that the instrument will meet the performance requirements listed in the Specification section. However, to fully ensure satisfactory performance, it is recommended that the Performance Check be performed after any adjustment is made.

Services Available

Tektronix, Inc. provides complete instrument repair and adjustment at local Field Service Centers and at the Factory Service Center. Contact your local Tektronix Field Office or representative for further information.

Test Equipment Required

The test equipment listed in Table 3-1, or equivalent, is required for adjustment of the AM 502. Specifications given for the test equipment are the minimum necessary for accurate adjustment and measurement. All test equipment is assumed to be correctly calibrated and operating within specification.

If other test equipment is substituted, control settings or calibration setup may need to be altered to meet the requirements of the equipment used.

A flexible plug-in extender, Tektronix Part No. 067-0645-03, is useful for troubleshooting or adjusting the AM 502; however, the complete Adjustment Procedure can be performed without use of the extender. Remove the power module cabinet to make adjustments to the AM 502 inside the power module.

Table 3-1
LIST OF TEST EQUIPMENT REQUIREMENTS

Description	Performance Requirements	Applications	Example
Oscilloscope	Bandwidth, dc to 2 MHz; minimum deflection factor, 20 mV/div; sweep rate, .2 ms/div to 1 μ s/div; accuracy, within 3%	Used throughout pro- cedure to provide display	TEKTRONIX 5110. 5A13N, 5B10N
Digital Voltmeter	Range, 0 to 50 V; accuracy, within 0.1%	Voltage measurements. Output voltage flatness check.	TEKTRONIX DM 501 Digita Multimeter.*
Power Module	Three compartments or more.	All tests	TEKTRONIX TM 503 or TM 504.
Calibration Gen- erator	Amplitude calibration, 5 mV to 10 V; accuracy, ±0.25% into 1 MΩ; output, square wave at approximately 1 kHz.	Amplifier gain check	TEKTRONIX PG 506 Calibration Generator.*
Function Generator	Waveforms, sine and square; voltage amplitude, 10 V p-p; frequency range, 1 kHz to 50 kHz; accuracy, within 3%.	Signal generation for cross neutral-ization and high frequency cmrr check.	TEKTRONIX FG 501 Function Generator.*
Autotransformer with ac voltmeter	Capable of supplying an output voltage from 90 to 132 V, ac; 120 watts of power at the upper limit.	Power supply check.	General Radio W10MTR3W Variac Autotransformer.
Input RC normalizer	Time constant, 1 MΩ x 47 pF; connectors, bnc.	Input capacitance check.	Tektronix Part No. 067-0541-00
Termination	Impedance, 50 Ω; connectors, bnc.	Output termination for signal generator.	Tektronix Part No. 011-0049-01.
Attenuator, 10X (4 required)	Impedance, 50 Ω ; connectors, bnc.	Output attenuation for signal gen-erator. Noise check.	Tektronix Part No. 011-0059-02
Coaxial cable (2 required)	Impedance, 50 Ω; length, 36 inches; connectors, bnc.	Provides signal interconnection.	Tektronix Part No. 012-0057-01.

^{*}Requires TM 500 Series power module

Preparation

- a. Remove the left and right side covers of the AM 502 to gain access to the component side of the circuit boards. Pull the rear end of the side cover outward from the side of the instrument (the cover snaps into place).
- b. Install the AM 502 into the left power module compartment, or if appropriate, connect the AM 502 to the power module by means of the flexible plug-in extender.
- c. Set the power module for the line voltage to be applied (see power module manual) and connect it to the variable autotransformer; connect the autotransformer to the line voltage source. Be sure that the power switch is off.
- d. Install the TM 500-series equipment, including the AM 502 into the power module.
- e. Connect all test equipment to a suitable line voltage source.
- f. Turn on all test equipment and allow at least 20 minutes for the equipment to warm up and stabilize.

Initial Control Settings

Set the following controls during warm-up time:

AM 502

STEP GAIN DC BAL	midrange
÷100	pushbutton out
DC OFFSET	
(FINE and COARSE)	midrange
GAIN	100
CAL	fully clockwise (cal)
HF -3dB	.1 kHz
LF -3dB	10 kHz
+ AC	pushbutton out
+ GND	pushbutton in
- AC	pushbutton out
– GND	pushbutton in

Oscilloscope

Intensity, Focus	Set for well-defined trace and normal brightness.

Vertical Amplifier

Volts/Div	1 V
+ Input	dc
Variable	fully clockwise (cal)

Time Base Plug-In

Time/Div	.5 ms
Variable	(cal in)
Triggering	
+ Slope	selected
Mode	peak-to-peak Auto
Coupling	ac
Source	internal
Position	Set so trace starts
	at left side of graticule.
Display Mode	main sweep
Magnifier	X1

ADJUSTMENT PROCEDURE

1. Check +15 Volt Power Supply

- a. Connect the digital voltmeter between the $\pm 15~\rm V$ test point on the Main circuit board, and chassis ground. See Fig. 3-1 for voltage test point location.
- b. Check—for a meter reading of +14.25 to +15.75 volts.
 - c. Disconnect the digital voltmeter.

2. Check -15 Volt Power Supply

- a. Connect the digital voltmeter between the $-15~\rm V$ test point on the Main circuit board, and chassis ground. See Fig. 3-1 for voltage test point location.
- b. Check—for a meter reading of -14.25 to -15.75 volts.
 - c. Disconnect the digital voltmeter.

3. Check +5 Volt Power Supply

- a. Connect the digital voltmeter between the ± 5 V test point on the Main circuit board, and chassis ground. See Fig. 3-1 for voltage test point location.
 - b. Check—for a meter reading of +4.50 to +5.50 volts.
- c. Adjust the autotransformer output voltage from the low limit to the high limit as indicated in Table 3-2. Meter reading should not vary more than $\pm 500\, \text{millivolts}$. Repeat this check for the ± 15 volt and $\pm 15\, \text{volt}$ supplies, except the meter reading should not vary more than $\pm 750\, \text{millivolts}$. Return the autotransformer to the nominal line voltage setting.
 - d. Disconnect the digital voltmeter.

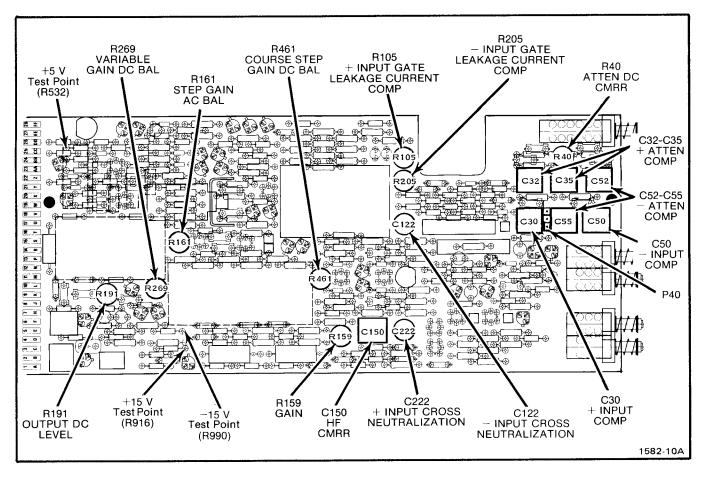


Fig. 3-1. Location of power supply test points, and all adjustments.

Table 3-2
POWER MODULE UNIVERSAL TRANSFORMER

Line Selector	Regulating Ranges		
Block Position	110-Volts Nominal	220-Volts Nominal	
L	90 V ac to 110 V ac	180 V ac to 220 V ac	
М	99 V ac to 121 V ac	198 V ac to 242 V ac	
Н	108 V ac to 132 V ac	216 V ac to 264 V ac	
Line Fuse			
Data	1.6 A slow-blow	0.8 A slow-blow	

4. Adjust Step Gain AC Balance

- a. Connect a 50-ohm cable from the AM 502 OUTPUT connector the vertical amplifier input.
- b. Check—for a trace shift of not more than three divisions vertically while switching the AM 502 GAIN switch bwtween the 100 and 100K positions.

- c. Adjust—Step Gain AC Bal, R161, for minimum trace shift while switching the AM 502 GAIN switch between the 100 and 100K positions. See Fig. 3-1 for adjustment location.
 - d. Return the AM 502 GAIN switch to the 100 position.

5. Adjust Variable Gain DC Balance

- a. Set the vertical amplifier deflection factor to .1 V/div.
- b. Check—for a trace shift of not more than one division vertically while rotating the AM 502 CAL control from fully clockwise to fully counterclockwise position.
- c. Adjust—Variable Gain DC BAL, R269, for minimum trace shift while rotating the AM 502 CAL control from fully clockwise to fully counterclockwise position. See Fig. 3-1 for adjustment location. Return control to the fully clockwise position and note the position of the trace.

6. Adjust Output DC Level

- a. Set the vertical amplifier deflection factor to 10 $\,$ mV/div.
- b. Without changing the AM 502 control positions from the previous step, note the position of the trace.
- c. Adjust—Output DC Level, R191, to place trace at same position as was established in Step 5 part c.

7. Adjust Coarse Step Gain DC Balance

- a. Set the vertical amplifier deflection factor to $2\,\text{V/div}$. Position the trace to the center of the graticule.
- b. Set the LF -3 dB switch to the DC position (not DC OFFSET), the GAIN switch to 10K, and STEP GAIN DC BAL control to midrange. Leave the other AM 502 controls as in the previous step.
- c. Adjust—Coarse Step Gain DC Bal, R461, to return trace to the center of the graticule. See Fig. 3-1 for adjustment location.
- d. Turn the AM 502 STEP GAIN DC BAL control fully clockwise and fully counterclockwise, and observe trace shift of at least two divisions above and below the graticule centerline, respectfully. Return the trace to the graticule center.

8. Adjust + Input Gate Leakage Current Compensation

- a. Connect a 50-ohm termination to the AM 502 \pm input connector.
- b. Set the vertical amplifier deflection factor to .1 V/div. Press in to release the \pm GND pushbutton.
- c. Adjust—the + Input Gate Leakage Current Comp, R105, while alternately depressing and releasing the +AC pushbutton, for minimum trace shift.
- d. Press in to release the +GND and +AC pushbuttons. Remove the termination from the + input connector.

9. Adjust — Input Gate Leakage Current Compensation

a. Connect a 50-ohm termination to the AM 502 - input connector.

- b. Press in to release the GND pushbutton.
- c. Adjust—the Input Gate Leakage Current Comp, R205, while alternately depressing and releasing the AC pushbutton, for minimum trace shift.
- d. Press in the GND pushbutton, and press in to release the -AC pushbutton. Remove the termination from the input connector.

NOTE

The oscilloscope vertical amplifier system gain, the 50-ohm termination, and 50-ohm attenuator are required to be calibrated within 0.5% accuracy prior to proceeding with Steps 10 and 11. The PG 506 calibration generator may be used to set the system to 0.5% accuracy.

10. Check Amplifier Gain.

- a. Set the vertical amplifier deflection factor to 1 V/div.
- b. Set the AM 502 GAIN control to 100, the LF -3dB switch to DC OFFSET, and the HF -3dB switch to 1 MHz.
- c. Connect the calibration generator to the vertical amplifier input, using a 50-ohm cable.
- d. Set the calibration generator for a five-division, square-wave display; set the Deflection Error control for zero per cent error. Disconnect the 50-ohm cable from the vertical amplifier input.
- e. Connect the calibration generator to the AM 502 + input connector through a 50-ohm termination, and a 50-ohm 10X attenuator, using a 50-ohm cable.
- f. Connect a 50-ohm cable from the AM 502 OUTPUT connector to the vertical amplifier input connector.
- g. Adjust the AM 502 DC OFFSET controls (FINE and COARSE) until the OVER RANGE indicator light goes out. Position the crt display to the center of the graticule with the vertical amplifier position control.
- h. Adjust-Gain, R159, for a five-division display amplitude.

Adjustment—AM 502

- i. Rotate the AM 502 CAL control fully counterclockwise and check for two divisions or less of display amplitude.
- j. Set the AM 502 CAL control fully clockwise and in the detent position.
- k. Check—using the AM 502 GAIN and calibration generator amplitude settings given in Table 3-3, check the vertical deflection within the given limits.

NOTE

The OVER RANGE indicator light must remain off during all switch settings. Adjust the DC OFFSET controls as required during the check procedure.

NOTE

Install both plug-in side covers and insert plugin into power module for checking gain below 5K setting to minimize noise on the display.

Table 3-3

AMPLIFIER GAIN ACCURACY

AM 502 GAIN Setting	HF -3 dB Switch Setting	Calibration Generator Amplitude Setting	Vertical Deflection in Divisions
100	1 MHz	1 V	4.9 - 5.1
200	1 MHz	.5 V	4.9 - 5.1
500	1 MHz	.2 V	4.9 - 5.1
1K	1 MHz	.1 V	4.9 - 5.1
2K	1 MHz	50 mV	4.9 - 5.1
5K	1 MHz	20 mV	4.9 - 5.1
10K	10 kHz	10 mV	4.9 - 5.1
20K	10 kHz	5 mV	4.9 - 5.1
50K	10 kHz	2 mV	4.9 - 5.1
100K	10 kHz	1 mV	4.9 - 5.1

11. Check \div 100 Amplifier Gain Accuracy. Gain Accuracy is within 2%

- a. Disconnect all cables and repeat Step 10 parts a through ${\bf f}$.
- b. Set the AM 502 GAIN switch to 10K and the HF -3dB switch to 1 MHz; push in the \div 100 pushbutton.
- c. Set the calibration generator amplitude control to 0.5 V. Adjust the AM 502 DC OFFSET control (FINE and COARSE) until the OVER RANGE indicator light goes out.

- d. Position the crt display to the center of the graticule with the vertical amplifier position control.
- e. Check—the crt display amplitude for 4.9 to 5.1 divisions.
- f. Disconnect the 10X attenuator from the AM 502 \pm input connector and connect it to the \pm input connector. Push in the \pm input GND pushbutton and push in to release the \pm input GND pushbutton.
 - g. Check-repeat parts d and e of this step.
- h. Disconnect the 50-ohm termination, 10X attenuator, and cable.

12. Adjust Input Cross Neutralization

- a. Set the AM 502 \div 100 pushbutton out, the + and input GND pushbuttons out, the + and AC pushbuttons in, the LF -3dB switch to DC, and the GAIN switch to 100.
- b. Connect a 50-ohm termination and a 10X attenuator to the AM 502 + input connector.
- c. Connect a 50-ohm cable from the fast-rise output of the calibration generator to the 10X attenuator on the + input of the AM 502.
- d. Set the calibration generator amplitude control to produce a five-division crt display (1 ms period). Set the time-base triggering controls for a stable, triggered display.
- e. Check—the crt display for roll-off or overshoot (upper corner of the leading edge) within 0.2 division.
- f. Adjust— + Input Cross Neutralization, C222, for minimum roll-off or overshoot (upper front corner of the leading edge).
- g. Disconnect the 50-ohm termination (with 10X attenuator and 50-ohm cable attached) from the + input connector, and connect it to the input connector.
 - h. Repeat part e of this step.
- i. Adjust— —Input Cross Neutralization, C122, for minimum roll-off or overshoot (upper corner of the leading edge).

3-6 REV. F MAR 1978

13. Adjust Input Capacitance

- a. Disconnect the 50-ohm termination from the input and connect a 47 pF Normalizer to the input connector; connect the 50-ohm termination to the other end of the 47 pF Normalizer.
- b. Increase the signal amplitude of the calibration generator to maintain five divisions of display.
- c. Check—the upper front corner of the leading edge for roll-off or overshoot, within 0.2 division.
- d. Adjust— Input Comp, C50, for the best upper corner on the leading edge of the displayed square wave.
- e. Disconnect the 47 pF Normalizer, with 50-ohm termination, 10X attenuator, and 50-ohm cable from the input connector and connect them to the + input connector.
 - f. Repeat part c of this step.
- g. Adjust— +Input Comp, C30, for best flat top (minimum roll-off or overshoot on upper front corner) on displayed waveform. See Fig. 3-1 for adjustment location.

14. Adjust Attenuator Compensation

- a. Press in the \div 100 pushbutton. Set the vertical amplifier deflection factor to .1 V/division.
- b. Disconnect the 10X attenuator that is connected between the 50-ohm termination and the 50-ohm cable; reconnect the 50-ohm cable to the termination.
- c. Readjust the calibration generator amplitude control to maintain a five-division display.
- d. Check—the crt display for roll-off or overshoot (upper corner of the leading edge) within 0.2 division.
- e. Adjust— + Atten Comp, C32-C35, for best flat top (minimum roll-off or overshoot on upper front corner) on displayed waveform. See Fig. 3-1 for adjustment location.
- f. Disconnect the 47 pF Normalizer, 50-ohm termination, and 50-ohm cable from the \pm input connector and connect them to the \pm input connector.
- g. Readjust the calibration generator amplitude control to maintain a five-division display.

- h. Check—the crt display for roll-off or overshoot (upper corner of the leading edge) within 0.2 division.
- i. Adjust— —Atten Comp, C52-C55, for best flat top (minimum roll-off or overshoot on upper front corner) on displayed waveform. See Fig. 3-1 for adjustment location.
- j. Disconnect the 47 pF Normalizer, 50-ohm termination and cable from the AM 502 input connector.

15. Adjust High Frequency Common Mode Rejection

- a. Connect a dual-input connector cable between the AM 502 + input and the input connector.
- b. Set the AM 502 \div 100 pushbutton out; the + input and input AC and GND pushbuttons out, and the GAIN switch to 100.
- c. Set the AM 502 HF -3 dB switch to .3 MHz and LF -3 dB switch to DC.
- d. Set the vertical amplifier deflection factor for 10 mV/div.
- e. Connect a 50-ohm cable from the function generator output to the center connector of the dual-input connector. Set the function generator for a 10 V, 50 kHz sine-wave output signal.
- f. Check—the crt display for one division of vertical deflection or less.
- g. Adjust—HF CMRR, C150, for minimum crt display amplitude within one division or less.

16. Adjust ÷ 100 Common Mode Rejection

- a. Set the AM 502 ÷ 100 pushbutton in.
- b. Set the vertical amplifier deflection factor for .1 V/div.
- c. Disconnect the 50-ohm cable from the function generator output and connect it to the calibration generator output. Set the calibration generator for a 100 volt output signal.
- d. Check—the crt display for one division of vertical deflection or less.

Adjustment—AM 502

- e. Adjust— Atten Comp, C52, and Atten DC CMRR, R40, for minimum crt display amplitude within one division or less. (Interaction between C52 and R40 will occur. Adjust for optimum display.)
- f. Disconnect the dual-input connector and 50-ohm cable from AM 502.

17. Check HF -3dB Bandwidth

- a. Set the AM 502 HF -3dB switch to 1 MHz, the \div 100 pushbutton out, the input pushbutton to GND.
- b. Set the vertical amplifier deflection factor to 1 V/div, and the time-base sweep rate to 1 ms/div.
- c. Set the function generator controls for a 50 mV, 1 kHz sine-wave output signal. (Use appropriate attenuation to eliminate input overdrive condition.)
- d. Connect a 50-ohm cable from the function generator output to the AM 502 + input connector. Adjust the vertical deflection amplitude for a five-division display.
- e. Set the function generator output frequency to 1 $\,\mathrm{MHz}.$
- f. Check—the amplitude of the crt display for 3.15 to 3.85 divisions.

NOTE

The specification in part f of this step must be met before proceeding with part g of this step.

g. Check—the remaining settings of the HF -3dB switch, using Table 3-4 as reference. (Change time-base sweep rate as needed for lower frequencies.)

Table 3-4

HF -3 dB BANDWIDTH ACCURACY

AM 502 HF -3 dB Switch Setting	Function Generator Output Frequency	Vertical Deflection in Divisions
.3 MHz	300 kHz	3.15 - 3.85
.1 MHz	100 kHz	3.15 - 3.85
30 kHz	30 kHz	3.15 - 3.85
10 kHz	10 kHz	3.15 - 3.85
3 kHz	3 kHz	3.15 - 3.85
1 kHz	1 kHz	3.15 - 3.85
.3 kHz	300 Hz	3.15 - 3.85
.1 kHz `	100 Hz	3.15 - 3.85

h. Set the HF -3 dB switch to 1 MHz.

18. Check LF -3 dB Bandwidth

a. Check—the LF -3 dB switch settings, using Table 3-5 as reference. (Adjust the time-base sweep rate as needed for lower frequency settings.)

Table 3-5

LF -3 dB BANDWIDTH ACCURACY

AM 502 LF -3 dB Switch Setting	Function Generator Output Frequency	Vertical Deflection in Divisions
10 kHz	10 kHz	3.15 - 3.85
1 kHz	1 kHz	3.15 - 3.85
.1 kHz	.1 kHz	3.15 - 3.85
10 Hz	10 Hz	3.15 - 3.85
1 Hz	1 Hz	3.15 - 3.85

NOTE

The components used in the .1 Hz position are also used in the other positions of the switch; therefore, the tolerance of the .1 Hz position is checked.

- b. Press in the AM 502 + input AC pushbutton. Set the function generator controls for a 2 kHz output signal.
- c. Check—the amplitude of the crt display for 3.15 to 3.85 divisions.
- d. Disconnect the 50-ohm cable from the AM 502 \pm input connector and connect it to the input connector.
- e. Press in the AM 502 \pm input GND pushbutton and \pm input AC pushbutton. Press to release the \pm input GND pushbutton.
- f. Check—the amplitude of the crt display for 3.15 to 3.85 divisions.
- g. Set the LF $-3\,\mathrm{dB}$ switch to the DC OFFSET position. Press in the AM 502 input GND pushbutton.
- h. Disconnect the 50-ohm cable from the AM 502 input connector.

19. Check Overall Noise (Tangentially Measured)

- a. Set the AM 502 GAIN control to 100K; press in to release the \pm input AC and GND pushbutton, and the \pm input AC pushbutton.
- b. Set the vertical amplifier deflection factor to 5 V/division. Set the time-base sweep rate to $10\,\mu\text{s}/\text{division}$ and the trigger source switch to external.
- c. Connect a 50-ohm termination to the AM 502 \pm input connector; connect four series-connected 10X attenuators to the 50-ohm termination.
- d. Connect a 50-ohm cable from the calibration generator fast-rise output connector to the end of the attenuator string. Set the pulse duration control for 1 ms.
- e. Adjust the AM 502 DC OFFSET controls (FINE and COARSE) until the OVER RANGE indicator light goes out.
- f. Adjust the calibration generator pulse amplitude control and observe two noise bands as shown in Fig. 3-2A (remove one attenuator if necessary, to produce the desired display).
- g. Decrease the calibration generator pulse amplitude until the noise bands just merge. See Fig. 3-2B.
- h. Remove three of the attenuators and connect the signal through the 50-ohm attenuator (including the 50-ohm termination), to the vertical amplifier input and measure the pulse amplitude. Calculate the tangentially measured display noise as follows:

Noise (in
$$\mu$$
V) =
$$\frac{\text{Signal level (measured in part h)}}{\text{Attenuation Removed}}$$

Typical figures are:

$$\frac{12 \text{ mV}}{10^{-3}}$$
 = 12 x 10⁻⁶ = 12 μ V of noise

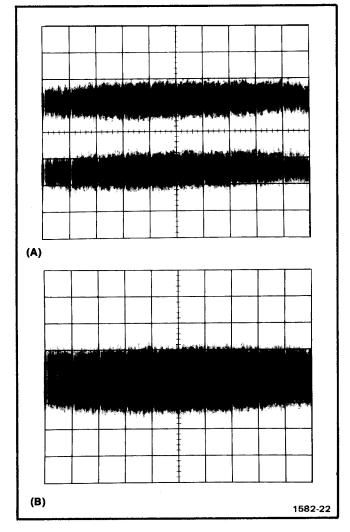


Fig. 3-2. Typical display of (A) two noise bands and (B) noise bands merged.

i. Disconnect all cables and equipment.

This completes the Adjustment procedure of the AM 502 Differential Amplifier.

MAINTENANCE AND INTERFACING INFORMATION

Preventive Maintenance

There are no special preventive maintenance procedures that apply to the AM 502. Refer to the power module instruction manual for general preventive maintenance procedures and instructions.

Corrective Maintenance

Refer to the power module instruction manual for general corrective maintenance procedures and instructions.

Troubleshooting

Use the Performance Check, Adjustment Procedure, and Circuit Description as aids to locate trouble in the event of equipment failure. The test equipment listed in the Performance Check and Adjustment Procedure will prove useful in troubleshooting the AM 502.

Functions Available at Rear Connector

A slot between pins 23 and 24 on the rear connector identifies the AM 502 as a member of the signal source family. Insert a barrier in the corresponding position of the power module jack to prevent other than signal source plug-ins from being used in that compartment; this protects the plug-in should specialized connections be made to that compartment. Consult the Building A System section of the power module manual for further information.

Signal outputs, or other specialized connections, are made to the rear interface connectors as shown in Fig. 4-1. The Signal In and Signal Ground are not factory wired. If Signal In connections are made, use coaxial cable. Connect one end to the pads as shown in Fig. 4-1 and the other end to the front panel bnc connector.

REPACKAGING FOR SHIPMENT

If the Tektronix instrument is to be shipped to a Tektronix Service Center for service or repair, attach a tag showing: owner (with address) and the name of an individual at your firm that can be contacted. Include complete instrument serial number and a description of the service required.

Save and re-use the package in which your instrument was shipped. If the original packaging is unfit for use or not available, repackage the instrument as follows:

Surround the instrument with polyethylene sheeting to protect the finish of the instrument. Obtain a carton of corrugated cardboard of the correct carton strength and having inside dimensions of no less than six inches more than the instrument dimensions. Cushion the instrument by tightly packing three inches of dunnage or urethane foam between carton and instrument, on all sides. Seal carton with shipping tape or industrial stapler.

The carton test strength for your instrument is 200 pounds.

4-1

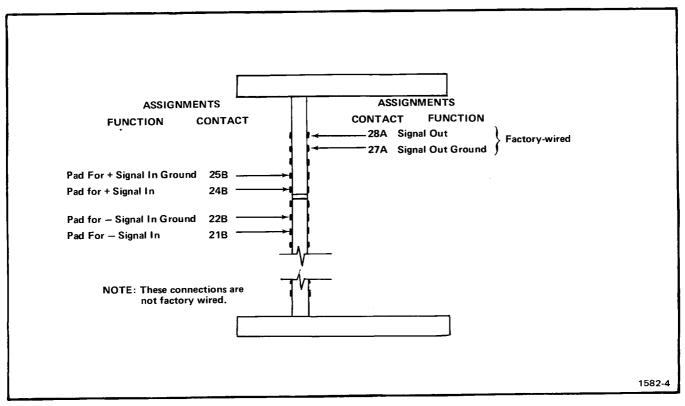


Fig. 4-1. Input/Output assignments at rear connector.

CIRCUIT DESCRIPTION

Introduction

This section of the manual contains a description of the circuitry used in the AM 502 Differential Amplifier. Individual descriptions are separated into the following parts; Input Coupling, Input Attenuators, Input Overdrive Protection, Gate Current Compensation, Preamp, Floating Power Supply, Common Mode Rejection, Cross Neutralization, DC Offset, LF —3 dB Frequency Selector, Gain Switching Amplifier, HF —3 dB Frequency Selector, Variable Gain Stage, Output Amplifier, Overrange Indicator, and Power Supplies. The circuit titles correspond to those listed in the Block Diagram. The numbered diamond by each title refers to the corresponding circuit diagram in the Diagrams section of this manual.

Input Coupling



A signal applied to the + or - input connector passes through the input-coupling selector switch to the input-attenuator circuit. The signals can be ac coupled, dc coupled or disconnected internally. The + and - inputs are identical except for circuit numbers. Except where needed for clarification, only the amplifier for the + input will be described in detail.

When the input-coupling switch is in the DC position, the input signal is coupled directly to the attenuator. In the AC position, the signal is coupled through capacitor C10. The capacitor prevents the dc component of the signal from passing to the attenuator. The GND position disconnects the input signal from the attenuators and connects it to ground through C10 and R10. This represents the same load to the signal source as was presented when the input switch was set in the AC mode. Thus, a ground reference to the input of the amplifier is provided without the need to remove the applied signal source. The PRECHG position limits the instantaneous current, caused by connecting a large signal to the input connector, by inserting R10 in series with C10.

Input Attenuators



The input attenuators for the + and - inputs are identical and are conventional RC type attenuators. The attenuators are frequency-compensated voltage dividers that provide a straight-through or NORM position, or a $\div 100$ position. Resistor R30 provides the 1 megohm input resistance, while C28, C30, and stray capacitance make up the 47 picofarad input capacitance in the NORM position. In the $\div 100$ position, R35 in series with R38 comprise the 1 megohm input resistance, and also form

the 100X divider for the dc component of the signal. The parallel combination of C32 and C34 in series with C38, comprise the ac 100X circuit. Each attenuator contains an adjustable capacitor to provide correct attenuation at high frequencies, and adjustable shunt capacitance to provide correct input capacitance.

Input Overdrive Protection



Fuses F100 and F200, combined with diode clamps CR108, CR109, CR208, and CR209 provide overdrive protection for Q121A and B at the ± 16.2 volt level. The ± 16.2 volt level is set by zeners VR438 and VR448, through isolation diodes CR438 and CR448, and the clamp diodes. If the signal amplitude exceeds approximately 16.2 volts, depending on the duration of the excess voltage, the fuses will open.

Gate Current Compensation



The leakage associated with the gates of the input Field-Effect Transistors (FET) may be as high as 100 picoamperes. This amount of leakage current (through 1 megohm, R30 or R50, to ground) will produce an offset of 100 microvolts, which at high gain settings, drives the output into overrange. To compensate for this effect, the gates of the input FET may be adjusted to zero volts by returning R107 and R106 through variable controls R105 and R205 to a slightly negative supply voltage.

Leakage current associated with the gates of the input FET and the overdrive protection diodes, increases rapidly with temperature, approximately doubling for every 10 degree Celsius. To compensate for this increase, a temperature sensitive input-current balancing network is included, using thermistors RT104 and RT203 as sensing elements.

As the gate current increases due to a temperature increase, the resistance of the thermistors decreases, thus maintaining a constant voltage difference as a function of temperature.

Preamp



From the input attenuators, the signal is coupled to the preamp. The preamp consists of two identical feedback amplifiers connected in a differential configuration. Except where needed for clarification, only the amplifier for the + input will be described in detail.

Circuit Description—AM 502

The supply voltages for the two amplifiers are obtained from a common power supply that is bootstrapped to the input to improve the common-mode rejection ratio of the preamp.

The feedback circuit for the + input consists of comparator Q121A, error amplifier Q131A and Q136, feedback modifier R125 and R126, and output load resistor R136. R126 is the gain-setting resistor for the amplifier. Diode CR131 connected between the base-emitter of Q131A protects the transistor against reverse-bias breakdown, and also provides negative feedback to stabilize the circuit at higher frequencies.

The last stage of the preamp consists of Q152A and Q152B. Zener diode VR252 provides the collector supply voltage for this stage. The gain of this stage is approximately 5.

STEP GAIN DC BAL control, R158, balances the dc level at the output and is adjusted for minimum dc level shift when the GAIN switch is rotated throughout its range.

Floating Power Supply



Refer to Fig. 5-1. The supply voltages for the preamp and the dc offset circuit are derived from a seriesconnected resistor chain consisting of R431, R433, R435, two current sources, plus a current-setting transistor, Q443. The positive current source consists of Q410A and B, and Q416. The negative current source consists of Q422A and B, and Q426.

Any common-mode changes that occur in the input amplifier, except at the output, are coupled to the power supply through Q451, which is connected as an emitter follower. It acts as a X1 bootstrap amplifier whose gain is maintained very close to unity by the minimum loading presented to its output by the high collector impedance of Q121A and B, thus achieving good bootstrap efficiency.

Quiescently, approximately 13 milliamperes of current flows through Q422A and B. R421 shunts about 4 milliamperes of this current from Q422A, leaving a balance of 9 milliamperes flowing through the positive current source, Q410A, Q410B, and Q416. Approximately 12 milliamperes of current from the negative current source also flows from Q422B and Q426, and to the dc offset circuit, Q141A and B, and Q145A and B. In any mode except the dc offset, the 12 milliamperes of current splits equally and becomes the negative current source for Q121A and B. About 4 milliamperes is shunted through Q136 and Q236. This leaves the balance of 9 milliamperes of current to flow through Q121A and B, and to the summing point at the collector of Q416.

Resistors R440, R441, R442, and R443 set the base voltage of Q443, with CR441 providing the temperature compensation. Current through R445 divides between Q443 and Q445. The current through Q445 is the negative current source for Q131A and B. After flowing through Q131A and B, it returns and passes through emitter-follower Q431, then to the summing point at the collector of Q416.

The 1.9 milliamperes of current flowing from Q443 and through resistor chain R435, R433, and R431, continues to the summing point at the collector of Q416. The voltage drop across each resistor, in reference to the emitter of Q451, sets the "floating" power supplies. The difference in the 9 milliamperes of current needed to satisfy the current demand of the positive current source, and what flows into the summing point from Q433, Q431, and Q121A and B, will come from Q451.

Common Mode Rejection



One of the primary functions of the preamp is to reject any common-mode component in the input signals, and amplify only the difference. In the extreme case of the inputs tied together and a common voltage applied, the output of the amplifier is ideally zero, and would actually be zero provided that the characteristics of all corresponding elements on the two sides of the amplifier (see Figs. 1-5 and 1-6) were exactly matched. In practice, any mismatch will cause a differential output. Even with perfect matching, there is still a common-mode output current resulting in an undesirable common-mode signal applied to the subsequent stages of the amplifier (common mode gain).

The floating power supply eliminates these difficulties and improves the common mode rejection. It can be seen that the entire power supply and amplifier moves an amount equal to the common-mode voltage, and that no change in voltage or current levels occur anywhere within the amplifier as a result of this common-mode voltage, except for Q136 and Q236 drain to gate voltage. Thus the only mismatch of importance is that of Q136, Q236 amplification factors. Being in the third stage of the amplifier at a point of relatively large differential signal level, this causes only a small degradation in common mode rejection.

At higher frequencies the stray capacitances from various points in the preamp to ground begin to inject significant current into the amplifier as a result of common mode signals. Differential capacitor C150, connected from a point in the floating power supply to the two output lines, is used to inject adjustable currents into the output to compensate for the change.

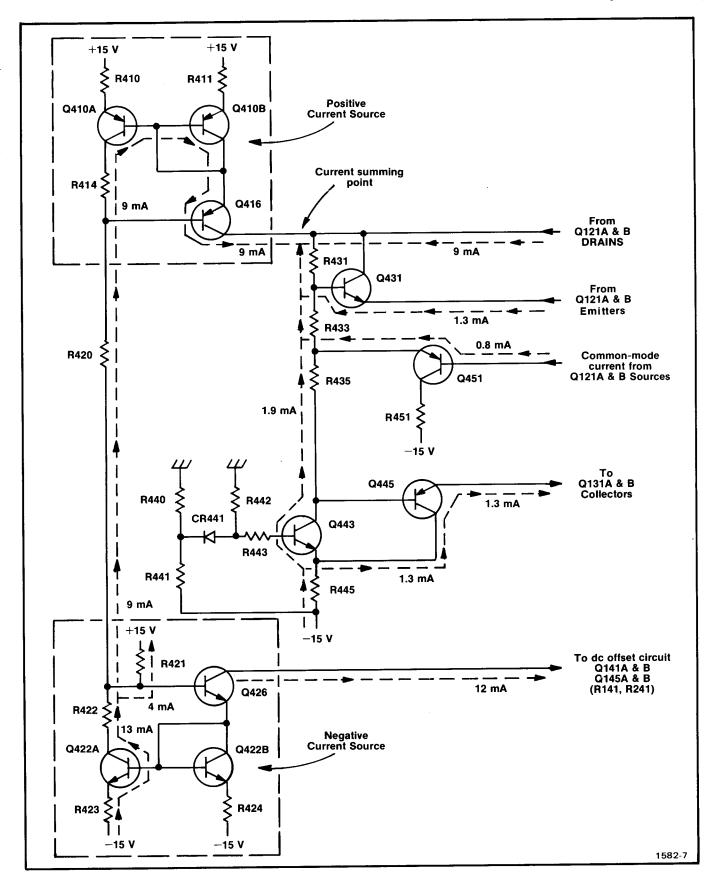


Fig. 5-1. Preamp floating power supply showing positive and negative current source.

Cross Neutralization



Refer to Fig. 5-2. The use of a common bootstrap power supply results in an undesirable capacitive coupling between the two inputs. Consider the effect of applying +1 volt to the + input of Q121A while keeping the - input of Q121B at zero volts.

The results are: an output current of i_o of 4 milliamperes, and a shift of all supply voltages and several other voltage levels by ± 0.5 volt due to the divider action of R126 and R226, operating into the bootstrapped power supply system. Specifically, the drain of Q121B also rises ± 0.5 volt and injects current i_a through the drain to gate capacitance, C_{dg} of Q121B, and into the \pm input. If there is any impedance between the \pm input and ground, i_a will develop a voltage across the impedance, and since it's applied to the \pm input, subtracts from the original \pm input, thereby causing an erroneous output (see Fig. 5-3).

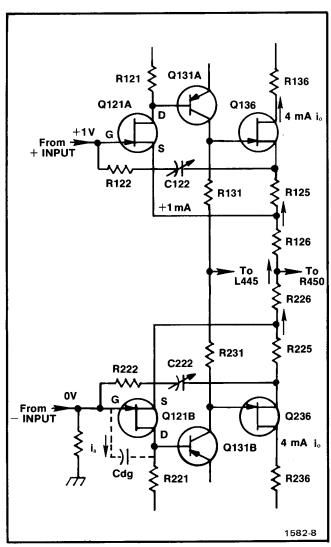


Fig. 5-2. Simplified preamp circuit showing cross neutralization circuit.

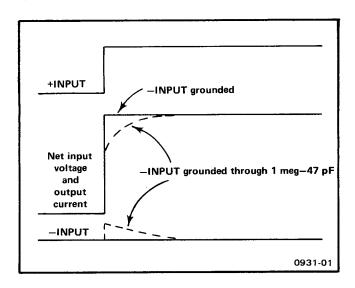


Fig. 5-3. Effects of high impedance to ground in the - input.

Capacitor C122, connected in the + input circuit, is adjusted to divert the undesirable capacitive current away from the input point, minimizing the unbalance. Capacitor C222 performs the same function for signals applied to the - input circuit.

DC Offset 1

The purpose of the dc offset system is to allow a differential slide-back measurement; that is, to buck out small dc components of input signals and allow the amplifier to amplify only the varying components, while keeping the differential capability. This means that when a dc voltage is applied across the inputs, some means must be found to balance out the resultant output current.

The dc offset system, with zero input signal conditions and resultant currents is illustrated in Fig. 5-4. Source and output currents remain at 2 milliamperes and 5 milliamperes respectively, thereby producing a balanced output. The amplifier and dc offset system illustrated in Fig. 5-5 shows the current conditions when a 0.25 volt dc signal voltage is applied to the + input. This voltage causes an additional 1 milliampere of current to flow through R126 and R226. However, if a 1 milliampere offset current is supplied by the offset generator, the net output current will remain at its zero signal value.

When the dc offset is not used, the dc offset controls are switched by the LF -3 dB switch, to a fixed reference voltage divider R460, R461, and R462.


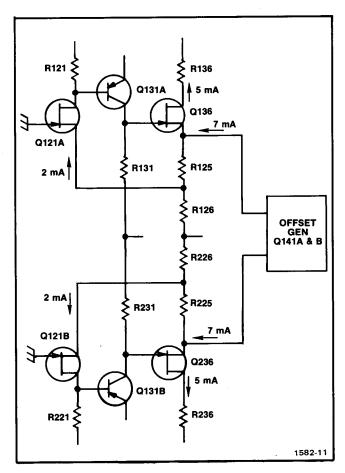


Fig. 5-4. Amplifier and dc offset system, showing balanced output conditions.

LF −3 dB Frequency Selector <

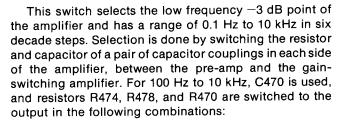

-3 dB frequency	100 Hz	1 kHz	10 kHz
output resistors	R470	R470 & R474	R470 & R478

Whenever R474 or R478 is not used on the output side of C258, the resistance (R474 or R478) is placed across the input to keep the high frequency load resistance, as seen by the pre-amp, constant.

Capacitor C472 is switched across C470 for the lower three ranges (0.1 Hz, 1 Hz, and 10 Hz) and C258 is shorted out for dc coupling. Resistors R474, R478, and R470 are switched in the same combination for the lower three positions as they are for the three higher positions.

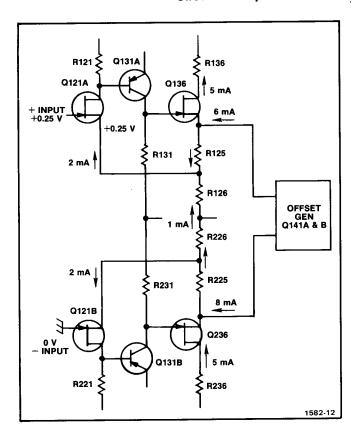


Fig. 5-5. Amplifier and dc offset system, showing a 0.25 volt dc offset and resultant current conditions.

Gain Switching Amplifier



The gain switching amplifier is a balanced differential configuration very similar to the preamp but with a fixed power supply. A gain of 0.48 to 480 is changed by switching different values of resistance with switch S480. Q162A, Q164A, and Q168 comprise the switching amplifier. Q167 and Q267 serve two functions; the base-emitter junctions serve as reverse bias protection for Q168 and Q268 and they drive the over range indicator circuitry to indicate a differential over range condition.

Step Gain AC Balance control, R161, is used to remove any initial unbalance in the gate-to-source voltages of Q162A and B, and sets the voltage across the gain-setting resistors when the differential input (gate-to-gate) voltage is zero.

Capacitors C164 and C264 are in the circuit for peaking at frequencies near 2 MHz. Inductors L164 and L264 reduce the gain at high frequencies, thus preventing oscillations.

HF -3 dB Frequency Selector



This switch selects the high frequency —3 dB point of the amplifier and has a range of 100 Hz through 1 MHz in nine steps of a 1-3 sequence. This is accomplished by simply switching different values of capacitance across the output of the gain-switched stage.

Variable Gain Stage



The variable gain stage (Q171A, Q171B and Q181, Q281) is a feedback amplifier having a gain range of approximately 2.5:1, determined by CAL control R178. Feedback resistors R173 and R273, plus gain setting resistors R177, R277, and R178 set the overall gain of the stage at 10.

The common mode signal at the output of the feedback amplifier is compared to ground by comparator Q195 and Q196. A dc feedback loop is provided to feed back the output of the comparator to Q171A and B via emitter follower Q269 to maintain the junction of R175 and R275 at a zero volt level.

Variable Gain DC Balance R269 adjusts the input levels of Q171A and B so their source voltages will be equal. With the source levels equal, rotating the front panel VAR control causes no shift of the dc output level.

Output DC Level control R191 adjusts for a zero-volt dc level at the output of the amplifier. R191 and R193 act as a current source or sink to set the voltage difference across R175 and R275. Their common junction is at zero volts due to feedback action described earlier. Q185 and Q285 act as the negative current source for the variable gain stage.

Diodes CR181 and CR182 act as peak detectors, providing a signal to the over-range detector circuit (Q501 emitter). A voltage swing greater than +6 volts on the collector of Q181 or Q281 will produce an over range indication, lighting indicator lamp DS50.

Output Amplifier



The output amplifier stage consists primarily of transistors Q561A and B and Q571A and B and provides a ± 5 volt differential voltage swing with a zero-volt dc level. The output R_o is 5 ohms or less. Transistors Q555 and Q558 comprise the positive and negative current sources respectively.

The output stage is short-circuit protected to ground. If the current in R571 causes the voltage at Q571B collector to fall below the level of Q555 emitter, CR555 conducts, taking current from Q555, thereby reducing the drive to Q571B, limiting the output current.

Over Range Indicator



When the LF -3 dB switch is in a position other than DC, there is no on-screen indication of the dc conditions in the preamp or output amplifier, and it may be driven into non-linearity or overload by a dc component, leading to erroneous displays. The over range indicator circuit detects this condition and so indicates by means of indicator lamp DS50.

The over range circuit consists of three input transistors Q501, Q511, and Q521, monostable multivibrator Q532 and Q541, indicator lamp driver Q540. When any of the three input transistors become turned on by an overrange signal from the respective amplifiers, a turn-on current is supplied to the base of Q532. This starts the multivibrator action which has a minimum on time of approximately 200 milliseconds; thus the indicator lamp is held on for a minimum of 200 milliseconds regardless of the duration of the over-range signal. The timing components are C536, R536, and R531.

Power Supplies



Regulation of the +15 volt and -15 volt supply is accomplished through the respective operational amplifier, U910 and U990. Zener diode VR910 sets the reference voltage for U910; the +15 volt supply sets the reference for U990. The operational amplifiers control the current to the AM 502 through series-pass transistors located in the power module. The inputs of both integrated circuits seek common levels by varying the current supplied to the unit. Should the power supply load increase, pin 2 of U990 goes more positive, causing pin 6 to go more negative. This action increases current flow in the series-pass transistor and restores the voltage to the initial level. The +15 volt supply operates in the same manner except the polarities are different.

Transistors Q910 and Q990 act as current sensing regulators. Should the current increase to a level that the voltage drop across R970 (R920) is great enough to cause Q990 (Q910) to conduct, current is taken from the seriespass transistor until it is turned off. This action limits the current available to the amplifier, should an overload or short circuit occur.

The ± 5 volt supply is comprised of Q940 and Q942. The base voltage level of Q942 is set by divider R945 and R946. The base-emitter drop of the transistors produce the ± 5 volt supply.

OPTIONS

No options are available at this time.

REPLACEABLE ELECTRICAL PARTS

PARTS ORDERING INFORMATION

Replacement parts are available from or through your local Tektronix, Inc. Field Office or representative.

Changes to Tektronix instruments are sometimes made to accommodate improved components as they become available, and to give you the benefit of the latest circuit improvements developed in our engineering department. It is therefore important, when ordering parts, to include the following information in your order: Part number, instrument type or number, serial number, and modification number if applicable.

If a part you have ordered has been replaced with a new or improved part, your local Tektronix, Inc. Field Office or representative will contact you concerning any change in part number.

Change information, if any, is located at the rear of this manual.

SPECIAL NOTES AND SYMBOLS

X000 Part first added at this serial number

00X Part removed after this serial number

ITEM NAME

In the Parts List, an Item Name is separated from the description by a colon (:). Because of space limitations, an Item Name may sometimes appear as incomplete. For further Item Name identification, the U.S. Federal Cataloging Handbook H6-1 can be utilized where possible.

ABBREVIATIONS

ACTR	ACTUATOR	PLSTC	PLASTIC
ASSY	ASSEMBLY	QTZ	QUARTZ
CAP	CAPACITOR	RECP	RECEPTACLE
CER	CERAMIC	RES	RESISTOR
CKT	CIRCUIT	RF	RADIO FREQUENCY
COMP	COMPOSITION	SEL	SELECTED
CONN	CONNECTOR	SEMICOND	SEMICONDUCTOR
ELCTLT	ELECTROLYTIC	SENS	SENSITIVE
ELEC	ELECTRICAL	VAR	VARIABLE
INCAND	INCANDESCENT	ww	WIREWOUND
LED	LIGHT EMITTING DIODE	XFMR	TRANSFORMER
NONWIR	NON WIREWOUND	XTAL	CRYSTAL

CROSS INDEX—MFR. CODE NUMBER TO MANUFACTURER

O0853 SANGAMO ELECTRIC CO., S. CAROLINA DIV. O BOX 128 PICKENS, SC 29671 O1002 GEMERAL ELECTRIC COMPANY, INDUSTRIAL AND POWER CAPACITOR PRODUCTS DEPARTMENT 1201 2ND STREET HUDSON FALLS, NY 12839 MILWAUKEE, WI 53204 O1121 ALLEN-BRADLEY COMPANY O1295 TEXAS INSTRUMENTS, INC., SEMICONDUCTOR GROUP O BOX 5012, 13500 N CENTRAL EXPRESSWAY O5091 TRI-ORDINATE CORPORATION O5091 TRI-ORDINATE CORPORATION O7100 TELEDYNE SEMICONDUCTOR O1100	Mfr. Code	Manufacturer	Address	City, State, Zip
ALLEN-BRADLEY COMPANY 1201 2ND STREET HUDSON FALLS, NY 12839			P O BOX 128	PICKENS, SC 29671
Olipho		AND POWER CAPACITOR PRODUCTS DEPARTMENT		
O BOX 5012, 13500 N CENTRAL EXPRESSWAY DALLAS, TX 75222 DALLAS TX SALES DALLAS TX TS 222			1201 2ND STREET SOUTH	MILWAUKEE, WI 33204
MOTOROLA, INC., SEMICONDUCTOR PROD. DIV. 5005 E MCDOWELL RD, PO BOX 20923 PHOENIX, AZ 85036 BERKELEY HEIGHTS, NJ 07922	01295		D O DOY EGIS 13500 N CENTRAL	
MOTOROLA, INC., SEMICONDUCTOR PROD. DIV. 5005 E MCDOWELL RD, PO BOX 20923 PHOENIX, AZ 85036 BERKELEY HEIGHTS, NJ 07922 PARCHILD SEMICONDUCTOR, A DIV. OF FAIRCHILD CAMERA AND INSTRUMENT CORP. 464 ELLIS STREET MOUNTAIN VIEW, CA 94042 HAWTHORNE, CA 90250 HAWTHORNE, CA 90261 HAWTHORNE, CA 90250 HAWTHORNE, CA 90		GROUP		DAT.T.AS TX 75222
TRI-ORDINATE CORPORATION 343 SNYDER AVENUE BERKELEY HEIGHTS, NJ 0/922		CONTROL DOOR DEED DEED DEED		
TRI-UNIAL COMPUTOR, A DIV. OF FAIRCHILD SEMICONDUCTOR, A DIV. OF FAIRCHILD CAMERA AND INSTRUMENT CORF. 464 ELLIS STREET MOUNTAIN VIEW, CA 94042 12515 CHADRON AVE. HAWTHORNE, CA 90250			3/3 CNVDED AVENUE	BERKELEY HEIGHTS, NJ 07922
FAIRCHILD CAMERA AND INSTRUMENT CORP. 464 ELLIS STREET MOUNTAIN VIEW, CA 94042 07910 TELEDYNE SEMICONDUCTOR 12515 CHADRON AVE. HAWTHORNE, CA 90250 08806 GENERAL ELECTRIC CO., MINIATURE LAMP PRODUCTS DEPARTMENT NELA PARK CLEVELAND, OH 44112 11237 CTS KEENE, INC. 3230 RIVERSIDE AVE. PASO ROBLES, CA 93446 15454 RODAN INDUSTRIES, INC. 2905 BLUE STAR ST. ANAHEIM, CA 92806 156289 SPRAGUE ELECTRIC CO. 2905 BLUE STAR ST. ANAHEIM, CA 92806 156280 SPRAGUE ELECTRIC CO. 1007 MCGRAW-EDISON CO. 2536 W. UNIVERSITY ST. ST. LOUIS, MO 63107 17440 CHICAGO MINIATURE LAMP WORKS 4431 RAVENSWOOD AVE. CHICAGO, IL 60640 17292 ERIE TECHNOLOGICAL PRODUCTS, INC. 644 W. 12TH ST. ERIE, PA 16512 173138 BECKMAN INSTRUMENTS, INC., HELIPOT DIV. 2500 HARBOR BLUD. FULLERTON, CA 92634 174970 JOHNSON, E. F., CO. 299 10TH AVE. S. W. WASECA, MN 56093 175042 TRW ELECTRONIC COMPONENTS, IRC FIXED 178488 STACKPOLE CARBON CO. 299 10TH AVE. S. W. WASECA, MN 56093 178488 STACKPOLE CARBON CO. 5135 MAGNOLIA AVE. PHILADELPHIA, PA 19108 178488 STACKPOLE CARBON CO. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 180099 TEXTRONIX, INC. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 1801087, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 1801087, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 1801087, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 1801087, INC. 3029 E WASHINGTON STREET 190201 MALLORY CAPACITOR CO., DIV. OF 191637 DALE ELECTRONICS, INC. 9. O. BOX 609 180201 MALLORY CAPACITOR CO., DIV. OF 191637 DALE ELECTRONICS, INC. 9. O. BOX 609 180203 DALE ELECTRONICS, INC. 9. O. BOX 609 180204 BOLLER AND CO., INC. 9. O. BOX 609 180205 DALE ELECTRONICS, INC. 9. O. BOX 609			343 SHIDER AVENUE	
THELEDYNE SEMICONDUCTOR 12515 CHADRON AVE. HAWTHORNE, CA 90250	07263		AGA FILIS STREET	MOUNTAIN VIEW, CA 94042
RELECTRIC CO., MINIATURE LAMP PRODUCTS DEPARTMENT NELA PARK CLEVELAND, OH 44112				
LAMP PRODUCTS DEPARTMENT NELA PARK CLEVELAND, OH 44112			12010 Cimbion May	•
11237 CTS KEENE, INC. 3230 RIVERSIDE AVE. PASO ROBLES, CA 93446 15454 RODAN INDUSTRIES, INC. 2905 BLUE STAR ST. ANAHEIM, CA 92806 56289 SPRAGUE ELECTRIC CO. NORTH ADAMS, MA 01247 71400 BUSSMAN MFG., DIVISION OF MCGRAW-EDISON CO. 2536 W. UNIVERSITY ST. ST. LOUIS, MO 63107 71744 CHICAGO MINIATURE LAMP WORKS 4433 RAVENSWOOD AVE. CHICAGO, IL 60640 72982 ERIE TECHNOLOGICAL PRODUCTS, INC. 644 W. 12TH ST. ERIE, PA 16512 73138 BECKMAN INSTRUMENTS, INC., HELIPOT DIV. 2500 HARBOR BLVD. FULLERTON, CA 92634 74970 JOHNSON, E. F., CO. 299 10TH AVE. S. W. WASECA, MN 56093 75042 TRW ELECTRONIC COMPONENTS, IRC FIXED REISTORS, PHILADELPHIA DIVISION 401 N. BROAD ST. PHILADELPHIA, PA 19108 78488 STACKPOLE CARBON CO. FOR BOX 500 BEAVERTON, OR 97077 80099 TEXTRONIX, INC. FOR BOX 500 BEAVERTON, OR 97077 80294 BOURNS, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. LOS ANGELES, CA 90069 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. P. O. BOX 609 COLUMBUS, NE 68601 91637 DALE ELECTRONICS, INC. P. O. BOX 609 COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBUS, NE 68601 TECHNOLOR OF THE COLUMBU	08806		NETA PARK	CLEVELAND, OH 44112
15454 RODAN INDUSTRIES, INC. 2905 BLUE STAR ST. ANAHEIM, CA 92806 56289 SPRAGUE ELECTRIC CO. 71400 BUSSMAN MFG., DIVISION OF MCGRAW- EDISON CO. 2536 W. UNIVERSITY ST. ST. LOUIS, MO 63107 71744 CHICAGO MINIATURE LAMP WORKS 4433 RAVENSWOOD AVE. CHICAGO, IL 60640 72982 ERIE TECHNOLOGICAL PRODUCTS, INC. 644 W. 12TH ST. ERIE, PA 16512 73138 BECKMAN INSTRUMENTS, INC., HELIPOT DIV. 2500 HARBOR BLVD. FULLERTON, CA 92634 74970 JOHNSON, E. F., CO. 299 10TH AVE. S. W. WASECA, MN 56093 75042 TRW ELECTRONIC COMPONENTS, IRC FIXED RESISTORS, PHILADELPHIA DIVISION 401 N. BROAD ST. PHILADELPHIA, PA 19108 878488 STACKPOLE CARBON CO. BEAVERTON, OR 97077 80294 BOURNS, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. LOS ANGELES, CA 90069 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 90 BOX 372 INDIANAPOLIS, IN 46206 91637 DALE ELECTRONICS, INC. P. O. BOX 609 COLUMBUS, NE 68601 91637 DALE ELECTRONICS, INC. P. O. BOX 609 COLUMBUS, NE 68601	11007			
SPRAGUE ELECTRIC CO. SPRAGUE ELECTRIC CO. BUSSMAN MFG., DIVISION OF MCGRAW- EDISON CO. 2536 W. UNIVERSITY ST. ST. LOUIS, MO 63107		•		ANAHEIM, CA 92806
### PROPRIES BUSSMAN MFG., DIVISION OF MCGRAW-EDISON CO. 2536 W. UNIVERSITY ST. ST. LOUIS, MO 63107		· · · · · · · · · · · · · · · · · · ·	,	NORTH ADAMS, MA 01247
EDISON CO. 2536 W. UNIVERSITY ST. CHICAGO MINIATURE LAMP WORKS 4433 RAVENSWOOD AVE. CHICAGO, IL 60640 ERIE, PA 16512 FULLERTON, CA 92634 FULLERTON, C				
71744 CHICAGO MINIATURE LAMP WORKS 72982 ERIE TECHNOLOGICAL PRODUCTS, INC. 73138 BECKMAN INSTRUMENTS, INC., HELIPOT DIV. 74970 JOHNSON, E. F., CO. 75042 TRW ELECTRONIC COMPONENTS, IRC FIXED RESISTORS, PHILADELPHIA DIVISION 78488 STACKPOLE CARBON CO. 80009 TEKTRONIX, INC. 80294 BOURNS, INC., INSTRUMENT DIV. 80294 BOURNS, INC., INSTRUMENT DIV. 81483 INTERNATIONAL RECTIFIER CORP. 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 91637 DALE ELECTRONICS, INC. 14433 RAVENSWOOD AVE. 644 W. 12TH ST. 6716 644 W. 12TH ST. 6716 644 W. 12TH ST. 6716	71400		2536 W. UNIVERSITY ST.	ST. LOUIS, MO 63107
72982 ERIE TECHNOLOGICAL PRODUCTS, INC. 644 W. 12TH ST. ERIE, PA 16512 73138 BECKMAN INSTRUMENTS, INC., HELIPOT DIV. 2500 HARBOR BLVD. 299 10TH AVE. S. W. WASECA, MN 56093 74970 JOHNSON, E. F., CO. 299 10TH AVE. S. W. WASECA, MN 56093 75042 TRW ELECTRONIC COMPONENTS, IRC FIXED RESISTORS, PHILADELPHIA DIVISION 401 N. BROAD ST. PHILADELPHIA, PA 19108 78488 STACKPOLE CARBON CO. ST. MARYS, PA 15857 80009 TEKTRONIX, INC. PO BOX 500 BEAVERTON, OR 97077 80294 BOURNS, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. LOS ANGELES, CA 90069 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 3029 E WASHINGTON STREET PO BOX 372 INDIANAPOLIS, IN 46206 91637 DALE ELECTRONICS, INC. 10. BOX 609 COLUMBUS, NE 68601	71744			CHICAGO, IL 60640
73138 BECKMAN INSTRUMENTS, INC., HELIPOT DIV. 2500 HARBOR BLVD. 299 10TH AVE. S. W. WASECA, MN 56093 74970 JOHNSON, E. F., CO. 299 10TH AVE. S. W. WASECA, MN 56093 75042 TRW ELECTRONIC COMPONENTS, IRC FIXED RESISTORS, PHILADELPHIA DIVISION 401 N. BROAD ST. PHILADELPHIA, PA 19108 78488 STACKPOLE CARBON CO. ST. MARYS, PA 15857 80009 TEXTRONIX, INC. PO BOX 500 BEAVERTON, OR 97077 80294 BOURNS, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. LOS ANGELES, CA 90069 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 3029 E WASHINGTON STREET PO BOX 372 INDIANAPOLIS, IN 46206 91637 DALE ELECTRONICS, INC. P. O. BOX 609 COLUMBUS, NE 68601	-			ERIE, PA 16512
74970 JOHNSON, E. F., CO. 75042 TRW ELECTRONIC COMPONENTS, IRC FIXED RESISTORS, PHILADELPHIA DIVISION 78488 STACKPOLE CARBON CO. 80009 TEKTRONIX, INC. 80294 BOURNS, INC., INSTRUMENT DIV. 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. 90201 MALLORY AND CO., INC. 90201 MALLORY AND CO., INC. 91637 DALE ELECTRONICS, INC. 299 10TH AVE. S. W. WASECA, MN 56093 PHILADELPHIA, PA 19108 ST. MARYS, PA 15857 BEAVERTON, OR 97077 6135 MAGNOLIA AVE. 9220 SUNSET BLVD. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 1NDIANAPOLIS, IN 46206 COLUMBUS, NE 68601 TICKANOF N. 1077				FULLERTON, CA 92634
TRW ELECTRONIC COMPONENTS, IRC FIXED RESISTORS, PHILADELPHIA DIVISION 78488 STACKPOLE CARBON CO. 80009 TEKTRONIX, INC. 80294 BOURNS, INC., INSTRUMENT DIV. 80294 BOURNS, INC., INSTRUMENT DIV. 81483 INTERNATIONAL RECTIFIER CORP. 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 90205 WASHINGTON STREET P O BOX 372 P. O. BOX 609 11 INDIANAPOLIS, IN 46206 COLUMBUS, NE 68601 TICKNANDE N. 1077				WASECA, MN 56093
RESISTORS, PHILADELPHIA DIVISION 78488 STACKPOLE CARBON CO. 80009 TEKTRONIX, INC. 80294 BOURNS, INC., INSTRUMENT DIV. 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. 9220 SUNSET BLVD. 9220 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 9080 372 91637 DALE ELECTRONICS, INC. 9090 INTERNATIONAL RECTIFIER CORP. 9100 BEAVERTON, OR 97077 8029 SUNSET BLVD. 9100 SANTA ANA, CA 92711 9100 BOX 372 P. O. BOX 609 9100 COLUMBUS, NE 68601 9100 STEVENDER NO 1007				
78488 STACKPOLE CARBON CO. 80009 TEKTRONIX, INC. 80294 BOURNS, INC., INSTRUMENT DIV. 81483 INTERNATIONAL RECTIFIER CORP. 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 91637 DALE ELECTRONICS, INC. ST. MARYS, PA 15857 BEAVERTON, OR 97077 RIVERSIDE, CA 92506 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 1NDIANAPOLIS, IN 46206 COLUMBUS, NE 68601 TICKANOF N. 10707	73042		401 N. BROAD ST.	PHILADELPHIA, PA 19108
80009 TEKTRONIX, INC. P O BOX 500 BEAVERTON, OR 97077 80294 BOURNS, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. LOS ANGELES, CA 90069 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 3029 E WASHINGTON STREET P O BOX 372 INDIANAPOLIS, IN 46206 91637 DALE ELECTRONICS, INC. P. O. BOX 609 COLUMBUS, NE 68601	78488			ST. MARYS, PA 15857
80294 BOURNS, INC., INSTRUMENT DIV. 6135 MAGNOLIA AVE. RIVERSIDE, CA 92506 81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. LOS ANGELES, CA 90069 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 3029 E WASHINGTON STREET P O BOX 372 INDIANAPOLIS, IN 46206 91637 DALE ELECTRONICS, INC. P. O. BOX 609 COLUMBUS, NE 68601			P O BOX 500	BEAVERTON, OR 97077
81483 INTERNATIONAL RECTIFIER CORP. 9220 SUNSET BLVD. LOS ANGELES, CA 90069 87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 3029 E WASHINGTON STREET P O BOX 372 INDIANAPOLIS, IN 46206 91637 DALE ELECTRONICS, INC. P. O. BOX 609 COLUMBUS, NE 68601			6135 MAGNOLIA AVE.	RIVERSIDE, CA 92506
87034 ILLUMINATED PRODUCTS INC., A SUB OF OAK INDUSTRIES, INC. 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 91637 DALE ELECTRONICS, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92711 1029 E WASHINGTON STREET P O BOX 372 INDIANAPOLIS, IN 46206 P. O. BOX 609 COLUMBUS, NE 68601			9220 SUNSET BLVD.	LOS ANGELES, CA 90069
OAK INDUSTRIES, INC. 90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 91637 DALE ELECTRONICS, INC. 2620 SUSAN ST, PO BOX 11930 SANTA ANA, CA 92/11 3029 E WASHINGTON STREET P O BOX 372 INDIANAPOLIS, IN 46206 P. O. BOX 609 COLUMBUS, NE 68601				
90201 MALLORY CAPACITOR CO., DIV. OF P. R. MALLORY AND CO., INC. 90201 P. R. MALLORY AND CO., INC. 90201 P. R. MALLORY AND CO., INC. 90202 WASHINGTON STREET P O BOX 372 INDIANAPOLIS, IN 46206 P. O. BOX 609 COLUMBUS, NE 68601 91637 DALE ELECTRONICS, INC. 91637 DALE ELECTRONICS, INC.	0,051		2620 SUSAN ST, PO BOX 11930	SANTA ANA, CA 92711
P. R. MALLORY AND CO., INC. 91637 DALE ELECTRONICS, INC. 3029 E WASHINGTON STREET P O BOX 372 P. O. BOX 609 COLUMBUS, NE 68601 THE PROPERTY 10707	90201			
91637 DALE ELECTRONICS, INC. P. O. BOX 609 COLUMBUS, NE 68601			3029 E WASHINGTON STREET	
9163/ DALLE ELECTRONICS, INC.		•	P O BOX 372	
MICHAUDE NV 10707	91637	DALE ELECTRONICS, INC.	P. O. BOX 609	
		· · · · · · · · · · · · · · · · · · ·	40 MARBLEDALE ROAD	TUCKAHOE, NY 10707

Ckt No.		Serial/Mod Eff	el No. Dscont	Name & Description	Mfr Code	Mfr Part Number
Al	670-2733-00	во10100	в051720	CKT BOARD ASSY:MAIN	80009	670-2733-00
Al.	670-2733-01	B051721		CKT BOARD ASSY:MAIN	80009	670-2733-01
A1	670-2733-02	в052350		CKT BOARD ASSY:MAIN	80009	670-2733-02
C10) C20)	295-0155-00			CAP.,SET,MTCHD:(2) 0.1UF,1%	80009	295-0155-00
C28	281-0663-00			CAP.,FXD,CER DI:10.4PF,1%,500V	72982	
C30	281-0081-00			CAP., VAR, AIR DI:1.8-13PF, 375VDC	74970	189-6-5
C32	281-0081-00			CAP., VAR, AIR DI:1.8-13PF, 375VDC	74970	189-6-5
C34	281-0645-00			CAP.,FXD,CER DI:8.2PF,+/-0.25PF,500V	• 72982	374-011C0H0829C
C35	281-0081-00			CAP., VAR, AIR DI:1.8-13PF, 375VDC	74970	
C38	283-0626-00			CAP., FXD, MICA D:1800PF, 5%, 500V		D195E182J0
C48	281-0663-00			CAP.,FXD,CER DI:10.4PF,1%,500V	72982	
C50	281-0081-00			CAP., VAR, AIR DI:1.8-13PF, 375VDC	74970	189-6-5
C52	281-0081-00			CAP., VAR, AIR DI:1.8-13PF, 375VDC	74970	189-6-5
C54	281-0645-00			CAP., FXD, CER DI:8.2PF, +/-0.25PF, 500V		374-011C0H0829C
C55	281-0081-00			CAP., VAR, AIR DI:1.8-13PF, 375VDC		189-6-5
C58	283-0626-00			CAP.,FXD,MICA D:1800PF,5%,500V	00853	D195E182J0
C107	283-0005-00			CAP.,FXD,CER DI:0.01UF,+100-0%,250V	72982	8131N300Z5U0103P
C122	281-0093-00	во10100	в039999	CAP., VAR, CER DI:5.5-18PF	72982	538-011A5.5-18
C122	281-0092-00	в040000		CAP., VAR, CER DI:9-35PF, 200V	72982	538-011 D9-35
C131	281-0612-00			CAP.,FXD,CER DI:5.6PF,+/-0.5PF,500V	72982	374-001C0H0569D
C146	283-0177-00			CAP.,FXD,CER DI:1UF,+80~20%,25V	72982	8131N039 E 105Z
C150	281-0114-00			CAP., VAR, AIR DI:1.3-5.4PF, 750V	74970	189-0352-075
C164	281-0523-00			CAP., FXD, CER DI:100PF, +/-20PF, 500V	72982	301-000U2M0101M
C176	281-0618-00			CAP., FXD, CER DI:4.7PF,+/-0.5PF,200V	72982	374-001C0H0479D
C195	281-0518-00			CAP.,FXD,CER DI:47PF,+/-9.4PF,500V	72982	301-000U2J0470M
C207	283-0005-00			CAP.,FXD,CER DI:0.01UF,+100-0%,250V	72982	8131N300Z5U0103P
C222	281-0093-00	в010100	в039999	CAP., VAR, CER DI:5.5-18PF	72982	538-011A5.5~18
C222	281-0092-00	в040000		CAP., VAR, CER DI:9-35PF, 200V	72982	538-011 D9-35
C231	281-0612-00			CAP.,FXD,CER DI:5.6PF,+/-0.5PF,500V	72982	374-001C0H0569D
C252	290-0512-00	XB052350		CAP.,FXD,ELCTLT:22UF,20%,15V	56289	196D226X0015KA1
C264	281-0523-00			CAP.,FXD,CER DI:100PF,+/-20PF,500V	72982	301-000U2M0101M
C276	281-0618-00			CAP.,FXD,CER DI:4.7PF,+/-0.5PF,200V	72982	374-001C0H0479D
C431	290-0527-00			CAP.,FXD,ELCTLT:15UF,20%,20V	90201	TDC156M020FL
C433	290-0527-00			CAP., FXD, ELCTLT: 15UF, 20%, 20V	90201	TDC156M020FL
C435	290-0527-00			CAP.,FXD,ELCTLT:15UF,20%,20V	90201	TDC156M020FL
C438	283-0177-00			CAP.,FXD,CER DI:1UF,+80-20%,25V	72982 90201	8131N039 E 105Z TDC156M020FL
C443	290-0527-00			CAP.,FXD,ELCTLT:15UF,20%,20V		
C448	283-0177-00			CAP., FXD, CER DI:1UF, +80-20%, 25V		8131N039 E 105Z
C450	281-0613-00	B010100	B051720	CAP.,FXD,CER DI:10PF,+/-1PF,200V	72982	374-001C0G0100F
C450	281-0629-00	B051721		CAP.,FXD,CER DI:33PF,5%,600V	72982	308-000C0G0330J
C470	283-0594-00			CAP., FXD, MICA D:0.001UF, 1%, 100V	00853	D151F102F0
C471	283-0594 - 00			CAP.,FXD,MICA D:0.001UF,1%,100V	00853	D151F102F0
C472	285-0809-00			CAP.,FXD,PLSTC:1UF,10%,50V	56289	LP66A1A105K
C473	285-0809-00			CAP., FXD, PLSTC:1UF, 10%, 50V	56289	LP66A1A105K
C490	281-0511-00			CAP., FXD, CER DI:22PF,+/-2.2PF,500V	72982	301-000C0G0220K 301-000P3K720J
C491	281-0721-00			CAP.,FXD,CER DI:72PF,5%,500V	72982 00853	D153F3040G0
C493	283-0604-00			CAP.,FXD,MICA D:304PF,2%,300V		
C494	283-0594-00			CAP., FXD, MICA D:0.001UF, 1%, 100V	00853	D151F102F0
C495	285-0627-00			CAP.,FXD,PLSTC:0.0033UF,5%,100V	56289	410P33251
C496	285-0598-00			CAP.,FXD,PLSTC:0.01UF,5%,100V	01002	61F10AC103
C497	285-0702-00			CAP.,FXD,PLSTC:0.033UF,5%,100V	56289	410P33351
C498	285-0703-00			CAP.,FXD,PLSTC:0.1UF,5%,100V	56289	410P112

Ckt No.	Tektronix Part No.	Serial/Model No. Eff Dscont	Name & Description	Mfr Code	Mfr Part Number
	005 0500 00		CAP., FXD, PLSTC: 0.22UF, 20%, 100V	56289	410P22491
C499	285-0633-00		· · · · · · · · · · · · · · · · · · ·	72982	8121-N088Z5U104M
C530	283-0111-00		CAP.,FXD,CER DI:0.1UF,20%,50V		
C536	290-0534-00		CAP., FXD, ELCTLT: 1UF, 20%, 35V	56289	196D105X0035HA1
C552	281-0511-00		CAP., FXD, CER DI:22PF, +/-2.2PF, 500V	72982	301-000C0G0220K
C560	283-0111-00		CAP.,FXD,CER DI:0.1UF,20%,50V	72982	8121-N088Z5U104M
C570	283-0111-00		CAP.,FXD,CER DI:0.1UF,20%,50V	72982	8121-N088Z5U104M
C910	290-0512-00		CAP.,FXD,ELCTLT:22UF,20%,15V	56289	196D226X0015KA1
	283-0111-00		CAP., FXD, CER DI:0.1UF, 20%, 50V	72982	8121-N088Z5U104M
C922			CAP., FXD, CER DI:0.1UF, 20%, 50V	72982	
C924	283-0111-00		CAP.,FAD,CER DI:0.10F,209,50V	72982	8121-N088Z5U104M
C926	283-0111-00		CAP.,FXD,CER DI:0.1UF,20%,50V	72302	Olai Moodaotom
C930	290-0527-00		CAP.,FXD,ELCTLT:15UF,20%,20V	90201	TDC156M020FL
C942	283-0111-00		CAP., FXD, CER DI:0.1UF, 20%, 50V	72982	8121-N088Z5U104M
			CAP., FXD, ELCTLT:15UF, 20%, 20V	90201	TDC156M020FL
C960	290-0527-00			72982	
C964	283-0111-00		CAP.,FXD,CER DI:0.1UF,20%,50V	72982	8121-N088Z5U104M
C966	283-0111-00		CAP.,FXD,CER DI:0.1UF,20%,50V	12902	0121-N000230104H
C968	283-0111-00		CAP.,FXD,CER DI:0.1UF,20%,50V	72982	8121-N088Z5U104M
gp3.00	152-0333-00		SEMICOND DEVICE:SILICON, 35V, 100MA	80009	152-0323-00
CR108	152-0323-00		SEMICOND DEVICE:SILICON, 35V, 100MA		152-0323-00
CR109	152-0323-00		SEMICOND DEVICE:SILICON, 35V, 100MA		1N4152
CR131	152-0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA		
CR150	152-0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA	07910	1N4152
CR181	152 - 0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA	07910	1N4152
CR208	152-0323-00		SEMICOND DEVICE:SILICON, 35V, 100MA	80009	152-0323-00
CR209	152-0323-00		SEMICOND DEVICE:SILICON, 35V, 100MA	80009	152-0323-00
	152-0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA	07910	1N4152
CR231			SEMICOND DEVICE:SILICON, 30V, 150MA	07910	1N4152
CR250 CR281	152-0141-02 152-0141-02		SEMICOND DEVICE:SILICON,30V,150MA	07910	1N4152
			THE SAME STATE OF THE SAME	07910	ln4152
CR438	152-0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA		
CR441	152-0141-02		SEMICOND DEVICE: SILICON, 30V, 150MA	07910	1N4152
CR448	152-0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA	07910	1N4152
CR531	152-0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA	07910	1N4152
CR555	152-0141-02		SEMICOND DEVICE: SILICON, 30V, 150MA	07910	lN4152
CDEEG	152-0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA	07910	1N4152
CR558	152-0141-02		SEMICOND DEVICE:SILICON,30V,150MA	07910	1N4152
CR560	152-0141-02			07910	1N4152
CR561	152-0141-02		SEMICOND DEVICE:SILICON,30V,150MA	07910	1N4152
CR570	152-0141-02		SEMICOND DEVICE:SILICON, 30V, 150MA		
CR571	152-0141-02		SEMICOND DEVICE:SILICON,30V,150MA	07910	1N4152
DOFO	150-0107-00	во10100 во52349	LAMP.INCAND:12V,0.04A,T-1	87034	227AS25
DS50			LAMP, INCAND: 5V, 0.06A, SEL	08806	683AS15
DS50	150-0048-01	B052350	TAMP, INCAMD: 5V, 0. OOA, SEE		227AS25
DS900	150-0107-00		LAMP, INCAND:12V, 0.04A, T-1		CM7220
DS900	150-0109-00	в052350	LAMP, INCAND: 18V, 26MA	/1/44	CH1/220
F100	150-0024-00	1	FUSE, CARTRIDGE: 3AG, 0.06A, 250V, FAST BLOW	71400	AGC 1/16
F100 F200	159-0024-00 159-0024-00		FUSE, CARTRIDGE: 3AG, 0.06A, 250V, FAST BLOW	71400	AGC 1/16
T1.0	131-0955-00	1	CONNECTOR, RCPT, : BNC, FEMALE, W/HARDWARE	05091	31-279
J10			POST, BDG, ELEC: ASSEMBLY	80009	129-0103-00
J15	129-0103-00		CONNECTOR, RCPT, :BNC, FEMALE, W/HARDWARE	05091	
J20 J550	131-0955-00 131-0818-00		CONNECTOR, RCPT, :BNC, FEMALE CONNECTOR, RCPT, :BNC, FEMALE	91836	
				80009	108-0409-00
L164	108-0409-00		COIL, RF:17.5UH	80009	
L264	108-0409-00)	COIL, RF:17.5UH		
L416	276-0507-00)	SHIELDING BEAD,: 0.6UH	78488	
L426	276-0507-00)	SHIELDING BEAD,:0.6UH	78488	
L445	276-0507-00		SHIELDING BEAD,: 0.6UH	78488	57-0180-7D 500B

	Tektronix	Serial/Mod	el No.		Mfr	
Ckt No.	Part No.	Eff	Dscont	Name & Description	Code	Mfr Part Number
Q121A,B	151-1027-00			TRANSISTOR:SILICON, JFE, N-CHAN	80009	151-1027-00
Q131A,B	151-0261-00			TRANSISTOR:SILICON, PNP, DUAL	80009	151-0261-00
0136	151-1025-00			TRANSISTOR:SILICON, JFE, N-CHANNEL	80009	151-1025-00
Q141A,B	151-0176-00			TRANSISTOR:SILICON, NPN, DUAL	04713	SD555
Q145A,B	151-0261-00			TRANSISTOR:SILICON, PNP, DUAL	80009	151-0261-00
Q152A,B	153-0586-00			TRANSISTOR:SILICON, PNP, PAIR	80009	153-0586-00
Q162A,B	151-1036-00			TRANSISTOR:SILICON, JFE, N-CHANNEL, DUAL	80009	151-1036-00
Q164A,B	153-0586-00			TRANSISTOR:SILICON, PNP, PAIR	80009	153-0586-00
Q167	151-0190-00			TRANSISTOR:SILICON, NPN	80009	151-0190-00
Q168	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q171A,B	151-1029-00			TRANSISTOR:SILICON, JFE, N-CHANNEL, DUAL	80009	151-1029-00
Q181	151-0188-00			TRANSISTOR:SILICON, PNP	01295	2N3906
Q185	151-0190-00			TRANSISTOR:SILICON, NPN	80009	151-0190-00
Q195	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q196	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q236	151-1025-00			TRANSISTOR:SILICON, JFE, N-CHANNEL	80009	151-1025-00
Q267	151-0190-00			TRANSISTOR:SILICON, NPN	80009	151-0190-00
Q268	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q269	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q281	151-0188-00			TRANSISTOR:SILICON, PNP	01295	2N3906
Q285	151-0190-00			TRANSISTOR:SILICON,NPN	80009	151-0190-00
Q410A,B	153-0586-00			TRANSISTOR:SILICON, PNP, PAIR	80009	153-0586-00
Q416	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q422A,B	153-0587-00			TRANSISTOR:SILICON, NPN, PAIR	80009	153-0587-00
Q426	151-0190-00			TRANSISTOR:SILICON, NPN	80009	151-0190-00
Q431	151-0190-00			TRANSISTOR:SILICON,NPN	80009	151-0190-00
Q443	151-0190-00			TRANSISTOR:SILICON, NPN	80009	151-0190-00
Q445	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q451	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q501	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q511	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q521	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Õ532	151-0190-00			TRANSISTOR:SILICON, NPN	80009	151-0190-00
Q540	151-0260-00			TRANSISTOR:SILICON, NPN	80009	151-0260-00
Q541	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
Q555	151-0220-00			TRANSISTOR:SILICON, PNP	80009	151-0220-00
ç558	151-0190-00			TRANSISTOR:SILICON, NPN	80009	151-0190-00
Q561A,B	151-0307-00			TRANSISTOR:SILICON, PNP, DUAL	07263	SP13404
Õ571A,B	151-0309-00			TRANSISTOR:SILICON, NPN, DUAL	07263	SP13378
Q910	151-0301-00	B010100	B040389	TRANSISTOR:SILICON, PNP	04713	2N2907A
0910	151-0302-00	в040390		TRANSISTOR:SILICON, NPN	04713	2N2222A
Q940	151-0334-00			TRANSISTOR:SILICON, NPN	80009	151-0334-00
0942	151-0190-00			TRANSISTOR:SILICON, NPN	80009	151-0190-00
õ990	151-0302-00	B010100	B040289	TRANSISTOR:SILICON, NPN	04713	2N2222A
Q990	151-0301-00			TRANSISTOR:SILICON, PNP	04713	2N2907A
R10	315-0105-00)		RES.,FXD,CMPSN:1M OHM,5%,0.25W	01121	CB1055
R20	315-0105-00			RES.,FXD,CMPSN:1M OHM,5%,0.25W	01121	CB1055
R30	322-0687-07			RES.,FXD,FILM:100.5M OHM,0.1%,0.25W	91637	MFF1421C10053B
R35	322-0624-07			RES., FXD, FILM: 990K OHM, 0.1%, 0.25W	75042	CEBTO-9903B
R38	321-0289-03			RES.,FXD,FILM:10K OHM,0.25%,0.125W	91637	MFF1816D10001C
R40	311-1566-00	во10100	в052349	RES., VAR, NONWIR: 200 OHM, 20%, 0.50W	73138	91A R200
R40	311-1568-00			RES., VAR, NONWIR:50 OHM, 20%, 0.50W	73138	91A R50
R50	322-0687-07			RES.,FXD,FILM:100.5M OHM,0.1%,0.25W	91637	MFF1421C10053B
	,					

	Taktroniy	Serial/Mod	al Ma		Mfr	
Ckt No.	Tektronix Part No.	Eff	Dscont	Name & Description	Code	Mfr Part Number
	322-0624-07			RES.,FXD,FILM:990K OHM,0.1%,0.25W	75042	CEBT0-9903B
R55 R58	321-0289-03			RES.,FXD,FILM:10K OHM,0.25%,0.125W	91637	MFF1816D10001C
R103	321-0289-03			RES., FXD, FILM:8.25K OHM, 1%, 0.125W	91637	MFF1816G82500F
R105	311-1559-00			RES., VAR, NONWIR: 10K OHM, 20%, 0.50W	73138	91A-10001M
R105	315-0107-00			RES.,FXD,CMPSN:100M OHM,5%,0.25W	01121	CB1075
RIOO	313-0107-00			,		
R107	315-0107-00			RES.,FXD,CMPSN:100M OHM,5%,0.25W	01121	CB1075
R120	315-0101-00			RES., FXD, CMPSN:100 OHM, 5%, 0.25W	01121	CB1015
R121	321-0222-00			RES.,FXD,FILM:2K OHM,1%,0.125W	91637	MFF1816G20000F
R122	315-0512-00			RES.,FXD,CMPSN:5.1K OHM,5%,0.25W	01121	CB5125
R125	321-0114-00			RES.,FXD,FILM:150 OHM,1%,0.125W	91637	MFF1816G150R0F
R126	321-0927-07			RES.,FXD,FILM:125 OHM,0.1%,0.125W	91637	MFF1816C125R0B
R131	321-0251-00			RES.,FXD,FILM:4.02K OHM,1%,0.125W	91637	MFF1816G40200F
R136	321-0199-06			RES.,FXD,FILM:1.15K OHM,0.25%,0.125W	91637	MFF1816C11500C
R137	315-0101-00			RES.,FXD,CMPSN:100 OHM,5%,0.25W	01121	CB1015
R141	321-0813-07			RES.,FXD,FILM:495 OHM,0.1%,0.125W	91637	MFF1816C495R0B
R145	321-0353-00			RES.,FXD,FILM:46.4K OHM,1%,0.125W	91637	MFF1816G46401F
R146	321-0272-00			RES., FXD, FILM: 6.65K OHM, 1%, 0.125W	91637	MFF1816G66500F
R152	321-0222-00			RES., FXD, FILM: 2K OHM, 1%, 0.125W	91637	MFF1816G20000F
R154	321-0225-06			RES.,FXD,FILM:2.15K OHM,0.25%,0.125	91637	MFF1816C21500C
R155	321-0154-00	B010100	в039999	RES.,FXD,FILM:392 OHM,1%,0.125W	91637	MFF1816G392R0F
R155	321-0150-00	в040000		RES.,FXD,FILM:357 OHM,1%,0.125W	91637	MFF1816G357R0F
R157	315-0564-00			RES., FXD, CMPSN: 560K OHM, 5%, 0.25W	01121	
R158	311-0467-00			RES., VAR, NONWIR: 100K OHM, 20%, 0.50W	11237	300SF-41334
R159	311-1567-00		в039999	RES., VAR, NONWIR: 100 OHM, 20%, 0.50W	73138	91A-100ROM
R159	311-1566-00	B040000		RES., VAR, NONWIR: 200 OHM, 20%, 0.50W	73138	91A R200
					01101	an 1 02 F
R160	315-0102-00			RES., FXD, CMPSN:1K OHM, 5%, 0.25W	01121	CB1025
R161	311-1564-00			RES., VAR, NONWIR:500 OHM, 20%, 0.50W	73138	91A R500
R162	321-0240-00			RES.,FXD,FILM:3.09K OHM,1%,0.125W	91637	
R164	321-0309-00			RES., FXD, FILM:16.2K OHM, 1%, 0.125W	91637 91637	MFF1816G16201F MFF1816G301R0F
R166	321-0143-00			RES.,FXD,FILM:301 OHM,1%,0.125W	91037	F# 1 1810G301R01
				DEG. EVD DILW.2 12 OUM 19 O 125W	91637	MFF1816G21000F
R168	321-0224-00			RES.,FXD,FILM:2.1K OHM,1%,0.125W RES.,FXD,CMPSN:1K OHM,5%,0.25W	01121	CB1025
R170	315-0102-00			RES.,FXD,FILM:7.5K OHM,1%,0.125W	91637	
R171	321-0277-00			RES.,FXD,FILM:3.01K OHM,0.1%,0.125W	91637	MFF1816C30100B
R173	321-0239-07			RES.,FXD,FILM:10K OHM,0.1%,0.125W	91637	MFF1816C10001B
R175	321-0289-07			RES., FAD, FILM: TOR OHM, 0:10, 0:125	2200.	
D176	221 1267 02			RES., FXD, FILM: 5.97K OHM, 0.25%, 0.125W	91637	MFF1816D59700C
R176	321-1267-03			RES., FXD, FILM:332 OHM, 1%, 0.125W	91637	MFF1816G332R0F
R177 R178	321-0147-00 311-1494-00			RES., VAR, NONWIR: 2K OHM, 10%, 1W	01121	11M187
R178 R181	321-0187-00			RES., FXD, FILM:866 OHM, 1%, 0.125W	91637	MFF1816G866R0F
	321-0107-00			RES., FXD, FILM:12.1K OHM, 1%, 0.125W	91637	MFF1816G12101F
R184	J21 J251-00					
R185	321-0236-00			RES.,FXD,FILM:2.8K OHM,1%,0.125W	91637	MFF1816G28000F
R191	311-1555-00			RES., VAR, NONWIR: 100K OHM, 20%, 0.5W	73138	91A R100K
R193	315-0474-00			RES., FXD, CMPSN:470K OHM, 5%, 0.25W	01121	CB4745
R195	315-0512-00			RES., FXD, CMPSN:5.1K OHM, 5%, 0.25W	01121	CB5125
R196	315-0682-00			RES.,FXD,CMPSN:6.8K OHM,5%,0.25W	01121	CB6825
10,50	010 0000 00					
R198	315-0472-00	, -		RES.,FXD,CMPSN:4.7K OHM,5%,0.25W	01121	
R205	311-1559-00			RES., VAR, NONWIR: 10K OHM, 20%, 0.50W	73138	
R206	315-0107-00			RES.,FXD,CMPSN:100M OHM,5%,0.25W		CB1075
R207	315-0107-00			RES.,FXD,CMPSN:100M OHM,5%,0.25W	01121	
R220	315-0101-00			RES.,FXD,CMPSN:100 OHM,5%,0.25W	01121	CB1015
R221	321-0222-00)		RES.,FXD,FILM:2K OHM,1%,0.125W	91637	
R222	315-0512-00			RES., FXD, CMPSN:5.1K OHM, 5%, 0.25W	01121	
R225	321-0114-00)		RES.,FXD,FILM:150 OHM,1%,0.125W	91637	MFF1816G150R0F

	T	0	-1 NI-		N.A.6.	
Ckt No.	Tektronix Part No.	Serial/Mod Eff	el No. Dscont	Name & Description	Mfr Code	Mfr Part Number
	203 2027 27			RES.,FXD,FILM:125 OHM,0.1%,0.125W	91637	MFF1816C125R0B
R226	321-0927-07			RES.,FXD,FILM:4.02K OHM,1%,0.125W	91637	MFF1816G40200F
R231	321-0251-00			RES.,FXD,FILM:1.15K OHM,0.25%,0.125W	91637	MFF1816C11500C
R236	321-0199-06 315-0101-00			RES.,FXD,CMPSN:100 OHM,5%,0.25W	01121	CB1015
R237				RES.,FXD,FILM:495 OHM,0.1%,0.125W	91637	MFF1816C495R0B
R241	321-0813-07			RES., FAD, FILM: 455 OHM, O. 18, O. 125h	31037	111 10100 1301103
R245	321-0353-00			RES., FXD, FILM: 46.4K OHM, 1%, 0.125W	91637	MFF1816G46401F
R246	321-0272-00			RES., FXD, FILM: 6.65K OHM, 1%, 0.125W	91637	MFF1816G66500F
R252	321-0222-00			RES.,FXD,FILM:2K OHM,1%,0.125W	91637	MFF1816G20000F
R254	321-0225-06			RES., FXD, FILM: 2.15K OHM, 0.25%, 0.125	91637	MFF1816C21500C
R255	321-0154-00	во10100	в039999	RES.,FXD,FILM:392 OHM,1%,0.125W	91637	MFF1816G392R0F
R255	321-0150-00	B040000		RES.,FXD,FILM:357 OHM,1%,0.125W	91637	MFF1816G357R0F
R260	315-0102-00			RES.,FXD,CMPSN:1K OHM,5%,0.25W	01121	CB1025
R261	321-0239-00			RES.,FXD,FILM:3.01K OHM,1%,0.125W	91637	MFF1816G30100F
R262	321-0240-00			RES.,FXD,FILM:3.09K OHM,1%,0.125W	91637	MFF1816G30900F
R264	321-0309-00			RES., FXD, FILM:16.2K OHM, 1%, 0.125W	91637	MFF1816G16201F
2000	221 0142 00			RES.,FXD,FILM:301 OHM,1%,0.125W	91637	MFF1816G301R0F
R266	321-0143-00			RES.,FXD,FILM:2.1K OHM,1%,0.125W	91637	MFF1816G21000F
R268	321-0224-00			RES., VAR, NONWIR:500 OHM, 20%, 0.50W	73138	91A R500
R269	311-1564-00				01121	CB1025
R270	315-0102-00			RES.,FXD,CMPSN:1K OHM,5%,0.25W RES.,FXD,FILM:7.5K OHM,1%,0.125W	91637	MFF1816G75000F
R271	321-0277-00			RES., PAD, PIEM. 1.3K OMA, 18, 0.123	3200	
R273	321-0239-07			RES.,FXD,FILM:3.01K OHM,0.1%,0.125W	91637	MFF1816C30100B
R275	321-0289-07			RES., FXD, FILM: 10K OHM, 0.1%, 0.125W	91637	MFF1816C10001B
R276	321-1267-03			RES.,FXD,FILM:5.97K OHM,0.25%,0.125W	91637	MFF1816D59700C
R277	321-0147-00	1		RES., FXD, FILM: 332 OHM, 1%, 0.125W	91637	MFF1816G332R0F
R284	321-0314-00			RES., FXD, FILM:18.2K OHM, 1%, 0.125W	91637	MFF1816G18201F
				RES.,FXD,FILM:2.8K OHM,1%,0.125W	91637	MFF1816G28000F
R285	321-0236-00			RES.,FXD,FILM:2.6R OHM,1%,0.125W RES.,FXD,FILM:49.9 OHM,1%,0.125W	91637	MFF1816G49R90F
R410	321-0068-00				91637	MFF1816G49R90F
R411	321-0068-00			RES., FXD, FILM: 49.9 OHM, 1%, 0.125W	01121	CB6805
R414	315-0680-00			RES.,FXD,CMPSN:68 OHM,5%,0.25W	75042	CEBT0-2941F
R420	322-0238-00)		RES.,FXD,FILM:2.94K OHM,1%,0.25W	73042	00010 00111
R421	321-0274-00)		RES.,FXD,FILM:6.98K OHM,1%,0.125W	91637	MFF1816G69800F
R422	315-0470-00			RES., FXD, CMPSN: 47 OHM, 5%, 0.25W	01121	CB4705
R423	321-0068-00			RES.,FXD,FILM:49.9 OHM,1%,0.125W	91637	MFF1816G49R90F
R424	321-0068-00			RES., FXD, FILM: 49.9 OHM, 1%, 0.125W	91637	MFF1816G49R90F
R431	321-0182-00			RES., FXD, FILM: 768 OHM, 1%, 0.125W	91637	MFF1816G768R0F
m 4 2 2	221 0170 02			RES.,FXD,FILM:715 OHM,1%,0.125W	91637	MFF1816G715R0F
R433	321-0179-00			RES.,FXD,FILM:2.61K OHM,1%,0.125W	91637	
R435	321-0233-00			RES.,FXD,CMPSN:4.7K OHM,5%,0.25W	01121	CB4725
R436	315-0472-00 315-0222-00			RES., FXD, CMPSN: 2.2K OHM, 5%, 0.25W	01121	CB2225
R438				RES.,FXD,FILM:6.04K OHM,1%,0.125W		MFF1816G60400F
R440	321-0268-00	,		IND. IL ADIT TERROLO IN CHAPTON CONTROL OF THE CONT		
R441	321-0247-00)		RES., FXD, FILM: 3.65K OHM, 1%, 0.125W	91637	MFF1816G36500F
R442	321-0400-00)		RES., FXD, FILM:143K OHM, 1%, 0.125W	91637	MFF1816G14302F
R443	315-0471-00)		RES.,FXD,CMPSN:470 OHM,5%,0.25W	01121	CB4715
R445	321-0218-00)		RES.,FXD,FILM:1.82K OHM,1%,0.125W	91637	MFF1816G18200F
R446	315-0472-00)		RES.,FXD,CMPSN:4.7K OHM,5%,0.25W	01121	CB4725
R448	315-0222-00)		RES.,FXD,CMPSN:2.2K OHM,5%,0.25W	01121	CB2225
R450	315-0222-00			RES., FXD, CMPSN: 2.2K OHM, 5%, 0.25W	01121	CB2225
R450 R451	315-0101-00			RES.,FXD,CMPSN:100 OHM,5%,0.25W	01121	CB1015
R451 R460	321-0229-00			RES.,FXD,FILM:2.37K OHM,1%,0.125W	91637	MFF1816G23700F
R460 R461	311-1566-00			RES., VAR, NONWIR: 200 OHM, 20%, 0.50W	73138	91A R200
V#OT	311-1300-00	•				
R462	321-0229-00)		RES.,FXD,FILM:2.37K OHM,1%,0.125W	91637	MFF1816G23700F
R464	321-0204-00)		RES.,FXD,FILM:1.3K OHM,1%,0.125W	91637	MFF1816G13000F
R466	315-0105-00)		RES., FXD, CMPSN:1M OHM, 5%, 0.25W	01121	CB1055

	Taktroniy	Serial/Model No.		Mfr	
Ckt No.	Tektronix Part N o.	Eff Dscont	Name & Description	Code	Mfr Part Number
OKT NO.	Tartito.	LII DOCUME	Hame & Bescription		
R467	311-1057-00		RES., VAR, NONWIR: 5K OHM, 3%, 0.01%	80294	3540S-420-502
R468	311-0887-00		RES., VAR, NONWIR: 50K OHM, 10%, 0.50W	80009	311-0887-00
R470	316-0155-00		RES.,FXD,CMPSN:1.5M OHM,10%,0.25W	01121	CB1551
R471	316-0155-00		RES.,FXD,CMPSN:1.5M OHM,10%,0.25W	01121	CB1551
R474	321-0408-00		RES.,FXD,FILM:174K OHM,1%,0.125W	91637	MFF1816G17402F
R475	321-0408-00		RES.,FXD,FILM:174K OHM,1%,0.125W	91637	MFF1816G17402F
R473 R478	321-0303-00		RES., FXD, FILM:14K OHM, 1%, 0.125W	91637	MFF1816G14001F
R479	321-0303-00		RES.,FXD,FILM:14K OHM,1%,0.125W	91637	MFF1816G14001F
R480	321-0289-03		RES., FXD, FILM:10K OHM, 0.25%, 0.125W	91637	MFF1816D10001C
R481	321-0816-07		RES.,FXD,FILM:5K OHM,0.1%,0.125W	91637	MFF1816C50000B
			• •		
R482	321-0222-03		RES., FXD, FILM: 2K OHM, 0.25%, 0.125W	91637	
R483	321-0193-03		RES., FXD, FILM: 1K OHM, 0.25%, 0.125W	91637	
R484	321-0612-03		RES.,FXD,FILM:500 OHM,0.25%,0.125W	91637	MFF1816D500R0C
R485	321-0126-03		RES.,FXD,FILM:200 OHM,0.25%,0.125W	91637	MFF1816D200R0C
R486	321-0097-03		RES.,FXD,FILM:100 OHM,0.25%,0.125W	91637	MFF1816D100R0C
D407	221_0751_06		RES.,FXD,FILM:50 OHM,0.25%,0.125W	91637	MFF1816C50R00C
R487	321-0751-06		RES.,FXD,FILM:30 OHM,0.23*,0.125W	91637	LFF18D2OR00D
R488	321-0030-02		RES.,FXD,FILM:10 OHM,0.5%,0.125W	91637	MFF1816G10R00D
R489	321-0001-01		RES., FXD, CMPSN: 2.2M OHM, 10%, 0.25W	01121	CB2251
R501 R503	316-0225-00 321-0452-00		RES.,FXD,FILM:499K OHM,1%,0.125W	91637	MFF1816G49902F
1003	321 0132 00				
R504	321-0423-00		RES.,FXD,FILM:249K OHM,1%,0.125W	91637	
R507	315-0473-00		RES.,FXD,CMPSN:47K OHM,5%,0.25W	01121	
R510	315-0104-00		RES.,FXD,CMPSN:100K OHM,5%,0.25W	01121	
R511	315-0104-00		RES.,FXD,CMPSN:100K OHM,5%,0.25W	01121	
R513	315-0104-00		RES.,FXD,CMPSN:100K OHM,5%,0.25W	01121	CB1045
DE 21	316-0276-00		RES.,FXD,CMPSN:27M OHM,10%,0.25W	01121	CB2761
R521 R523	321-0423-00		RES.,FXD,FILM:249K OHM,1%,0.125W	91637	MFF1816G24902F
R523	321-0452-00		RES.,FXD,FILM:499K OHM,1%,0.125W	91637	MFF1816G49902F
R524	315-0104-00		RES.,FXD,CMPSN:100K OHM,5%,0.25W	01121	CB1045
R531	315-0224-00		RES.,FXD,CMPSN:220K OHM,5%,0.25W	01121	CB2245
,	011 011 00				
R532	315-0103-00		RES., FXD, CMPSN:10K OHM, 5%, 0.25W	01121	CB1035
R533	315-0103-00		RES., FXD, CMPSN:10K OHM, 5%, 0.25W	01121	CB1035
R536	315-0104-00		RES.,FXD,CMPSN:100K OHM,5%,0.25W	01121	CB1045
R541	315-0681-00		RES., FXD, CMPSN:680 OHM, 5%, 0.25W	01121	
R542	315-0103-00		RES.,FXD,CMPSN:10K OHM,5%,0.25W	01121	CB1035
2550	215 0692 00		RES.,FXD,CMPSN:6.8K OHM,5%,0.25W	01121	CB6825
R550	315-0682-00		RES., FXD, CMPSN:15K OHM, 5%, 0.25W		CB1535
R551 R552	315-0153-00		RES.,FXD,CMPSN:1K OHM,5%,0.25W		CB1025
R552 R553	315-0102-00 315-0682-00		RES.,FXD,CMPSN:6.8K OHM,5%,0.25W	01121	CB6825
R555	315-0302-00		RES.,FXD,CMPSN:3K OHM,5%,0.25W	01121	CB3025
1			,		
R558	315-0302-00		RES.,FXD,CMPSN:3K OHM,5%,0.25W		CB3025
R 5 59	315-0202-00		RES., FXD, CMPSN: 2K OHM, 5%, 0.25W	01121	
R560	315-0101-00		RES.,FXD,CMPSN:100 OHM,5%,0.25W	01121	
R561	315-0181-00		RES.,FXD,CMPSN:180 OHM,5%,0.25W	01121	
R562	307-0104-00		RES.,FXD,CMPSN:3.3 OHM,5%,0.25W	01121	CB33G5
R570	315-0101-00	,	RES.,FXD,CMPSN:100 OHM,5%,0.25W	01121	CB1015
R570 R571	315-0101-00		RES.,FXD,CMPSN:180 OHM,5%,0.25W	01121	CB1815
R571 R572	307-0104-00		RES., FXD, CMPSN:3.3 OHM, 5%, 0.25W		CB33G5
R900	304-0102-00		RES., FXD, CMPSN:1K OHM, 10%, 1W	01121	
R900	303-0511-00		RES., FXD, CMPSN:510 OHM, 5%, 1W	01121	GB5115
1,500	555 5511 00				
R900	303-0751 - 00		RES.,FXD,CMPSN:750 OHM,5%,1W	01121	GB7515
R910	315-0821-00		RES.,FXD,CMPSN:820 OHM,5%,0.25W	01121	CB8215
R915	321-0285-00)	RES.,FXD,FILM:9.09K OHM,1%,0.125W	91637	MFF1816G90900F

	Tektronix	Serial/Model No.		Mfr	
Ckt No.	Part No.	Eff Dscont	Name & Description	Code	Mfr Part Number
R916	321-0268-00)	RES., FXD, FILM: 6.04K OHM, 1%, 0.125W	91637	MFF1816G60400F
R920	307-0105-00		RES., FXD, CMPSN:3.9 OHM, 5%, 0.25W	01121	CB39G5
R922	315-0100-00		RES., FXD, CMPSN:10 OHM, 5%, 0.25W	01121	CB1005
R924	315-0100-00		RES., FXD, CMPSN:10 OHM, 5%, 0.25W	01121	CB1005
R926	315-0100-00)	RES.,FXD,CMPSN:10 OHM,5%,0.25W	01121	CB1005
R940	307-0106-00)	RES.,FXD,CMPSN:4.7 OHM,5%,0.25W	01121	CB47G5
R942	315-0102-00)	RES., FXD, CMPSN:1K OHM, 5%, 0.25W		CB1025
R945	315-0912-00)	RES.,FXD,CMPSN:9.1K OHM,5%,0.25W		CB9125
R946	315-0622-00)	RES.,FXD,CMPSN:6.2K OHM,5%,0.25W		CB6225
R964	315-0100-00)	RES.,FXD,CMPSN:10 OHM,5%,0.25W	01121	CB1005
R966	315-0100-00)	RES.,FXD,CMPSN:10 OHM,5%,0.25W	01121	
R968	315-0100-00)	RES., FXD, CMPSN:10 OHM, 5%, 0.25W	01121	
R970	307-0105-00)	RES.,FXD,CMPSN:3.9 OHM,5%,0.25W		CB39G5
R989	315-0561-00)	RES., FXD, CMPSN: 560 OHM, 5%, 0.25W	01121	
R990	315-0102-00)	RES.,FXD,CMPSN:1K OHM,5%,0.25W	01121	CB1025
R994	321-0306-00)	RES.,FXD,FILM:15K OHM,1%,0.125W	91637	
R995	321-0306-00)	RES.,FXD,FILM:15K OHM,1%,0.125W	91637	MFF1816G15001F
RT104	307-0181-00)	RES.,THERMAL:100K OHM,10%,4MW/DEG C	15454	
RT203	307-0181-00)	RES.,THERMAL:100K OHM,10%,4MW/DEG C	15454	1DE104-K-220EC
S10A,B	260-1207-00)	SWITCH, PUSH: 2 MODULE	80009	
S20A,B	260-1207-00		SWITCH, PUSH: 2 MODULE	80009	
S30 S178	260-1209-00)	SWITCH, PUSH: 4PDT	80009	260-1209-00
S470	263-1004-00)	DRUM ASSY, CAM S:LF-3DB POINT	80009	263-1004-00
S480	263-1005-00)	DRUM ASSY, CAM S:GAIN	80009	
S490	263-1003-00)	DRUM ASSY, CAM S:HF-3DB POINT	80009	263-1003-00
U910	156-0067-00)	MICROCIRCUIT, LI: OPERATIONAL AMPLIFIER	80009	
บ990	156-0067-00)	MICROCIRCUIT, LI: OPERATIONAL AMPLIFIER	80009	156-0067-00
VR252	152-0168-00)	SEMICOND DEVICE: ZENER, 0.4W, 12V, 5%	04713	1N963B
VR438	152-0243-00		SEMICOND DEVICE: ZENER, 0.4W, 15V, 5%	81483	
VR448	152-0243-00)	SEMICOND DEVICE: ZENER, 0.4W, 15V, 5%	81483	
VR910	152-0123-00)	SEMICOND DEVICE: ZENER, 0.5W, 9V, 5%	80009	152-0123-00

and Number

DIAGRAMS AND CIRCUIT BOARD ILLUSTRATIONS

Symbols and Reference Designators

Electrical components shown on the diagrams are in the following units unless noted otherwise:

Values one or greater are in picofarads (pF). Capacitors =

Values less than one are in microfarads (μ F).

Resistors = Ohms (Ω) .

Graphic symbols and class designation letters are based on ANSI Standard Y32.2-1975.

Logic symbology is based on ANSI Y32.14-1973 in terms of positive logic. Logic symbols depict the logic function performed and may differ from the manufacturer's data.

The overline on a signal name indicates that the signal performs its intended function when it goes to the low state. Abbreviations are based on ANSI Y1.1-1972.

Other ANSI standards that are used in the preparation of diagrams by Tektronix, Inc. are:


Drafting Practices. Y14.15, 1966


Line Conventions and Lettering. Y14.2, 1973

Y10.5, 1968 Letter Symbols for Quantities Used in Electrical Science and

Electrical Engineering.


The following prefix letters are used as reference designators to identify components or assemblies on the diagrams.

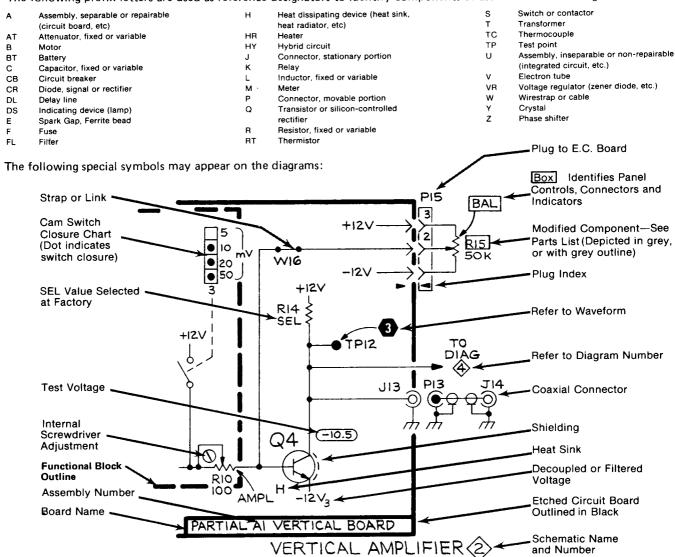

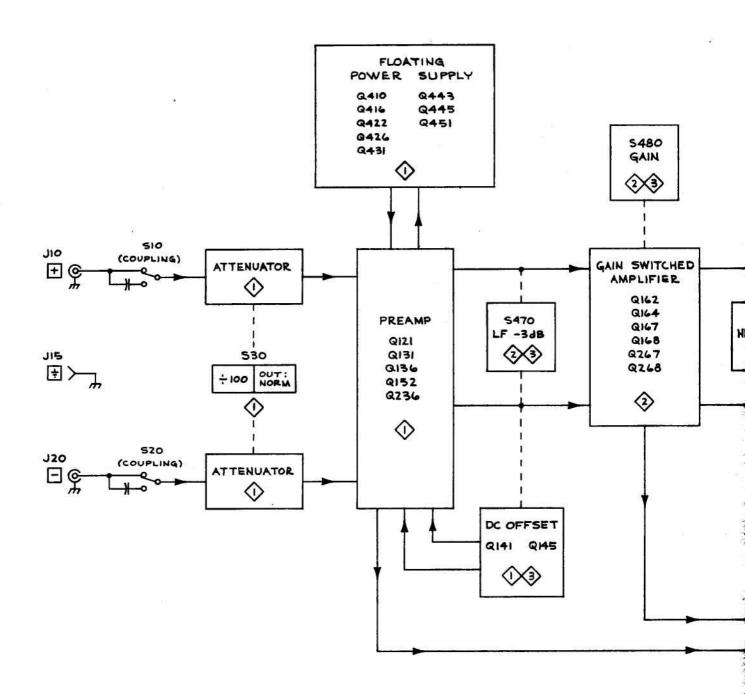




Digitally signed by http://www.aa4df.com

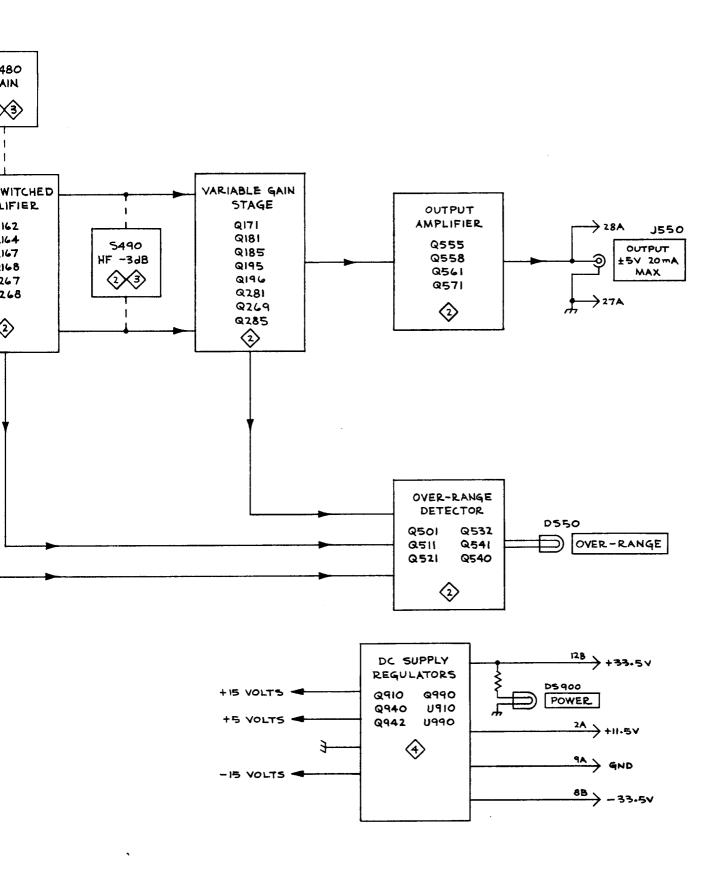


CKT NO	GRID LOC	CKT NO	GRID LOC	CKT NO	GRID LOC	CKT NO	GRID LOC	CKT NO	GRID LOC	CKT NO	GRID LOC
C28	КЗ	CR108	12	Q501	E1	R185	C2	R451	J4	R926	Н6
C30	К3	CR109	J3	Q511	F1	R191	В4	R460	F4	R940	C6
C32	K2	CR141	H4	Q521	E1	R193	C4	R461	G4	R942	E6
C34	K2	CR150	G4	Q532	C1	R195	C5	R462	F4	R945	E6
C35	L2	CR181	D2	Q540	В1	R196	B5	R464	H1	R946	F6
C38	K1	CR208	13	Q541	C1	R198	A5	R466	H1	R964	H6
C48	L3	CR209	J3	Q555	В3	R205	12	R470	F5	R966	Н6
C50	L3	CR231	H4	Q558	A2	R206	13	R471	E5	R968	E5
C54	M2	CR250	G4	Q561	C2	R207	13	R474	G5	R970	B5
C55	L3	CR281	D2	Q571	СЗ	R220	14	R475	F5	R989	B5
C58	L1	CR438	J2	Q910	В6	R221	14	R478	G6	R990	D5
C107	12	CR441	K4	Q940	В6	R222	H5	R479	F6	R994	C5
C122	13	CR448	К3	Q942	D6	R225	H5	R480	G2	R995	C5
C131	H4	CR531	D1	Q990	B5	R226	15	R481	G1	RT104	H2
C146	16	CR555	В3			R231	15	R482	G1	RT203	H2
C150	H5	CR558	B2	R10	К4	R236	G5	R483	G1		
C164	E3	CR560	B2	R20	K5	R237	G5	R484	G1	S10A	L4
C176	E2	CR561	B3	R30	K3	R241	J6	R485	G2	S10B	L4
C176	C5	CR570	B2	R35	K2	R245	J6	R486	G2 G2	S20A	L5
C207	12	CR570	B2	R38	K2	R246	Н6	R487	F1	S20B	L6
C222	15	C1137 1	02	R40	L2	R252	G5	R488	F1	S30	L1
C221	H5	F100	J3	R50	L2	R254	H5	R489	F2	S178	C3
	G5	F200	13	R55	L2	R255	G5	N409	12	S470 †	E5
C252		1 200	13	R58	L2	R260	D4	R501	D1	S480	G3
C264	E4	L164	F3	R103	H1	R261	E3	R503	D1	S490 †	B4
C276	E3	L264	F4	R105	12	R262	D4	R504	D1	3430 1	D-7
C431	J4	L416	J5	R106	12	R264	E4	R504	D1	U910	A6
C433	J4	L416	K5	R107	13	R266	E4	R510	F2	U990	A5
C435	K4	L445	K3	R120	14	R268	E4	R510	E1	0330	73
C438	J2	P40	K3	R121	14	R269	C4	R513	F1	VR252	G5
C443	L4	Q121	14	R122	Н4	R270	F3	R521	G3	VR438	J2
C448	К3	Q131	H4	R125	Н4	R271	E2	R523	E1	VR448	J3
C450	K4	Q136	H4	R126	14	R273	E2	R524	D1	VR910	C6
C470	F5	Q141	15	R131	13	R275	D3	R527	E1	V11310	CO
C471	E5	Q145	16	R136	G4	R276	D3	R532	В1		
C472	F6	Q152	G4	R137	G4	R277	D3	R533	C1		
C473	E6	Q162	D3	R141	J6	R284	D2	R536	C1		
C490	B4	Q164	E3	R145	J5	R285	C2	R541	В1		
C491	A4	Q167	F3	R146	H6	R410	J5	R542	A1		
C493 C494	B4 A4	Q168	E3	R152	G3	R411	J6	R550	B2		
C495 t	C2	Q171	F2	R154	нз	R414	J5	R551	A2		
C495 1	C2	0181	E2	R155	G4	R420	K5	R552	D2		
C497 †	C2	Q185	D2	R157	F1	R421	K5	R553	B2		
C498 †	B2	Q195	C4	R159	G5	R422	К5	R555	В3		
C499 †	B2	Q196	C5	R160	D3	R423	К5	R558	B2		
C530	F1	Q236	H4	R161	D3	R424	К6	R559	B2		
C536	C1	Q267	F3	R162	D3	R431	J4	R560	C2	1	
C552	C2	Q268	F3	R164	E3	R433	J4	R561	A2		
C560	B2	Q269	C5	R166	E3	R435	J4	R562	C2		
C570	B2	Q281	E2	R168	D4	R436	J3	R570	C2	1	
C910	B6	Q285	D2	R170	F2	R438	J2	R571	В3		
C922	D5	Q410	J5	R171	E2	R440	K4	R572	C3		
C922	H6	Q416	J4	R173	E2	R441	K4	R900	Α4		
C924	H6	Q422	J5	R175	D2	R442	K4	R910	D6		
C920	B5	Q426	K5	R176	D3	R443	K4	R915	C5		
C930	D6	Q431	J4	R177	D2	R445	К3	R916	D6		
	D6 С5	Q443	L4	R178	E3	R446	J3	R920	C5		
C960 C964	G6	Q445	K4	R181	E2	R448	J2	R922	E 5		
1.3904	GU			R184	C1	R450	J4	R924	Н6	1	
C966	G6	Q451	J4	IN 104	C I	111750		11327		1	

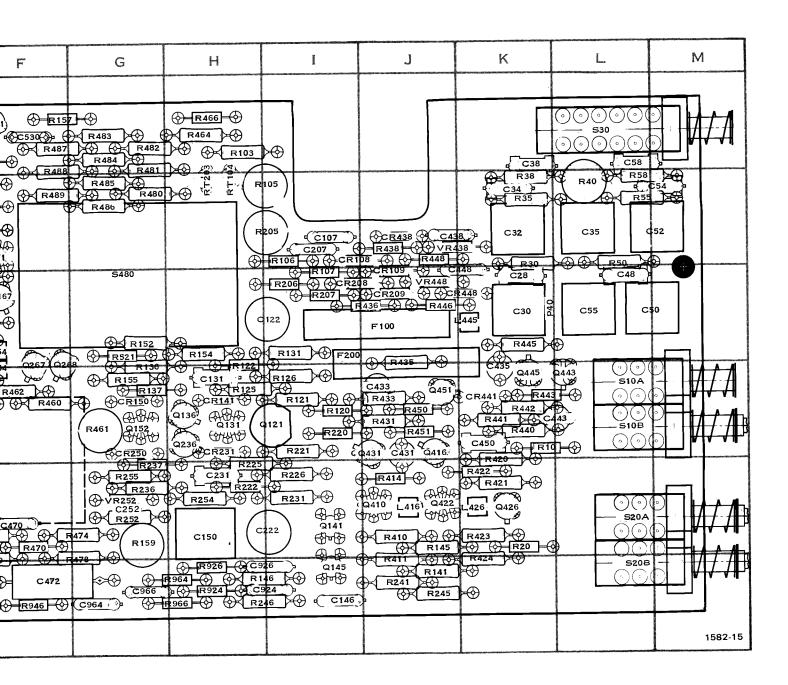

A1 Main Circuit Board 'A' side.

t Back of board.



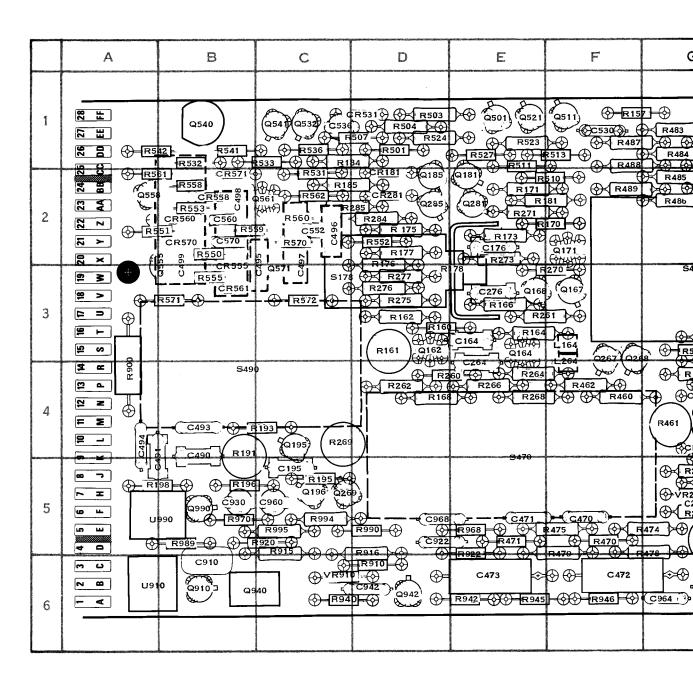

A1 Main Circuit



VOLTAGE AND WAVEFORM CONDITIONS

WARNING

Dangerous potentials exist at several points throughout this instrument. When the instrument is operated with the covers removed, do not touch exposed connections or components. Some transistors have voltages present on their cases. Disconnect the power source before replacing parts.


The voltages and waveforms shown on the diagram 1 and 2 were taken with the AM 502 front panel controls set as follows:

WAVEFORMS^a **VOLTAGES** 100 GAIN 100 **GAIN** pushbutton out ÷100 pushbutton out ÷100 HF -3 dB switch 1 MHz HF -3 dB switch 1 MHz LF -3 dB switch DC DC LF -3 dB switch pushbutton out + input AC + input AC pushbutton out pushbutton out **GND GND** pushbutton in pushbutton out pushbutton out - input AC - input AC pushbutton in **GND** pushbutton in **GND** STEP GAIN centered STEP GAIN centered DC BAL DC BAL DC OFFSET centered DC OFFSET centered (COURSE-FINE) (COURSE - FINE)

The waveforms shown were taken with a 50 mV square-wave input signal applied to the + input connector of the AM 502.

Voltage Conditions. The voltages shown on the diagram were obtained using a digital multimeter with a 10 megohm input impedance (TEKTRONIX DM 501 Digital Multimeter or TEKTRONIX 7D13 Digital Multimeter used with readout equipped, 7000-series oscilloscope).


Waveform Conditions. The waveforms shown are actual waveform photographs taken with a Tektronix Oscilloscope Camera System and Projected Graticule. Vertical deflection factor shown on the waveform is the actual deflection factor from the probe tip. Voltages and waveforms on the diagrams are not absolute and may vary between instruments because of component tolerances, internal calibration, or front-panel settings. Readouts are simulated in larger-than-normal type.

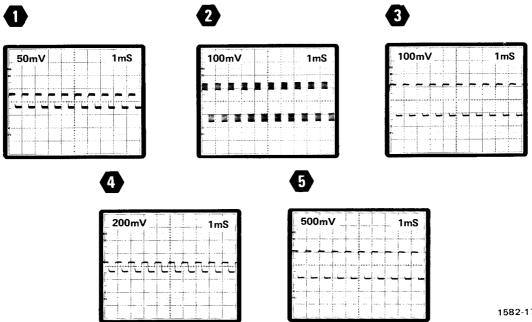



^{*}Ground reference: center horizontal graticule line


VOLTAGE & WAVEFORM CONDITIONS

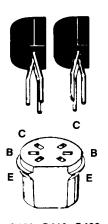

Q152, Q410, Q422





Q416, Q426, Q431 Q443, Q445, Q451

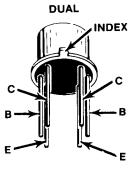

Q121, Q131, Q141, Q145



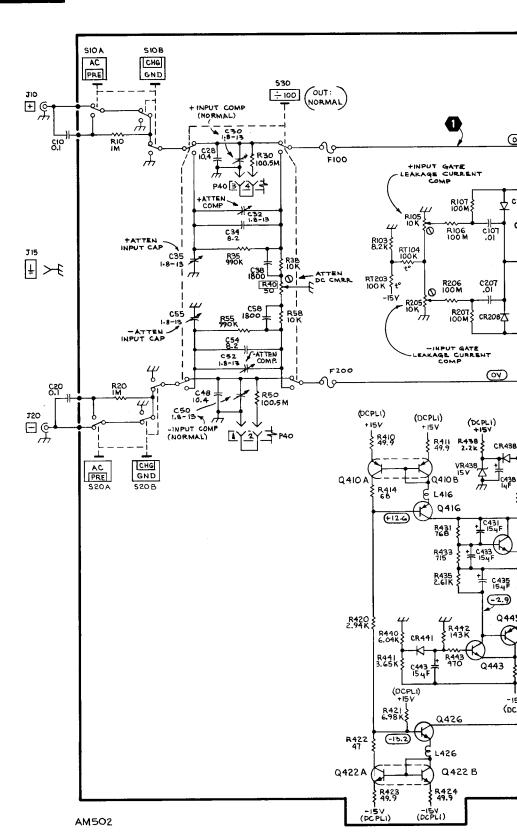

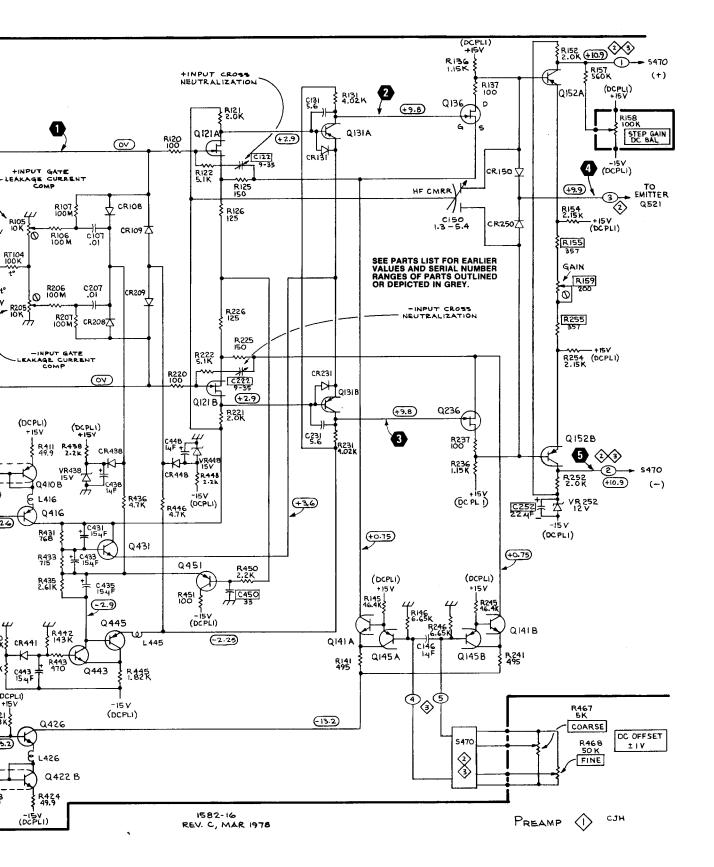

VOLTAGE AND WAVEFORM CONDITIONS

WARNING

Dangerous potentials exist at several points throughout this instrument. When the instrument is operated with the covers removed, do not touch exposed connections or components. Some transistors have voltages present on their cases. Disconnect the power source before replacing parts.

The voltages and waveforms shown on the diagram 1 and 2 were taken with the AM 502 front panel controls set as follows:

VOLIAGE		17777 = 1 0 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1 1				
GAIN	100	GAIN	100			
÷100	pushbutton out	÷100	pushbutton out			
HF -3 dB switch	1 MHz	HF -3 dB switch	1 MHz			
LF -3 dB switch	DC	LF -3 dB switch	DC			
+ input AC	pushbutton out	+ input AC	pushbutton out			
GND	pushbutton in	GND	pushbutton out			
- input AC	pushbutton out	input AC	pushbutton out			
GND	pushbutton in	GND	pushbutton in			
STEP GAIN	centered	STEP GAIN	centered			
DC BAL		DC BAL				
DC OFFSET	centered	DC OFFSET	centered			
(COURSE - FINE)		(COURSE- FINE)				


WAVEFORMS^a

VOLTAGES

The waveforms shown were taken with a 50 mV square-wave input signal applied to the \pm input connector of the AM 502.

Voltage Conditions. The voltages shown on the diagram were obtained using a digital multimeter with a 10 megohm input impedance (TEKTRONIX DM 501 Digital Multimeter or TEKTRONIX 7D13 Digital Multimeter used with readout equipped, 7000-series oscilloscope).


Waveform Conditions. The waveforms shown are actual waveform photographs taken with a Tektronix Oscilloscope Camera System and Projected Graticule. Vertical deflection factor shown on the waveform is the actual deflection factor from the probe tip. Voltages and waveforms on the diagrams are not absolute and may vary between instruments because of component tolerances, internal calibration, or front-panel settings. Readouts are simulated in larger-than-normal type.

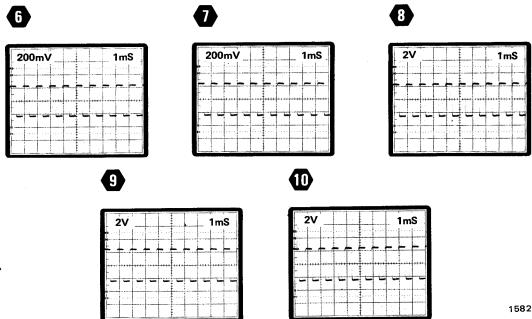





^aGround reference: center horizontal graticule line

VOLTAGE & WAVEFORM CONDITIONS

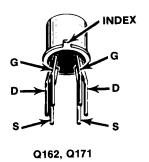

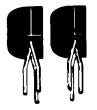


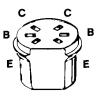



Q164, Q264




Q167, Q168, Q181, Q185 Q195, Q196, Q267, Q268 Q269, Q281, Q285

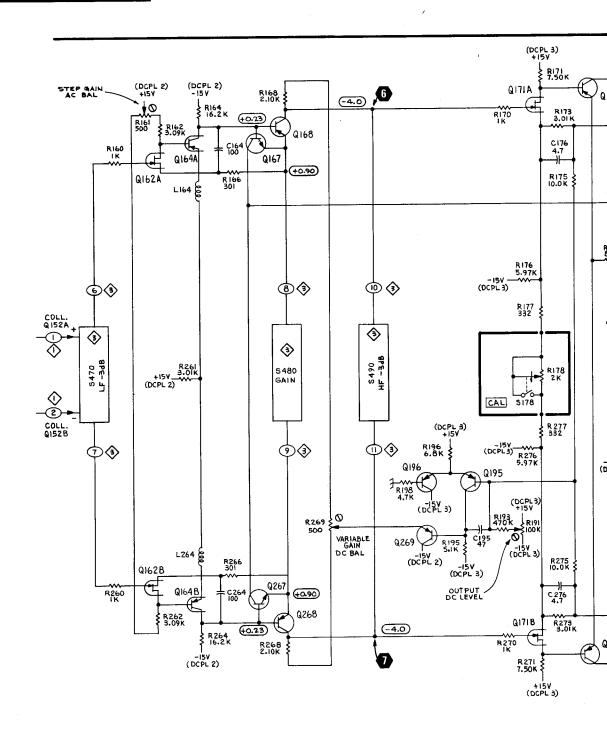

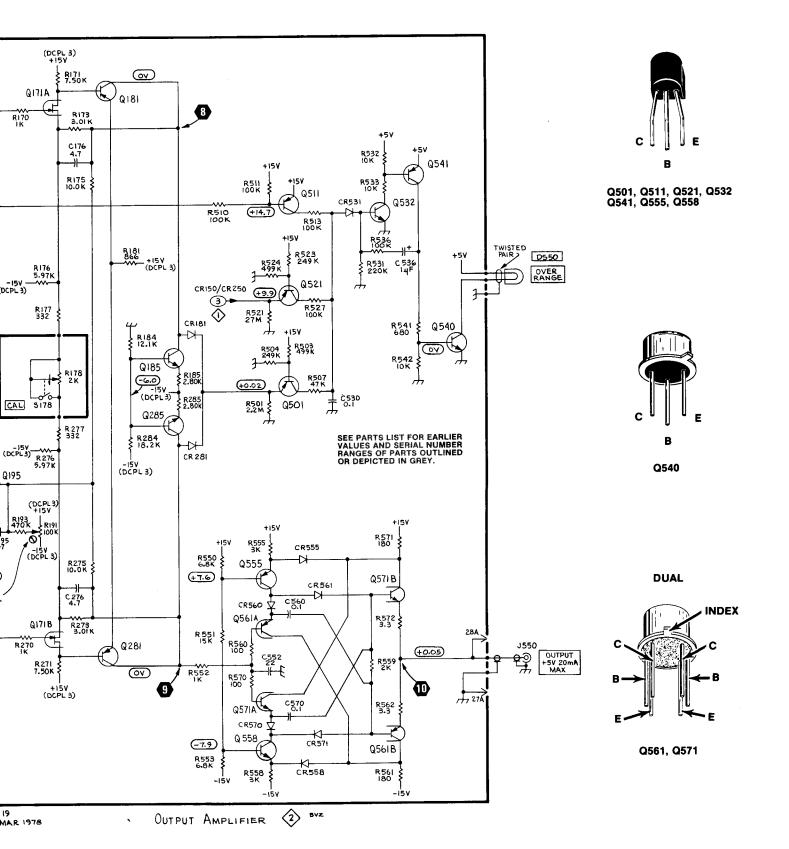

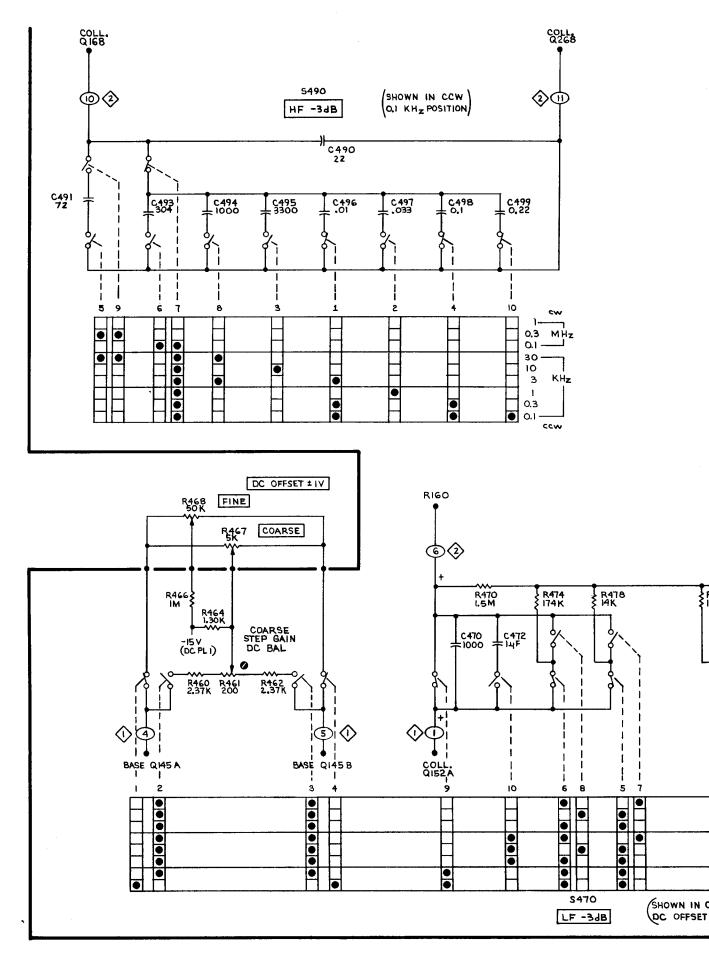


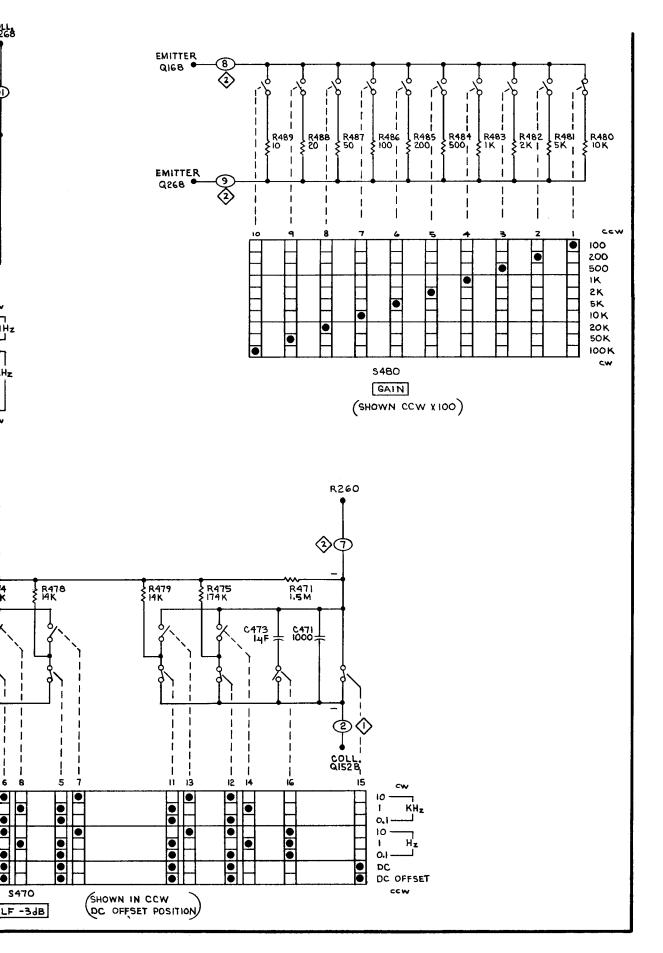


0 AR 1978 SWITCH DETAILS (3) CJH




VOLTAGE CONDITIONS

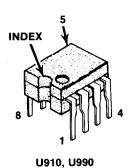


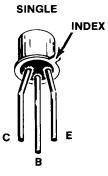

Q910, Q990

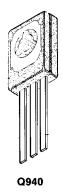




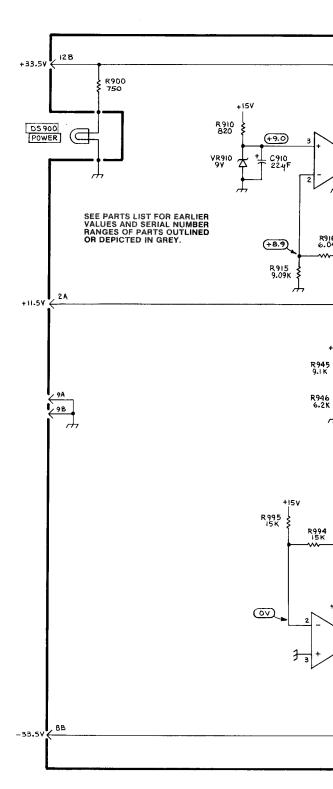

AM 502



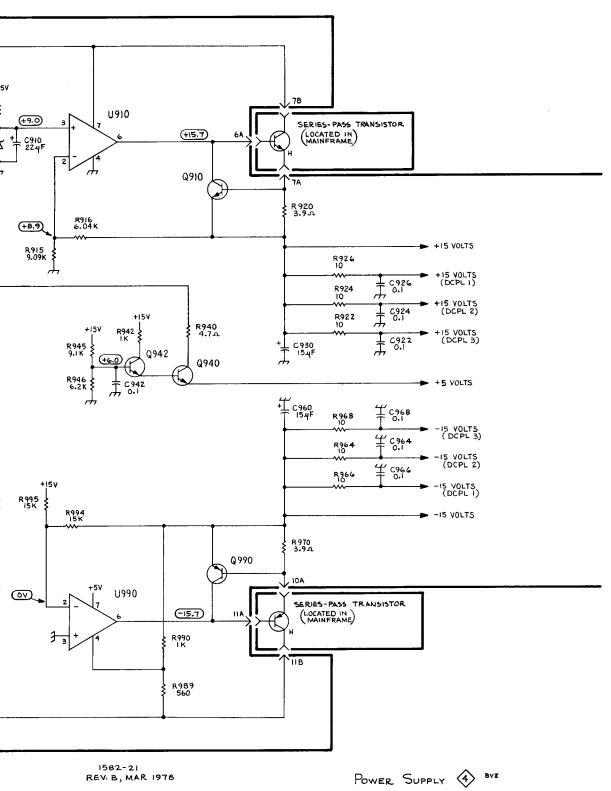

REPLACEABLE MECHANICAL PARTS

PARTS ORDERING INFORMATION

Replacement parts are available from or through your local Tektronix, Inc. Field Office or representative.

Changes to Tektronix instruments are sometimes made to accommodate improved components as they become available, and to give you the benefit of the latest circuit improvements developed in our engineering department. It is therefore important, when ordering parts, to include the following information in your order: Part number, instrument type or number, serial number, and modification number if applicable.

If a part you have ordered has been replaced with a new or improved part, your local Tektronix, Inc. Field Office or representative will contact you concerning any change in part number.

Change information, if any, is located at the rear of this manual.

SPECIAL NOTES AND SYMBOLS

X000 Part first added at this serial number

00X Part removed after this serial number

FIGURE AND INDEX NUMBERS

Items in this section are referenced by figure and index numbers to the illustrations.

INDENTATION SYSTEM

This mechanical parts list is indented to indicate item relationships. Following is an example of the indentation system used in the description column.

1 2 3 4 5

Name & Description

Assembly and/or Component
Attaching parts for Assembly and/or Component

Detail Part of Assembly and/or Component Attaching parts for Detail Part

---*---

Parts of Detail Part Attaching parts for Parts of Detail Part

Attaching Parts always appear in the same indentation as the item it mounts, while the detail parts are indented to the right. Indented items are part of, and included with, the next higher indentation. The separation symbol - - - * - - - indicates the end of attaching parts.

Attaching parts must be purchased separately, unless otherwise specified.

ITEM NAME

In the Parts List, an Item Name is separated from the description by a colon (:). Because of space limitations, an Item Name may sometimes appear as incomplete. For further Item Name identification, the U.S. Federal Cataloging Handbook H6-1 can be utilized where possible.

ABBREVIATIONS

# ACTR ADPTR ALIGN AL ASSEM ASSY ATTEN AWG	INCH NUMBER SIZE ACTUATOR ADAPTER ALIGNMENT ALUMINUM ASSEMBLED ASSEMBLY ATTENUATOR AMERICAN WIRE GAGE	ELCTRN ELEC ELCTLT ELEM EPL EQPT EXT FIL FLEX FLH	ELECTRON ELECTRICAL ELECTROLYTIC ELEMENT ELECTRICAL PARTS LIST EQUIPMENT EXTERNAL FILLISTER HEAD FLEXIBLE FLAT HEAD		INCH INCANDESCENT INSULATOR INTERNAL LAMPHOLDER MACHINE MECHANICAL MOUNTING NIPPLE NOT WIRE WOUND	SHLD SHLDR SKT SL SLFLKG SLVG SPR	SINGLE END SECTION SEMICONDUCTOR SHIELD SHOULDERED SOCKET SLIDE SELF-LOCKING SPRING SPRING
BD BRKT BRS BRZ BSHG CAB CER CHAS CKT COMP CONN COV CPLG CRT	BOARD BRACKET BRASS BRONZE BUSHING CABINET CAPACITOR CERAMIC CHASSIS CIRCUIT COMPOSITION CONNECTOR COVER COUPLING CATHODE RAY TUBE	FLTR FR FSTNR FT FXD GSKT HDL HEX HEX HD HEX SOC HLCPS HLEXT HV IC	FILTER FRAME OF FRONT FASTENER FOOT FIXED GASKET HANDLE HEXAGONAL HEAD HEXAGONAL SOCKET HELICAL COMPRESSION HELICAL EXTENSION HIGH VOLTAGE INTEGRATED CIRCUIT INSIDE DIAMETER	OBD OD OVH PH BRZ PL PLSTC PN PNH PWR RES RGD RLF RTNR SCH	ORDER BY DESCRIPTION OUTSIDE DIAMETER OVAL HEAD PHOSPHOR BRONZE PLAIN OF PLATE PLASTIC PART NUMBER PAN HEAD POWER RECEPTACLE RESISTOR RIGID RELIEF RETAINER SOCKET HEAD	SQ SSTL SW T TERM THD THK TNSN TPG TRH V VAR W/ WSHR	SQUARE STAINLESS STEEL STEEL SWITCH TUBE TERMINAL THREAD THICK TENSION TAPPING TRUSS HEAD VOLTAGE VARIABLE WITH WASHER
DEG DWR	DEGREE DRAWER	IDENT IMPLR	IDENTIFICATION IMPELLER	SCOPE SCR	OSCILLOSCOPE SCREW	XFMR XSTR	TRANSFORMER TRANSISTOR


CROSS INDEX—MFR. CODE NUMBER TO MANUFACTURER


Mfr. Code	Manufacturer	Address	City, State, Zip
02107 05091 08261 12327 22526 45722 70276	SPARTA MANUFACTURING COMPANY TRI-ORDINATE CORPORATION SPECTRA-STRIP CORP. FREEWAY CORPORATION BERG ELECTRONICS, INC. USM CORP., PARKER-KALON FASTENER DIV. ALLEN MFG. CO.	ROUTE NO. 2, BOX 128 343 SNYDER AVENUE 7100 LAMPSON AVE. 9301 ALLEN DRIVE YOUK EXPRESSWAY P. O. DRAWER 570	DOVER, OH 44622 BERRELEY HEIGHTS, NJ 07922 GARDEN GROVE, CA 92642 CLEVELAND, OH 44125 NEW CUMBERLAND, PA 17070 CAMPBELLSVILLE, KY 42718 HARTFORD, CT 06101
71159 71785 73743	BRISTOL SOCKET SCREW, DIV. OF AMERICAN CHAIN AND CABLE CO., INC. TRW, CINCH CONNECTORS FISCHER SPECIAL MFG. CO.	P O BOX 2244, 40 BRISTOL ST. 1501 MORSE AVENUE 446 MORGAN ST.	WATERBURY, CT 06720 ELK GROVE VILLAGE, IL 60007 CINCINNATI, OH 45206
73803 74445	TEXAS INSTRUMENTS, INC., METALLURGICAL MATERIALS DIV. HOLO-KROME CO.	34 FOREST STREET 31 BROOK ST. WEST	ATTLEBORO, MA 02703 HARTFORD, CT 06110
78189 78471	ILLINOIS TOOL WORKS, INC. SHAKEPROOF DIVISION TILLEY MFG. CO.	ST. CHARLES ROAD 900 INDUSTRIAL RD.	ELGIN, IL 60120 SAN CARLOS, CA 94070
79136 79807 80009 83385	WALDES, KOHINOOR, INC. WROUGHT WASHER MFG. CO. TEXTRONIX, INC. CENTRAL SCREW CO. SEASTROM MFG. COMPANY, INC.	47-16 AUSTEL PLACE 2100 S. O BAY ST. P O BOX 500 2530 CRESCENT DR. 701 SONORA AVENUE	LONG ISLAND CITY, NY 11101 MILWAUKEE, WI 53207 BEAVERTON, OR 97077 BROADVIEW, IL 60153 GLENDALE, CA 91201
86928 91836 97464	KINGS ELECTRONICS CO., INC. INDUSTRIAL RETAINING RING CO.	40 MARBLEDALE ROAD 57 CORDIER ST.	TUCKAHOE, NY 10707 IRVINGTON, NJ 07111

Fig. & Index No.		Serial/Mod Eff	del No. Dscont	Qty	1 2 3 4 5	Name & Description	Mfr Code	Mfr Part Number
		-010100	D150504		SHLD, ELECTRICAL:S	IDE	80009	337-1399-01
1-1	337-1399-01		B150524	2 1	SHLD, ELECTRICAL:S		80009	337-1399-01
	337-1399-01	B130525			SHLD, ELECTRICAL:S		80009	337-1399-00
-2	337-1399-00	в130525		i	KNOB:GRAY		80009	366-0494-00
- 3	366-0494-00			1		0.125 INCH, HEX SOC STL	74445	OBD
	213-0153-00 366-1518-00	D010100	B040700	1	KNOB:GRAY	0.110	80009	366-1518-00
-4	213-0153-00	B010100	B040799	ī		0.125 INCH, HEX SOC STL	74445	OBD
	366-1084-00	B010100	B040799	1		0,225 2,	80009	366-1084-00
	213-0725-00	B040800	B049999	ī	SETSCREW:3-48 X	0.095 INCH, HEX SOC STL	74445	OBD
	366-1023-00	B040600	B043333	1	KNOB:GRAY	, , , , , , , , , , , , , , , , , , , ,	80009	366-1023-00
				ī		0.093 INCH L, HEX SOC	71159	OBD
-	213-0246-00	возоооо		ī	KNOB:GRAY		80009	366-1101-00
- 5	366-1101-00			ì		0.125 INCH, HEX SOC STL	74445	OBD
_	213-0153-00			ī	KNOB:RED		80009	366-1317-00
- 6	366-1317-00			1	SETSCREW-5-40 X	0.125 INCH, HEX SOC STL	74445	OBD
-	213-0153-00			î	KNOB: GRAY, RF ATTE		80009	366-1001-00
- 7	366-1001-00			2	SETSCREW:5-40 X	0.125 INCH, HEX SOC STL	74445	OBD
•	213-0153-00			2	KNOB:LIGHT GRAY		80009	366-1163-00
-8	366-1163-00 213-0153-00			1		0.125 INCH, HEX SOC STL	74445	OBD
	366-1257-11				PUSH BUTTON: GRAY-		80009	366-1257-11
- 9					PUSH BUTTON: GRAY-		80009	366-1257-12
-10	366-1257-12				PUSH BUTTON: GRAY-		80009	366-1489-31
-11	366-1489-31				FR, PUSH BUTTON: GR		80009	426-0681-00
-12	426-0681-00	D010100	D0E2E70V		DIN KNOB SECRG-0	094 OD X 0.120 INCH LONG	80009	214-1840-00
	214-1840-00	B010100	B053575A	1		05. 05 0 unun =	80009	366-1422-01
-13	366-1422-01			1			80009	
	366-1690-00	в053580		1	POST, BDG, ELEC: ASS	FMRT.V	80009	129-0103-00
-14	129-0103-00				POSI, BDG, EDEC: ASS	(ATTACHING PARTS)		
	030 0455 00			1	MITT DIATM HEY .O	25-28 X 0.375 INCH, BRASS	73743	3089-402
-15	210-0455-00			1	TERMINAL, LUG:0.25		86928	A313-136
-16	210-0223-00			_	121411111111111111111111111111111111111	*		
				- .	POST ASSEMBLY I	NCLUDES:		
	200-0103-00			1	NUT.PLAIN.KNURL	:0.25-28 X 0.375" OD,BRASS	80009	200-0103-00
	129-0077-00			. ī	. STUD SHOULDERED	:0.938 INCH LONG, BRASS	80009	129-0077-00
- 17	131-0955-00			2	CONNECTOR . RCPT . : B	NC, FEMALE, W/HARDWARE	05091	31-279
-17 -18	131-0933-00				CONNECTOR, RCPT, : B		91836	KC19-153BNC
-19	384-1221-00		B049999	1	EXTENSION SHAFT:0	.081 DIA X 3.0 INCH L,SST	80009	384-1221-00
-19	384-0415-01			ī	EXTENSION SHAFT:3	.18 L X 0.08 OD STL	80009	384-0415-01
-20		200000		1		:(SEE R468 EPL)		
20						(ATTACHING PARTS)		
-21	210-0583-00			1	NUT, PLAIN, HEX.: 0.	25-32 X 0.312 INCH,BRS	73743	2X20224-402
-22	210-0940-00			1	WASHER, FLAT: 0.25	ID X 0.375 INCH OD,STL	79807	OBD
						*		
-23				1	RESISTOR, VARIABLE			
						(ATTACHING PARTS)	72742	2X20224-402
-24	210-0583-00	ı		1	NUT, PLAIN, HEX.: 0.	25-32 X 0.312 INCH,BRS		
-25	210-0940-00			1	WASHER, FLAT: 0.25	ID X 0.375 INCH OD, STL	79807	OBD
						*	00000	250 0270 00
-26	358-0378-00	ı		2	BUSHING, SLEEVE: PF	ESS MOUNT	80009	358-0378-00
-27	358-0029-00	1		1	BSHG, MACH. THD: HEX	x,0.375-32 x 0.438"LONG	80009	358-0029-00
						(ATTACHING PARTS)	72742	2X28269-402
-28	210-0590-00)		1		375 X 0.438 INCH,STL	73743	
-29	210-0978-00	1		1	WASHER, FLAT: 0.375	D ID X 0.50 INCH OD, STL	78471	OBD
						* = -	80009	333-1711-00
-30	333-1711-00)		1		<u>2</u>		
-31	214-1513-01	B010100	в053579	1	LCH.PLUG-IN RET:		80009	
	105-0719-00	в053579)	1	LATCH, RETAINING:		80009	103-0719-00
				_		(ATTACHING PARTS)	45722	OBD
-32	213-0254-00			1		2-32 X 0.250,100 DEG,FLH	80009	105-0718-00
-33	105-0718-00	В053580)	1	RELEASE, LATCH:	* *	60009	703-0110-00
				_	DAGE TAMBUOTEEN		80009	200-0935-00
-34	200-0935-00			2	LENS, LIGHT: AMBER	0.29 OD X 0.19 CASE	80009	378-0602-01
- 35	378-0602-01	-		1	HEND, LIGHT : AMDER		2222	

Fig. & Index No.		Serial/Mo Eff	del No. Dscont	Qty	1 2 3 4 5	Name & Description	Mfr Code	Mfr Part Number
1 26	378-0602-00			1	LENS, LIGHT: GREEN		80009	378-0602-00
1 - 36	352-0157-00				LAMPHOLDER:WHITE	PLASTIC	80009	352-0157-00
-37	386-2529-00	B010100	BU19999	_	SUBPANEL, FRONT: PI		80009	386-2529-00
-38		B050000	B043333		SUBPANEL, FRONT: PI		80009	386-2529-01
	386-2529-01	возоооо		_	DODI AMELI, I NOMI 11.	(ATTACHING PARTS)		
-39	213-0229-00			4	SCR, TPG, THD FOR:	6-20 X0.375"100 DEG,FLH STL	83385	OBD
-40	337-1782-00	B010100	B049999	1	SHLD, ELECTRICAL:	FRONT SUBPANEL	80009	337-1782-00
-40	337-1782-01			1	SHLD, ELECTRICAL:		80009	337-1782-01
-41	214-1061-00	D030000		1	SPRING, GROUND:FL		80009	214-1061-00
-41 -42	426-0725-00			1	FR SECT, PLUG-IN:		80009	426-0725-00
-42	420-0725-00			_		(ATTACHING PARTS)		
-43	213-0146-00			2	SCR, TPG, THD FOR:	6-20 X 0.313 INCH,PNH STL	83385	
-44	386-3657-00	XB053580		2	SUPPORT, PLUG-IN:		80009	386-3657-00
-45	210-1270-00			2	WASHER, FLAT: 0.14	1 ID X 0.04 THK,AL,0.219 OD	80009	210-1270-00
-45 -46	426-0724-00	7000000		1	FR SECT, PLUG-IN:		80009	426-0724-00
-40	-20-0/24-00			-	,	(ATTACHING PARTS)		
-47	213-0146-00			1	SCR.TPG.THD FOR:	6-20 X 0.313 INCH, PNH STL	83385	OBD
-47 -48	211-0012-00			ī	SCREW, MACHINE: 4-	40 X 0.375 INCH, PNH STL	83385	OBD
-49	210-0406-00			ī	NUT.PLAIN.HEX.:4	-40 X 0.188 INCH,BRS	73743	2X12161-402
-49	210-0400-00					· = *		
-50	376-0039-00			1	ADPT.SHAFT.CPLG:	0.128 AND 0.082"DIA SHAFT	80009	376-0039-00
-50	213-0075-00			2	. SETSCREW:4-40	X 0.094 INCH, HEX SOC STL	70276	OBD
						-MAIN(SEE Al EPL)		
-51				_	. CKT BOARD ASSY			
				1	TRANSISTOR: (SE	E 0940 EPL)		
- 52				_	TRANSISTOR I			
				1	WASHER LOCK:	0.228 ID X 0.375 INCH OD, STL	78189	4704-04-02
- 53	210-1122-00			1	WASHER NONME	TAL:0.188 ID X 0.313" OD, TEFLON	02107	OBD
-54	210-0910-00			36	CONTACT, ELEC: C	KT CD SW.SPR	80009	131-0604-00
-55	131-0604-00	D010100	в030294	72	SOCKET DIN TER	M:0.188 INCH LONG	22526	75060
-56	136-0252-04			56	SOCKET PIN TER	M:0.188 INCH LONG	22526	75060
	136-0252-04			2	SOCKET PLUG IN	:MICROCIRCUIT,8 CONTACT	73803	C950802
	136-0514-00	B030292	,	1	COV TRANSISTOR	R:0.438 DIA X 0.47 INCH H	80009	200-0687-01
- 57	200-0687-01			1	SOCKET PLUG-TN	1:6 CONTACT, ROUND	71785	133-96-12-062
- 58	136-0235-00			1		ABLE: (SEE R467 EPL)		
- 59						(ATTACHING PARTS)		
-60	210-0583-00			1	. NUT, PLAIN, HEX.	:0.25-32 X 0.312 INCH,BRS	73743	
-61	210-0046-00			1	. WASHER, LOCK: IN	TL,0.26 ID X 0.40" OD,STL	78189	1214-05-00 - 0541C
-62	210-1025-00			2	. WASHER, FLAT: 0.	312 ID X 0.469 OD,BRS	12327	OBD
02						*		
- 63	407-1337-00			1	•	(ATTACHING PARTS)	80009	407-1337-00
-64	211-0008-00			2	. SCREW, MACHINE:	:4-40 X 0.25 INCH, PNH STL	83385	
- 65	210-0551-00					.:4-40 X 0.25 INCH,STL	83385	
-66	337-1883-00				. SHLD, ELECTRICA		80009	337-1883-00 260-1207-00
-67	260-1207-00				. SWITCH, PUSH:2		80009 80009	361-0385-00
-68	361-0385-00		во52349		. SPACER, PB SW:			361-0384 - 00
	361-0384-00				. SPACER, PB SW:		80009	
- 69	260-1209-00			1	. SWITCH, PUSH:41	PDT	80009	260-1209-00
- 70	361-0383-00			2	. SPACER, PB SW:	CHARCOAL, 0.33 INCH LONG	80009	361-0383-00
- 71	344-0154-00			4	. CLIP, ELECTRICA	AL:FOR 0.25 INCH DIA FUSE	80009	344-0154-00
-72	384-0173-00			1	EXTENSION SHAFT	:0.125 DIA X 6.438 INCH L	80009	384-0173-00
-73	384-0348-01			1		:0.125 DIA X 9.161 INCH L	80009	384-0348-01
, ,	263-1004-00			1		(ATTACHING PARTS)	80009	263-1004-00
-74	211-0116-00)		4		4-40 X 0.312 INCH, PNH BRS	83385	OBD
		•			. ACTUATOR ASSY		00000	200-1598-00
- 75	200-1598-00)		1	. COVER, CAM SW:	16 ELEMENTS	80009	5103-25-MD-R
- 76	354-0219-00			1	. RING, RETAININ	G:FOR 0.25 INCH SHAFT	19130	2102-52-km−k

Fig. & Index No.	Tektronix Part No.	Serial Eff	/Model No. Dscont	Ωtv		1 2 3 4 5	Name & Descripti	on	Mfr Code	Mfr Part Number
			5000111							
1-77	401-0155-0					. BEARING, CAM S			80009	401-0155-00
-7 8	131-0840-0			1		. CONTACT, ELEC:			80009	131-0840-00
- 79	214-1704-0					•	AM SW DETENT, 0.008		80009	214-1704-01
-80	214-1127-0	0		2			:0.125 DIA X 0.125	INCH L	80009	214-1127-00
-81	105-0498-0	0				 DRUM, CAM SWIT 			80009	
-82	401-0156-0	0		1		. BEARING, CAM S			80009	
- 83	210-0406-0	0		4			.:4-40 X 0.188 INCH	,BRS	73743	
	263-1003-0	0		1	Ι	DRUM ASSY,CAM S	:HF-3DB POINT (ATTACHING PARTS	5)	80009	263-1003-00
-84	211-0116-0	0		4	5	SCR, ASSEM WSHR:	4-40 X 0.312 INCH,PN	NH BRS	83385	OBD
		-		-		ACTUATOR ASSY	INCLUDES:			
-85	200-1597-0	0		1		. COVER, CAM SW:	10 ELEMENTS		80009	200-1597-00
-86	354-0219-0	0		1		. RING, RETAININ	G:FOR 0.25 INCH SHAP	T	79136	5103-25-MD-R
-87	401-0155-0	0		1		. BEARING, CAM S	W:FRONT		80009	401-0155-00
~ 88	131-1248-0	0		1		. CONTACT, ELEC:	GROUNDING		80009	131-1248-00
-89	214-1704-0	0					AM SW DETENT, 0.006		80009	214-1704-00
- 90	214-1704-0	1					AM SW DETENT, 0.008		80009	214-1704-01
- 91	214-1127-0	0		2		. ROLLER, DETENT	:0.125 DIA X 0.125]	NCH L	80009	214-1127-00
-92	105-0502-0	0		1		. DRUM, CAM SWIT	CH:		80009	105-0502-00
- 93	401-0156-0	0		1		. BEARING, CAM S			80009	401-0156-00
-94	210-0406-0	0		4		. NUT, PLAIN, HEX	.:4-40 X 0.188 INCH	BRS	73743	2X12161-402
	263-1005-0	0		1		DRUM ASSY,CAM S			80009	263-1005-00
- 95	211-0116-0	0		4	5	SCR, ASSEM WSHR:	4-40 X 0.312 INCH,PN		83385	OBD
		_		_		ACTUATOR ASSY	INCLUDES:			
- 96	200-1596-0			1			10 ELEMENT, 0.83 INCH (ATTACHING PARTS		80009	200-1596-00
-97	211-0008-0	0		4		. SCREW.MACHINE	:4-40 X 0.25 INCH,PN		83385	OBD
-98	210-0004-0						NTL,0.12 ID X 0.26"C			1204-00-00-0541C
-99	210-0406-0						.:4-40 X 0.188 INCH,		73743	2X12161-402
-100	131-0963-0			1		. CONTACT, ELEC:			80009	131-0963-00
100	131 0300 0	•		_		·	*			
-101	354-0391-0	0		1		. RING, RETAININ	G:0.395"FREE ID X 0.	.025" STL	97464	3100-43-CD
-102	401-0081-0	2		1		. BEARING, CAM S	W:FRONT		80009	401-0081 - 02
-103	214-1139-0	3		1		. SPRING, FLAT: R	ED COLORED		80009	214-1139 - 03
-104	214-1139-0	2		1.		. SPRING, FLAT: G	REEN COLORED		80009	214-1139-02
-105	214-1127-0	0		2		. ROLLER, DETENT	:0.125 DIA X 0.125 I	NCH L	80009	214-1127-00
-106	105-0503-0	0		1		. DRUM, CAM SWIT	CH:		80009	105-0503-00
-107	401-0115-0	0		1		. BEARING, CAM S			80009	401-0115-00
	210-0406-0	0		3		. NUT, PLAIN, HEX	.:4-40 X 0.188 INCH,	BRS	73743	2X12161-402
	131-0963-0	0		1		. CONTACT, ELEC:	GROUNDING		80009	131-0963-00
-108	384-1220-0	0		1	E	EXTENSION SHAFT	:0.125 DIA X 8.038 I	NCH L	80009	384-1220-00
-109	376-0029-0	0		1	C	CPLG, SHAFT, RGD:	0.128 ID X 0.312 OD	X 0.5"L	80009	376-0029-00
	213-0075-0	0		2		SETSCREW:4-40	X 0.094 INCH, HEX SO	C STL	70276	OBD
-110				1	F	RESISTOR, VARIAB	LE:(SEE R178/S178 EF (ATTACHING PARTS			
-111	210-0583-0	0		1	1	UT.PLAIN.HEX.:	0.25-32 X 0.312 INCH		73743	2X20224-402
-112	210-0046-0			1			L,0.26 ID X 0.40" OF		78189	1214-05-00-0541C
-113	407-0803-0			ī		BRACKET, RES MTG			80009	407-0803-00
•••	185 0005 0	0		7 15		ATOM DIRONDICTI			08261	OBD
	175-0825-0					WIRE, ELECTRICAL				
-115						WIRE, ELECTRICAL		· DE		175-0826-00
-116				4			:0.48" L,22-26AWG WI	KE	22526	
-117				1		CONN BODY, PL, EL			80009	352-0162-02
-118	179-2028-0	0		1	V	WIRING HARNESS,	:MAIN		80009	159-2028-00

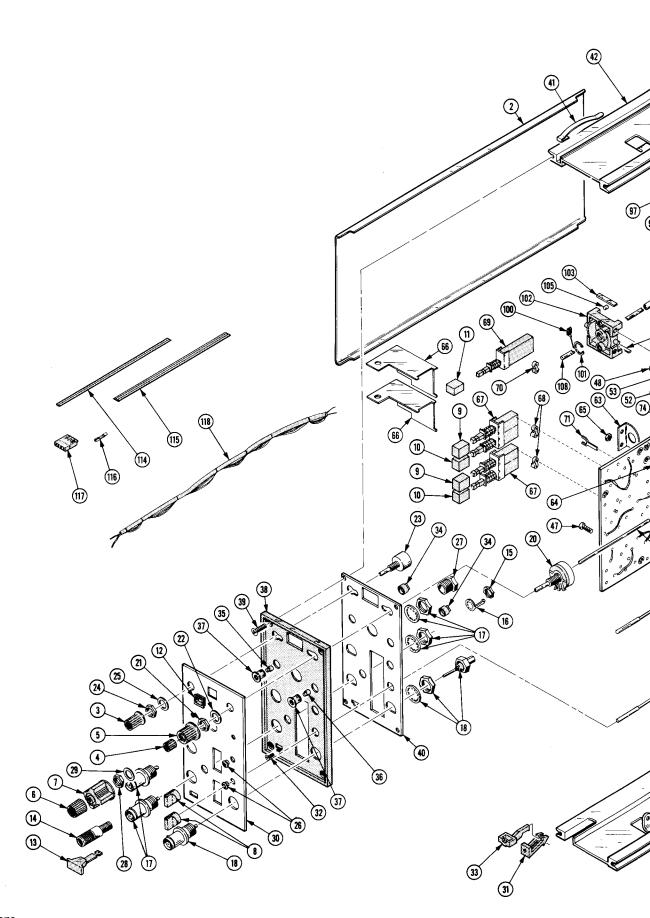

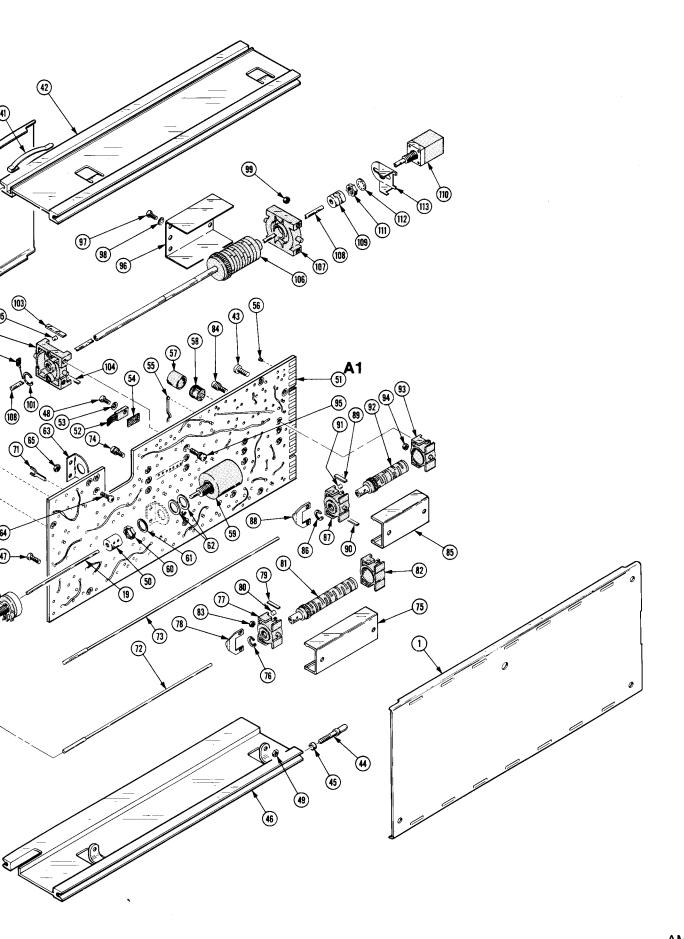

ACCESSORIES

STANDARD ACCESSORIES

Fig. & Index No.	Tektronix Part No.	Serial/N Eff	Model No. Dscont	Qty	1 2 3 4 5	Name & Description	Mfr Code	Mfr Part Number
	070-1582-0	1		1	MANUAL, TECH: IN	STRUCTION	80009	070-1582-01

MANUAL CHANGE INFORMATION

At Tektronix, we continually strive to keep up with latest electronic developments by adding circuit and component improvements to our instruments as soon as they are developed and tested.

Sometimes, due to printing and shipping requirements, we can't get these changes immediately into printed manuals. Hence, your manual may contain new change information on following pages.

A single change may affect several sections. Since the change information sheets are carried in the manual until all changes are permanently entered, some duplication may occur. If no such change pages appear following this page, your manual is correct as printed.

SERVICE NOTE

Because of the universal parts procurement problem, some electrical parts in your instrument may be different from those described in the Replaceable Electrical Parts List. The parts used will in no way alter or compromise the performance or reliability of this instrument. They are installed when necessary to ensure prompt delivery to the customer. Order replacement parts from the Replaceable Electrical Parts List.

CALIBRATION TEST EQUIPMENT REPLACEMENT

Calibration Test Equipment Chart

This chart compares TM 500 product performance to that of older Tektronix equipment. Only those characteristics where significant specification differences occur, are listed. In some cases the new instrument may not be a total functional replacement. Additional support instrumentation may be needed or a change in calibration procedure may be necessary.

Comparison of Main Characteristics

	Comparison of Main Character	
DM 501 replaces 7D13		
PG 501 replaces 107	PG 501 - Risetime less than 3.5 ns into 50 Ω.	107 - Risetime less than 3.0 ns into 50 Ω.
108	PG 501 - 5 V output pulse; 3.5 ns Risetime.	108 - 10 V output pulse; 1 ns Risetime.
111	PG 501 - Risetime less than 3.5 ns; 8 ns Pretrigger pulse delay.	111 - Risetime 0.5 ns; 30 to 250 ns Pretrigger Pulse delay.
114	PG 501 - ±5 V output.	114 - ±10 V output. Short proof output.
115	PG 501 - Does not have Paired, Burst, Gated,	115 - Paired, Burst, Gated, and Delayed
	or Delayed pulse mode; ±5 V dc	pulse mode; ±10 V output.
	Offset. Has ±5 V output.	Short-proof output.
PG 502 replaces 107		108 - 10 V output.
108 111	PG 502 - 5 V output PG 502 - Risetime less than 1 ns; 10 ns	111 - Risetime 0.5 ns; 30 to 250 ns
Uni	Pretrigger pulse delay.	Pretrigger pulse delay.
114	PG 502 - ±5 V output	114 - ±10 V output, Short proof output.
115	PG 502 - Does not have Paired, Burst, Gated,	115 - Paired, Burst, Gated, Delayed & Un-
	Delayed & Undelayed pulse mode;	delayed pulse mode; ±10 V output.
2121	Has ±5 V output.	Short-proof output. 2101 - Paired and Delayed pulse; 10 V
2101	PG 502 - Does not have Paired or Delayed pulse. Has ±5 V output.	output.
PG 506 replaces 106	PG 506 - Positive-going trigger output signal	106 - Positive and Negative-going trigger
, a ooo replaces 100	at least 1 V; High Amplitude out-	output signal, 50 ns and 1 V; High
	put, 60 V.	Amplitude output, 100 V.
067-0502-01	PG 506 - Does not have chopped feature.	0502-01 - Comparator output can be alter-
		nately chopped to a reference voltage.
		voltage.
SG 503 replaces 190,	CO FOO Amelitude reace 5 mV to 5 5 V a a	190B - Amplitude range 40 mV to 10 V p-p.
190A, 190B	SG 503 - Amplitude range 5 mV to 5.5 V p-p. SG 503 - Frequency range 250 kHz to 250 MHz.	191 - Frequency range 350 kHz to 100 MHz.
067-0532-01	SG 503 - Frequency range 250 kHz to 250 MHz.	0532-01 - Frequency range 65 MHz to 500 MHz.
TG 501 replaces 180,		
180A	TG 501 - Marker outputs, 5 sec to 1 ns.	180A - Marker outputs, 5 sec to 1 μs.
	Sinewave available at 5, 2, and 1 ns.	Sinewave available at 20, 10,
	Trigger output - slaved to marker	and 2 ns. Trigger pulses 1, 10, 100 Hz; 1, 10, and 100 kHz.
	output from 5 sec through 100 ns. One time-mark can be generated at a	Multiple time-marks can be
	time.	generated simultaneously.
181	TG 501 - Marker outputs, 5 sec to 1 ns. Sine-	181 - Marker outputs, 1, 10, 100, 1000,
404	wave available at 5, 2, and 1 ns. TG 501 - Marker outputs, 5 sec to 1 ns. Sine-	and 10,000 μs, plus 10 ns sinewave. 184 - Marker outputs, 5 sec to 2 ns. Sine-
184	wave available at 5, 2, and 1 ns.	wave available at 50, 20, 10, 5,
	Trigger output - slaved to marker	and 2 ns. Separate trigger pulses
	output from 5 sec through 100 ns.	of 1 and .1 sec; 10, 1, and .1 ms;
	One time-mark can be generated at	10 and 1 μs. Marker amplifier pro-
	a time.	vides positive or negative time marks of 25 V min. Marker
		intervals of 1 and .1 sec; 10, 1,
		and .1 ms; 10 and 1 μ s.
2901	TG 501 - Marker outputs, 5 sec to 1 ns. Sine-	2901 - Marker outputs, 5 sec to 0.1 μs.
	wave available at 5, 2, and 1 ns.	Sinewave available to 50, 10,
	Trigger output - slaved to marker	and 5 ns. Separate trigger pulses, from 5 sec to 0.1 μ s.
S	output from 5 sec through 100 ns. One time-mark can be generated at	Multiple time-marks can be gene-
	a time.	rated simultaneously.
	Name and American	

NOTE: All TM 500 generator outputs are short-proof. All TM 500 plug-in instruments require TM 500-Series Power Module.